

Commission Regulation (EEC) No 2568/91 of 11 July 1991 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis

COMMISSION REGULATION (EEC) No 2568/91

of 11 July 1991

on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis

THE COMMISSION OF THE EUROPEAN COMMUNITIES,

Having regard to the Treaty establishing the European Economic Community,

Having regard to Council Regulation No 136/66/EEC of 22 September 1966 on the establishment of a common organization of the market in oils and fats⁽¹⁾, as last amended by Regulation (EEC) No 3577/90⁽²⁾, and in particular Article 35a thereof,

Whereas the Annex to Regulation No 136/66/EEC contains the descriptions and definitions of olive oil and olive-residue oil marketed within each Member State, in intra-Community trade and in trade with third countries;

Whereas, for the purpose of differentiating between the various types of oil, the physical and chemical characteristics of each of them and the organoleptic characteristics of virgin oil should be defined, in order to guarantee the purity and quality of the products concerned, without prejudice to other existing provisions;

Whereas the presence of the characteristics of the different types of oil should be determined uniformly throughout the Community; whereas, to that end, Community methods of chemical analysis and organoleptic evaluation should be established; whereas the use should be permitted, for a transitional period, of other methods of analysis applied in the Member States provided that where there is a difference in the results, those obtained using the common method will be decisive;

Whereas the definition of the physical and chemical characteristics of olive oil and of the methods of analysis entails the amendment of the additional notes to Chapter 15 of the combined nomenclature;

Whereas the method of evaluating the organoleptic characteristics of virgin oil includes the setting up of panels of selected and trained tasters; whereas the period necessary for establishing such a structure should therefore be fixed; whereas in view of the difficulties that some Member States will encounter in setting up panels of tasters, the use of panels in other Member States should be authorized;

Whereas, in order to ensure that the system of levies applicable to imports of olive residues functions correctly, a single method for the determination of the oil content of these products should be laid down;

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Whereas, in order not to harm trade, provision should be made for oil packaged prior to the entry into force of this Regulation to be disposed of during a limited period;

Whereas it is necessary to repeal Commission Regulation (EEC) No 1058/77⁽³⁾, as last amended by Regulation (EEC) No 1858/88⁽⁴⁾;

Whereas the Management Committee for Oils and Fats has not delivered an opinion within the time limit set by its chairman,

HAS ADOPTED THIS REGULATION:

F¹Article 1

1 Oils, the characteristics of which comply with those set out in points 1 and 2 of Annex I to this Regulation, shall be deemed to be virgin olive oils within the meaning of point 1(a) and (b) of the Annex to Regulation No 136/66/EEC.

2 Oil, the characteristics of which comply with those set out in point 3 of Annex I to this Regulation, shall be deemed to be lampante olive oil within the meaning of point 1(c) of the Annex to Regulation No 136/66/EEC.

3 Oil, the characteristics of which comply with those set out in point 4 of Annex I to this Regulation, shall be deemed to be refined olive oil within the meaning of point 2 of the Annex to Regulation No 136/66/EEC.

4 Oil, the characteristics of which comply with those set out in point 5 of Annex I to this Regulation, shall be deemed to be olive oil composed of refined olive oils and virgin olive oils within the meaning of point 3 of the Annex to Regulation No 136/66/EEC.

5 Oil, the characteristics of which comply with those set out in point 6 of Annex I to this Regulation, shall be deemed to be crude olive-pomace oil within the meaning of point 4 of the Annex to Regulation No 136/66/EEC.

6 Oil, the characteristics of which comply with those set out in point 7 of Annex I to this Regulation, shall be deemed to be refined olive-pomace oil within the meaning of point 5 of the Annex to Regulation No 136/66/EEC.

7 Oil, the characteristics of which comply with those set out in point 8 of Annex I to this Regulation, shall be deemed to be olive-pomace oil within the meaning of point 6 of the Annex to Regulation No 136/66/EEC.]

Textual Amendments

- F1** Substituted by [Commission Regulation \(EC\) No 1989/2003 of 6 November 2003 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis.](#)

F²Article 2

1 The characteristics of oils laid down in Annex I shall be determined in accordance with the following methods of analysis:

- a for the determination of the free fatty acids, expressed as the percentage of oleic acid, the method set out in Annex II;
- b for the determination of the peroxide index, the method set out in Annex III;
- c for determination of the wax content, the method set out in Annex IV;

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- d for the determination of the composition and content of sterols and triterpene dialcohols by capillary-column gas chromatography, the method set out in Annex V;
- e for the determination of the percentage of 2- glyceryl monopalmitate, the method set out in Annex VII;
- f for spectrophotometric analysis, the method set out in Annex IX;
- g for the determination of the fatty acid composition, the method set out in Annex X A and X B;
- h for the determination of the volatile halogenated solvents, the method set out in Annex XI;
- i for the evaluation of the organoleptic characteristics of virgin olive oil, the method set out in Annex XII;
- j for the determination of stigmastadienes, the method set out in Annex XVII;
- k for determining the content of triglycerides with ECN42, the method set out in Annex XVIII;
- l for the determination of the aliphatic alcohol content, the method set out in Annex XIX;
- m for the determination of the content of waxes, fatty acid methyl esters and fatty acid ethyl esters, the method set out in Annex XX.

In order to detect the presence of extraneous vegetable oils in olive oils, the method of analysis set out in Annex XXa shall be applied.

2 Verification by national authorities or their representatives of the organoleptic characteristics of virgin oils shall be effected by tasting panels approved by the Member States.

The organoleptic characteristics of an oil as referred to in the first subparagraph shall be deemed consonant with the category declared if a panel approved by the Member State confirms the grading.

Should the panel not confirm the category declared as regards the organoleptic characteristics, at the interested party's request, the national authorities or their representatives shall have carried out without delay two counter-assessments by other approved panels, at least one by a panel approved by the producer Member State concerned. The characteristics concerned shall be deemed consonant with the characteristics declared if at least two of the counter-assessments confirm the declared grade. If that is not the case, the interested party shall be responsible for the cost of the counter-assessments.

3 When the national authorities or their representatives verify the characteristics of the oil as provided for in paragraph 1, samples shall be taken in accordance with international standards EN ISO 661 on the preparation of test samples and EN ISO 5555 on sampling. However, notwithstanding point 6.8 of standard EN ISO 5555, in case of batches of such oils in immediate packaging, the sample shall be taken in accordance with Annex Ia to this Regulation. In case of bulk oils for which the sampling cannot be performed according to EN ISO 5555, the sampling shall be performed in accordance with instructions provided by the competent authority of the Member State.

Without prejudice to standard EN ISO 5555 and Chapter 6 of standard EN ISO 661, the samples taken shall be put in a dark place away from strong heat as quickly as possible and sent to the laboratory for analysis no later than the fifth working day after they are taken, otherwise the samples shall be kept in such a way that they will not be degraded or damaged during transport or storage before being sent to the laboratory.

4 For the purposes of the verification provided for in paragraph 3, the analyses referred to in Annexes II, III, IX, XII and XX and, where applicable, any counter-analyses required

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

under national law, shall be carried out before the minimum durability date in case of packaged products. In case of sampling of bulk oils, those analyses shall be carried out no later than the sixth month after the month in which the sample was taken.

No time limit shall apply to the other analyses provided for in this Regulation.

Unless the sample was taken less than two months before the minimum durability date, if the results of the analyses do not match the characteristics of the category of olive oil or olive-pomace oil declared, the party concerned shall be notified no later than one month before the end of the period laid down in the first subparagraph.

5 For the purpose of determining the characteristics of olive oils by the methods provided for in the first subparagraph of paragraph 1, the analysis results shall be directly compared with the limits laid down in this Regulation.]

Textual Amendments

- F2** Substituted by [Commission Implementing Regulation \(EU\) No 1348/2013 of 16 December 2013 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)

[^{F3}Article 2a

1 For the purpose of this Article, ‘olive oil marketed’ means total quantity of olive oil and olive pomace oil of a relevant Member State that is consumed in that Member State or exported from that Member State.

2 Member States shall ensure that conformity checks are carried out selectively, based on a risk analysis, and with appropriate frequency, so as to ensure that the olive oil marketed is consistent with the category declared.

3 The criteria to assess the risk may include:

- a the category of oil, the period of production, the price of oils in relation to other vegetable oils, the blending and packing operations, the storage facilities and conditions, the country of origin, the country of destination, the means of transport or the volume of the lot;
- b the position of the operators in the marketing chain, the volume and/or value marketed by them, the range of oil categories they market, the type of business carried out such as milling, storage, refining, blending, packaging or retail sale;
- c findings made during previous checks including the number and type of defects found, the usual quality of oils marketed, the performance of technical equipment used;
- d the reliability of operators’ quality assurance systems or self-checking systems related to the conformity to marketing standards;
- e the place where the check is carried out, in particular if it is the first point of entry into the Union, the last point of exit from the Union or the place where the oils are produced, packaged, loaded or sold to the final consumer;
- f any other information that might indicate a risk of non-compliance.

4 Member States shall lay down in advance:

- a the criteria for assessing the risk of non-conformity of lots;
- b on the basis of a risk analysis for each risk category, the minimum number of operators or lots and/or quantities which will be subject to a conformity check.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

At least one conformity check per thousand tonnes of olive oil marketed in the Member State shall be carried out per year.

- 5 Member States shall verify compliance by:
- a carrying out, in any order, the analyses set out in Annex I; or
 - b following the order set out in Annex Ib on the decision tree until one of the decisions appearing in the decision tree is reached.]

Textual Amendments

- F3** Substituted by [Commission Implementing Regulation \(EU\) No 299/2013 of 26 March 2013 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)

^{F4}Article 3

Textual Amendments

- F4** Deleted by [Commission Regulation \(EC\) No 796/2002 of 6 May 2002 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis and the additional notes in the Annex to Council Regulation \(EEC\) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.](#)

^{F4}Article 3a

Textual Amendments

- F4** Deleted by [Commission Regulation \(EC\) No 796/2002 of 6 May 2002 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis and the additional notes in the Annex to Council Regulation \(EEC\) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.](#)

^{F3}Article 3

Where it is found that an oil does not correspond to its category description, the Member State concerned shall, without prejudice to any other penalties, apply effective, proportionate and dissuasive penalties to be determined in the light of the seriousness of the irregularity detected.

Where checks reveal significant irregularities, Member States shall increase the frequency of checks in relation to marketing stage, oil category, origin, or other criteria.]

Textual Amendments

- F3** Substituted by [Commission Implementing Regulation \(EU\) No 299/2013 of 26 March 2013 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

[^{F5} Article 4

[^{F61} The Member States may approve assessment panels so that national authorities or their representatives can assess and verify organoleptic characteristics.

The terms of approval shall be set by Member States and ensure that:

- the requirements of Annex XII.4 are met,
- the panel head is given training recognised for this purpose by the Member State,
- continued approval depends on performance in annual checks arranged by the Member State.

Member States shall notify to the Commission a list of approved panels and the action taken under this paragraph.]

2 Where Member States encounter difficulties in setting up tasting panels in their territory, they may call on a tasting panel approved in another Member State.

3 Each Member State draw up a list of tasting panels set up by professional or inter-branch organizations in accordance with the conditions laid down in paragraph 1 and shall ensure that those conditions are complied with.]

Textual Amendments

- F5** Substituted by [Commission Regulation \(EEC\) No 3288/92 of 12 November 1992 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and of the relevant methods of analysis.](#)
- F6** Substituted by [Commission Regulation \(EC\) No 796/2002 of 6 May 2002 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis and the additional notes in the Annex to Council Regulation \(EEC\) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.](#)

^{F4} Article 5

Textual Amendments

- F4** Deleted by [Commission Regulation \(EC\) No 796/2002 of 6 May 2002 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis and the additional notes in the Annex to Council Regulation \(EEC\) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.](#)

Article 6

1 The oil content of oil cake and other residues resulting from the extraction of olive oil (CN codes 2306 90 11 and 2306 90 19) shall be determined using the method set out in Annex XV.

2 The oil content referred to in paragraph 1 shall be expressed as a percentage of the weight of oil to the weight of dry matter.

[^{F1} Article 7

The Community provisions concerning the presence of contaminants shall apply.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

As regards halogenated solvents, the limits for all categories of olive oils are as follows:

- maximum content of each halogenated solvent detected: 0,1 mg/kg,
- maximum total content of halogenated solvents detected: 0,2 mg/kg.]

Textual Amendments

- F1** Substituted by [Commission Regulation \(EC\) No 1989/2003 of 6 November 2003 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis.](#)

[^{F7} Article 7a

Natural or legal persons and groups of persons who hold olive oil and olive pomace oil from the extraction at the mill up to the bottling stage included, for whatever professional or commercial purposes, shall be required to keep entry and withdrawal registers for each category of such oils.

Member State shall ensure that the obligation laid down in the first paragraph is duly complied with.]

Textual Amendments

- F7** Inserted by [Commission Implementing Regulation \(EU\) No 299/2013 of 26 March 2013 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)

[^{F3} Article 8

1 Member States shall notify the Commission of the measures implementing this Regulation. They shall inform the Commission of any subsequent amendments.]

2 Member States shall send the Commission, at the beginning of each half-year, a statement of the analytical data relating to the tests carried out during the previous half-year.

The results shall be considered by the Management Committee for Oils and Fats in accordance with the procedure laid down in Article 39 of Regulation No 136/66/EEC.

[^{F33} The notifications referred to in this Regulation shall be made in accordance with Commission Regulation (EC) No 792/2009⁽⁵⁾.]

Textual Amendments

- F3** Substituted by [Commission Implementing Regulation \(EU\) No 299/2013 of 26 March 2013 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)

Article 9

Regulation (EEC) No 1058/77 is hereby repealed.

Article 10

1 This Regulation shall enter into force on the third day following its publication in the *Official Journal of the European Communities*.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

However, the method set out in Annex XII shall apply from [^{F8}1 November 1992], *except in so far as operations relating to the intervention system are concerned.*

[^{F9}That method shall not apply to virgin olive oil prepared for the market prior to 1 November 1992.]

2 This Regulation shall not apply to olive oil and olive-residue oil packaged before the entry into force of this Regulation and marketed up to 31 October 1992.

Textual Amendments

- F8** Substituted by [Commission Regulation \(EEC\) No 3682/91 of 17 December 1991 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)
- F9** Inserted by [Commission Regulation \(EEC\) No 3288/92 of 12 November 1992 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and of the relevant methods of analysis.](#)

This Regulation shall be binding in its entirety and directly applicable in all Member States.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

[^{X1}ANNEXES

Summary

Editorial Information

X1 Substituted by [Corrigendum to Commission Regulation \(EEC\) No 2568/91 of 11 July 1991 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis \(Official Journal of the European Communities No L 248 of 5 September 1991\)](#).

[^{F2} Annex I:	Olive oil characteristics]
[^{F2} Annex Ia:	Sampling of olive oil or olive-pomace oil delivered in immediate packaging]
[^{F2} Annexe Ib:	Decision tree for verifying whether an olive oil sample is consistent with the category declared]
Annex II:	[^{F10} Determination of free fatty acids, cold method]
Annex III:	Determination of the peroxide value
Annex IV:	[^{F11} Determination of wax content by capillary column gas-liquid chromatography]
[^{F2} Annex V:	Determination of the composition and content of sterols and triterpenes dialcohols by capillary-column gas chromatography]
[^{F12}	
Annex VII:	[^{F10} Determination of the percentage of 2-glyceryl monopalmitate]
[^{F13}	
[^{F3} Annex IX:	Spectrophotometric investigation in the ultraviolet]
Annex XA:	Analysis by gas chromatography of methyl esters of fatty acids
Annex XB:	Preparation of methyl esters of fatty acids
Annex XI:	Determination of the volatile halogenated solvents of olive oil
[^{F2} Annex XII:	The international olive council's method for the organoleptic assessment of virgin olive oil]
[^{F13}	
[^{F4}	
Annex XV:	Oil content of olive residue

*Status: Point in time view as at 01/03/2014.**Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)*

Annex XVI:	Determination of iodine value
Annex XVII:	Determination of stigmastadienes in vegetable oils
[^{F3} Annex XVIII:	Determination of the difference between actual and theoretical content of triacylglycerols with ECN 42]
Annex XIX:	Method for determining aliphatic alcohol content
[^{F14} Annex XX:	Method for the determination of the content of waxes, fatty acid methyl esters and fatty acid ethyl esters by capillary gas chromatography]
[^{F15} Annex XXa:	Method for the detection of extraneous oils in olive oils]
[^{F7} ANNEX XXI:	Results of conformity checks carried out on olive oils referred to in Article 8(2)]

Textual Amendments

- F10** Substituted by Commission Regulation (EC) No 702/2007 of 21 June 2007 amending Commission Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F11** Substituted by Commission Regulation (EEC) No 183/93 of 29 January 1993 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F12** Deleted by Commission Implementing Regulation (EU) No 1348/2013 of 16 December 2013 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F13** Deleted by Commission Regulation (EC) No 1989/2003 of 6 November 2003 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis.
- F14** Inserted by Commission Regulation (EU) No 61/2011 of 24 January 2011 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F15** Inserted by Commission Implementing Regulation (EU) No 1348/2013 of 16 December 2013 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

[F2] ANNEX I

OLIVE OIL CHARACTERISTICS

Category	Fatty acid ethyl esters (FAEEs) mg/kg (*)	Acidity (%) (*)	Peroxide index mEq O ₂ /kg (*)	Waxes mg/kg (**)	2-glycerin/monopalmitic (%)	Stigmasterol mg/kg	Differs ECN42 ^a and ECN42 ^b (HPLC) (theoretical calculation)	K ₂₆₈ or K ₂₇₀ (*)	Delta-K (*)	Organoleptic evaluation of defect (Md) (*)	Organoleptic evaluation of median (Mf) (*)	Fruity
1.	FAEEs ≤ 40 (2013-2014 olive oil) FAEEs ≤ 35 (2014-2015 crop year) FAEEs ≤ 30 (after 2015 crop years)	≤ 0,8	≤ 20	C ₄₂ + C ₄₄ + C ₄₆ ≤ 350	if total palmitic acid % ≤ 14 % ≤ 1,0 if total palmitic acid % > 14 %	≤ 0,05	≤ 2,50	≤ 0,22	≤ 0,01	Md = 0	Mf > 0	
2.	— Virgin olive oil	≤ 2,0	≤ 20	C ₄₂ + C ₄₄ + C ₄₆ ≤ 350	if total palmitic acid % ≤ 14 % ≤ 1,0 if total palmitic acid % > 14 %	≤ 0,05	≤ 2,60	≤ 0,25	≤ 0,01	Md ≤ 3,5	Mf > 0	

a Total isomers which could (or could not) be separated by capillary column.

b The olive oil has to be in conformity with the method set out in Annex XXa.

c This limit applies to olive oils produced as from 1st March 2014

d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.

e Or where the median of defect is above 3,5 or the median of defect is less than or equal to 3,5 and the fruity median is equal to 0.

f Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

3.	Lampante olive oil	> 2,0	—	$C_{40} + C_{42} + C_{44} + C_{46} \leq 300$ ^a	% > 14 % if total palmitic acid % ≤ 14 % ≤ 1,1 if total palmitic acid % > 14 %	$\leq 0,50$	$\leq 0,3$	—	—	—	Md > 3,5 ^e	—
4.	Refined olive oil	≤ 0,3	≤ 5	$C_{40} + C_{42} + C_{44} + C_{46} \leq 350$	% > 14 % if total palmitic acid % ≤ 14 % ≤ 1,1 if total palmitic acid % > 14 %	$\leq 0,3$	—	≤ 1,10	≤ 0,16	—	—	—
5.	Olive oil composed of refined and virgin olive oils	≤ 1,0	≤ 15	$C_{40} + C_{42} + C_{44} + C_{46} \leq 350$	% > 14 % if total palmitic acid % ≤ 14 % ≤ 1,0 if total	$\leq 0,3$	—	≤ 0,90	≤ 0,15	—	—	—

a Total isomers which could (or could not) be separated by capillary column.

b The olive oil has to be in conformity with the method set out in Annex XXa.

c This limit applies to olive oils produced as from 1st March 2014

d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.

e Or where the median of defect is above 3,5 or the median of defect is less than or equal to 3,5 and the fruity median is equal to 0.

f Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

				palmitic acid % > 14 %						
6.	Crude olive-pomace	—	—	$C_{40} + C_{42} \leq 1,4 + C_{44} + C_{46} > 350$	$\leq 0,6$	—	—	—	—	—
7.	Refined olive-pomace oil	$\leq 0,3$	≤ 5	$C_{40} + C_{42} \leq 1,4 + C_{44} + C_{46} > 350$	$\leq 0,5$	—	$\leq 2,00$	$\leq 0,20$	—	—
8.	Olive-pomace oil	$\leq 1,0$	≤ 15	$C_{40} + C_{42} \leq 1,4 + C_{44} + C_{46} > 350$	$\leq 0,5$	—	$\leq 1,70$	$\leq 0,18$	—	—

- a Total isomers which could (or could not) be separated by capillary column.
- b The olive oil has to be in conformity with the method set out in Annex XXa.
- c This limit applies to olive oils produced as from 1st March 2014
- d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.
- e Or where the median of defect is above 3,5 or the median of defect is less than or equal to 3,5 and the fruity median is equal to 0.
- f Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

Category	Fatty acid composition ^a						Total Sterols composition						Total Erythrodiol			
	Myristic (%)	Palmitic (%)	Stearic (%)	Oleic (%)	Linoleic (%)	Other (%)	trans isomers (%)	linolenic isomers (%)	cholesterol (%)	stigmasterol (%)	β-sitosterol (%) ^c	sterols (mg/kg)	erythrodiol (%)	uvaol (%) ^(**)		
1. Extra virgin olive oil	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,20$	$\leq 0,20$	$\leq 0,05$	$\leq 0,05$	$\leq 0,5$	$\leq 0,1$	$\leq 4,0$	$< \text{Camp}$	$\geq 3,0$	$\leq 0,5$	≥ 1000	$\leq 4,5$
2. Virgin olive oil	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,20$	$\leq 0,20$	$\leq 0,05$	$\leq 0,05$	$\leq 0,5$	$\leq 0,1$	$\leq 4,0$	$< \text{Camp}$	$\geq 3,0$	$\leq 0,5$	≥ 1000	$\leq 4,5$

- a Other fatty acids content (%): palmitic: 7,50-20,00; palmitoleic: 0,30-3,50; heptadecanoic: $\leq 0,30$; heptadecenoic: $\leq 0,30$; stearic: 0,50-5,00; oleic: 55,00-83,00; linoleic: 3,50-21,00.
- b See the Appendix to this Annex.
- c App β-sitosterol: Delta-5,23-stigmastadienol+chlerosterol+beta-sitosterol+sitostanol+delta-5-avenasterol+delta-5,24-stigmastadienol.
- d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.
- e Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

3.	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,20$	$\leq 0,20$	$\leq 0,10$	$\leq 0,10$	$\leq 0,5$	$\leq 0,1$	$\leq 4,0$	—	$\geq 93,0$	$\leq 0,5$	≥ 1000	$\leq 4,5^d$
	Lampante olive oil															
4.	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,20$	$\leq 0,20$	$\leq 0,20$	$\leq 0,30$	$\leq 0,5$	$\leq 0,1$	$\leq 4,0$	< Camp	$\geq 93,0$	$\leq 0,5$	≥ 1000	$\leq 4,5$
	Refined olive oil															
5.	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,20$	$\leq 0,20$	$\leq 0,20$	$\leq 0,30$	$\leq 0,5$	$\leq 0,1$	$\leq 4,0$	< Camp	$\geq 93,0$	$\leq 0,5$	≥ 1000	$\leq 4,5$
	Olive oil composed of refined and virgin olive oils															
6.	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,30$	$\leq 0,20$	$\leq 0,20$	$\leq 0,10$	$\leq 0,5$	$\leq 0,2$	$\leq 4,0$	—	$\geq 93,0$	$\leq 0,5$	≥ 2500	$> 4,5^e$
	Crude olive-pomace oil															
7.	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,30$	$\leq 0,20$	$\leq 0,40$	$\leq 0,35$	$\leq 0,5$	$\leq 0,2$	$\leq 4,0$	< Camp	$\geq 93,0$	$\leq 0,5$	≥ 800	$> 4,5$
	Refined olive-pomace oil															
8.	$\leq 0,03$	$\leq 1,00$	$\leq 0,60$	$\leq 0,40$	$\leq 0,30$	$\leq 0,20$	$\leq 0,40$	$\leq 0,35$	$\leq 0,5$	$\leq 0,2$	$\leq 4,0$	< Camp	$\geq 93,0$	$\leq 0,5$	≥ 600	$> 4,5$
	Olive pomace oil															

a Other fatty acids content (%): palmitic: 7,50-20,00; palmitoleic: 0,30-3,50; heptadecanoic: $\leq 0,30$; heptadecenoic: $\leq 0,30$; stearic: 0,50-5,00; oleic: 55,00-83,00; linoleic: 3,50-21,00.

b See the Appendix to this Annex.

c App β -sitosterol: Delta-5,23-stigmastadienol+chlerosterol+beta-sitosterol+sitostanol+delta-5-avenasterol+delta-5,24-stigmastadienol.

d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.

e Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

Notes:

- (a) The results of the analyses must be expressed to the same number of decimal places as used for each characteristic. The last digit must be increased by one unit if the following digit is greater than 4.
- (b) If just a single characteristic does not match the values stated, the category of an oil can be changed or the oil declared impure for the purposes of this Regulation.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- (c) If a characteristic is marked with an asterisk (*), referring to the quality of the oil, this means the following: - for lampante olive oil, it is possible for both the relevant limits to be different from the stated values at the same time, - for virgin olive oils, if at least one of these limits is different from the stated values, the category of the oil will be changed, although they will still be classified in one of the categories of virgin olive oil.
- (d) If a characteristic is marked with two asterisks (**), this means that for all types of olive-pomace oil, it is possible for both the relevant limits to be different from the stated values at the same time.


Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix

Decision tree

Campesterol decision tree for virgin and extra virgin olive oils:

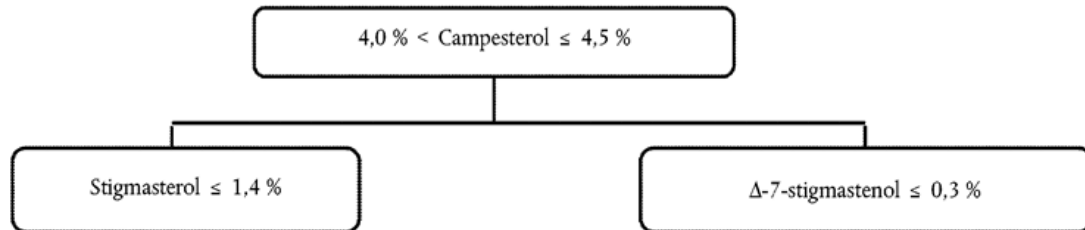

The other parameters shall comply with the limits fixed in this Regulation.

Delta-7-stigmasterol decision tree for:

— Extra virgin and virgin olive oils

The other parameters shall comply with the limits fixed in this Regulation.

— Olive-pomace oils (crude and refined]

[^{F2}ANNEX Ia

SAMPLING OF OLIVE OIL OR OLIVE-POMACE OIL DELIVERED IN IMMEDIATE PACKAGING

This method of sampling is applied to batches of olive oil or olive-pomace oil put up in immediate packaging. Different sampling methods apply, depending on whether the immediate packaging exceeds 5 litres or not.

‘Batch’ shall mean a set of sales units which are produced, manufactured and packed in circumstances such that the oil contained in each sales unit is considered to be homogenous in terms of all analytical characteristics. The individuation of a batch must be done in accordance with Directive 2011/91/EU of the European Parliament and of the Council⁽⁶⁾.

‘Increment’ shall mean the quantity of oil contained in an immediate package and taken from a random point of the batch.

1. CONTENT OF PRIMARY SAMPLE

1.1. Immediate packaging not exceeding 5 litres

‘Primary Sample’ for immediate packaging not exceeding 5 litres shall mean the number of increments taken from a batch and in agreement with Table 1.

TABLE 1

Primary sample minimum size must comprise the following

Where the immediate packaging has a capacity of	The primary sample must comprise the oil from
(a) 1 litre or more	(a) 1 immediate pack;

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

(b) less than 1 litre	(b) the minimum number of packs with a total capacity of at least 1,0 litre
-----------------------	---

The number of packs referred to in Table 1, which shall constitute a primary sample, can be increased by each Member State, according to their own needs (for example organoleptic assessment by a different laboratory from that which performed the chemical analyses, counter-analysis, etc.).

1.2. Immediate packaging exceeding 5 litres

‘Primary Sample’ for immediate packaging exceeding 5 litres shall mean a representative part of the total increments, obtained by a process of reduction and in agreement with Table 2. The primary sample must be composed of various examples.

‘Example’ of a primary sample shall mean each of the packages making up the primary sample.

TABLE 2

Minimum number of increments to be selected

Number of packages in the lot	Minimum number of increments to be selected
Up to 10	1
From ... 11 to 150	2
From ... 151 to 500	3
From ... 501 to 1 500	4
From ... 1 501 to 2 500	5
> 2 500 per 1 000 packages	1 extra increment

In order to reduce the volume of the sampling immediate packs, the content of the sampling increments is homogenised for the preparation of the primary sample. The portions of the different increments are poured into a common container for homogenisation by stirring, so that it will be best protected from air.

The content of the primary sample must be poured into a series of packages of the minimum capacity of 1,0 liter, each one of which constitutes an example of the primary sample.

The number of primary samples can be increased by each Member State, according to their own necessity (for example organoleptic assessment by a different laboratory from the one that performed the chemical analyses, counter-analysis, etc).

Each package must be filled in a way to minimise the air layer on top and then suitably closed and sealed to ensure the product is tamper-proof.

These examples must be labeled to ensure correct identification.

2. ANALYSES AND RESULTS

2.1. Each primary sample must be subdivided into laboratory samples, in accordance with point 2.5 of standard EN ISO 5555, and analysed according to the order shown in the decision tree set out in Annex Ib or in any other random order.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 2.2. Where all the results of the analyses comply with the characteristics of the category of oil declared, the whole batch is to be declared to comply.

If a single result of the analyses does not comply with the characteristics of the category of oil declared, the whole batch is to be declared non compliant.

3. VERIFICATION OF THE CATEGORY OF BATCH

- 3.1. In order to verify the batch category, the competent authority may increase the number of primary samples taken at different points of the batch according to the following table:

Table 3

Number of primary samples determined by the size of batch

Size of batch (litres)	Number of primary samples
Less than 7 500	2
From 7 500 to less than 25 000	3
From 25 000 to less than 75 000	4
From 75 000 to less than 125 000	5
Equal to and more than 125 000	6 + 1 each 50 000 litres more

Each increment constituting a primary sample must be taken from a continuous place in the batch; it is necessary to take note of the location of each primary sample and to identify it unambiguously.

The formation of each primary sample must be carried out according to the procedures referred to in points 1.1 and 1.2.

Each primary sample is then subjected to the analyses referred to in Article 2(1).

- 3.2. When one of the results of the analyses referred to in Article 2(1) of at least one primary sample does not comply with the characteristics of the declared category of oil, the whole sampling batch shall be declared non compliant.]

[^{F2}ANNEX Ib

DECISION TREE FOR VERIFYING WHETHER AN OLIVE OIL SAMPLE IS CONSISTENT WITH THE CATEGORY DECLARED

Table 1

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

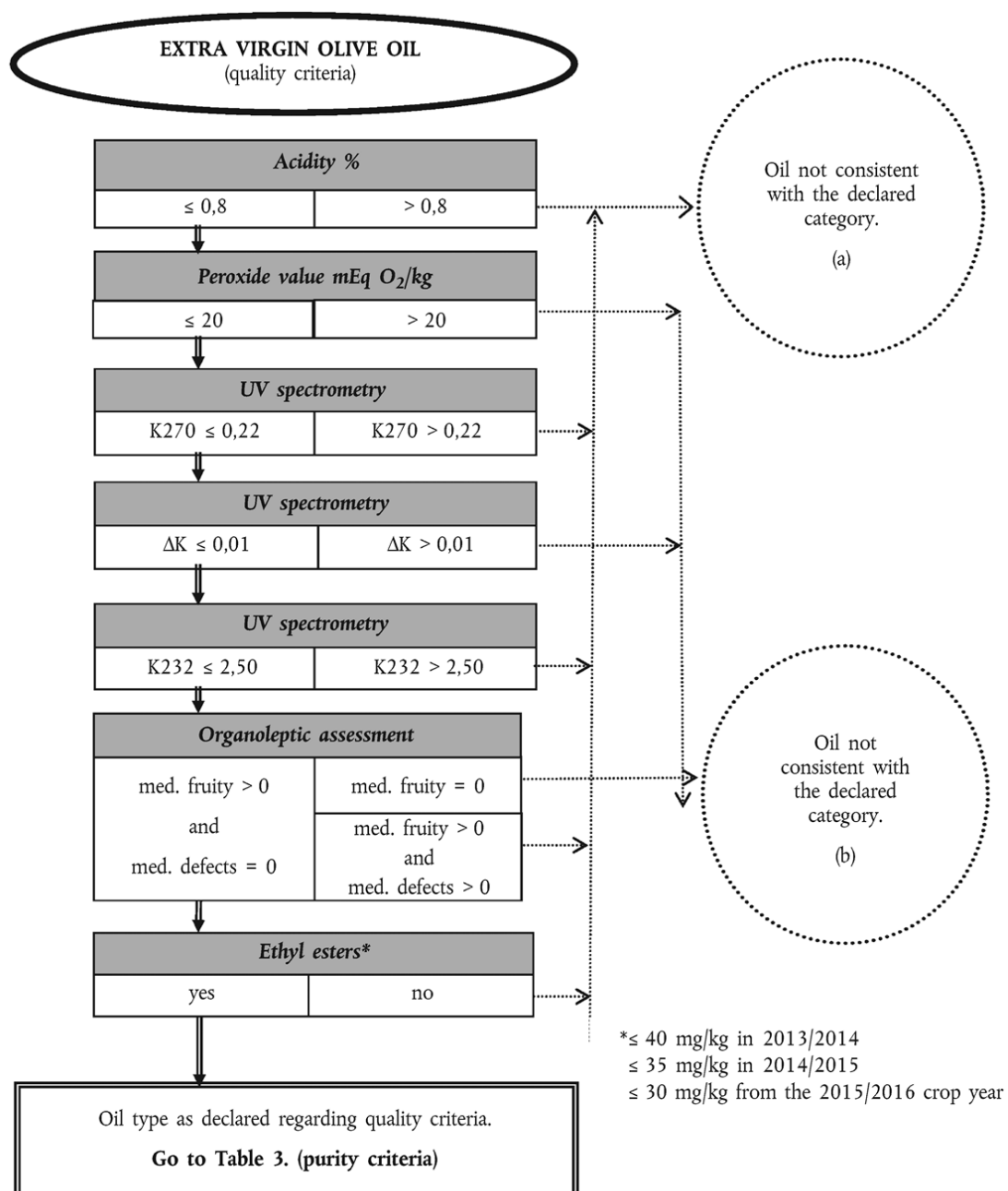


Table 2

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)


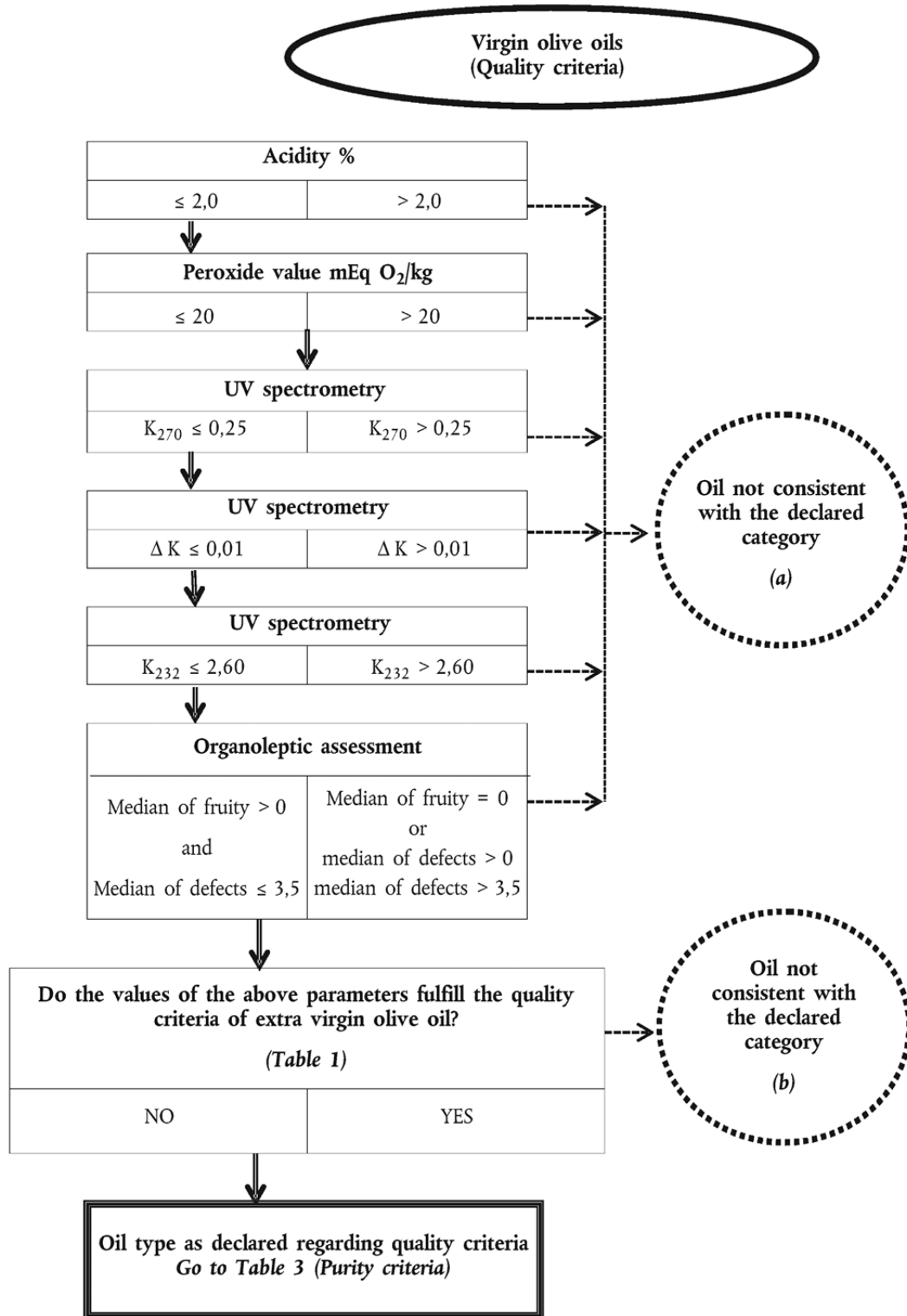


Table 3

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

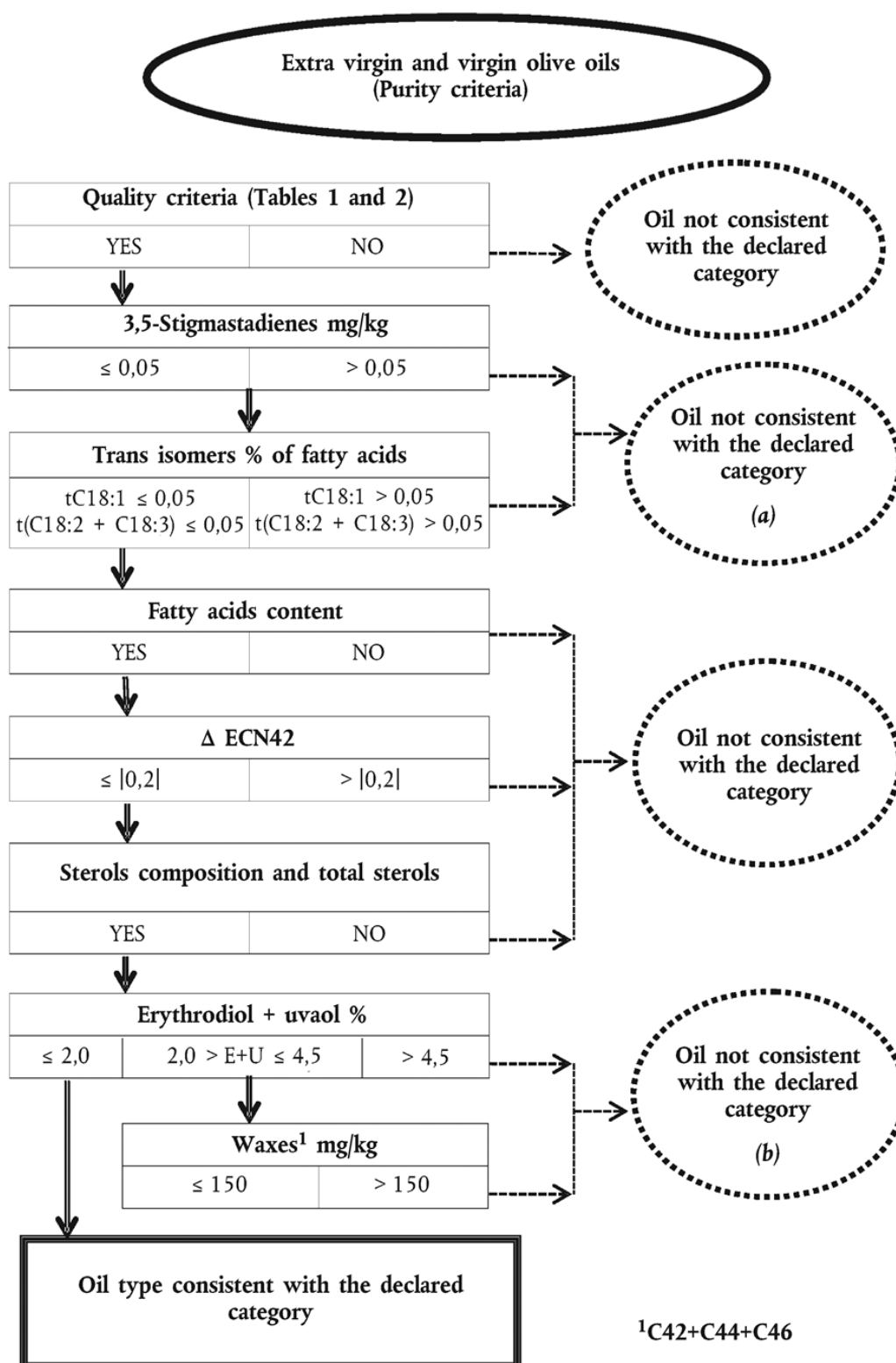

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix 1

TABLE OF EQUIVALENCE BETWEEN THE ANNEXES TO THIS REGULATION AND THE ANALYSES SPECIFIED IN THE DECISION TREE

—	Acidity	Annex II	Determination of free fatty acids, cold method
—	Peroxide value	Annex III	Determination of peroxide value
—	UV spectrometry	Annex IX	Spectrophotometric analysis
—	Organoleptic assessment	Annex XII	Organoleptic assessment of virgin olive oil
—	Ethyl esters	Annex XX	Method for the determination of the content of waxes, fatty acid methyl esters and fatty acid ethyl esters by capillary gas chromatography
—	3,5-Stigmastadienes	Annex XVII	Method of determining stigmastadienes in vegetable oils
—	Trans isomers of fatty acids	Annex X A and	Analysis by gas chromatography of methyl esters of fatty acids
		Annex X B	Preparation of methyl esters of fatty acids
—	Fatty acids content	Annex X A and	Analysis by gas chromatography of methyl esters of fatty acids
		Annex X B	Preparation of methyl esters of fatty acids
—	Δ ECN42	Annex XVIII	Determination of the composition of triglycerides with ECN42 (difference between the HPLC data and theoretical content)
—	Sterols composition and total sterols Erythrodiol and Uvaol	Annex V	Determination of the composition and content of sterols and triterpene dialcohols by capillary-column gas chromatography
—			
—	Waxes	Annex IV	Determination of wax content by capillary-column gas chromatography

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

—	Aliphatic alcohols	Annex XIX	Determination of aliphatic alcohols content by capillary-column gas chromatography
—	Saturated fatty acids in position 2	Annex VII	Determination of the percentage of 2-glyceryl monopalmitate]

ANNEX II

[^{F10}DETERMINATION OF FREE FATTY ACIDS, COLD METHOD]

1. DETERMINATION OF ACIDITY

The determination of free fatty acids in olive oils. The content of free fatty acids is expressed as acidity calculated conventionally.

1.1. Principle

A sample is dissolved in a mixture of solvents and the free fatty acids present titrated using an ethanolic solution of potassium hydroxide.

1.2. Reagents

All the reagents should be of recognized analytical quality and the water used either distilled or of equivalent purity.

1.2.1. [^{X1}Diethyl ether]; 95 % ethanol (v/v), mixture of equal parts by volume.

Note:

[^{X1}Diethyl ether] is highly inflammable and may form explosive peroxides. Special care should be taken in its use.

Neutralize precisely at the moment of use with the potassium hydroxide solution (1.2.2), with the addition of 0,3 ml of the phenolphthalein solution (1.2.3) per 100 ml of mixture.

Note:

If it is not possible to use [^{X1}diethyl ether], a mixture of solvents containing ethanol and toluene may be used. If necessary, ethanol may be replaced by propanol-2.

1.2.2. Potassium hydroxide, titrated ethanolic solution, c(KOH) about 0,1 mol/l or, if necessary, c(KOH) about 0,5 mol/l.

The exact concentration of the ethanolic solution of potassium hydroxide must be known and checked immediately prior to use. Use a solution prepared at least five days before use and decanted into a brown glass bottle with a rubber stopper. The solution should be colourless or straw coloured.

Note:

A stable colourless solution of potassium hydroxide may be prepared as follows. Bring to the boil 1 000 ml of ethanol with 8 g of potassium hydroxide and 0,5 g of aluminium shavings and continue boiling under reflux for one hour. Distill immediately. Dissolve in the distillate the required quantity of potassium hydroxide. Leave for several days and decant the clear supernatant liquid from the precipitate of potassium carbonate.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The solution may also be prepared without distillation as follows: to 1 000 ml of ethanol add 4 ml of aluminium butylate and leave the mixture for several days. Decant the supernatant liquid and dissolve the required quantity of potassium hydroxide. The solution is ready for use.

1.2.3. Phenolphthalein, 10 g/l solution in 95 to 96 % ethanol (v/v) or alkaline blue, (in the case of strongly coloured fats) 20 g/l solution in 95 to 96 % ethanol (v/v).

1.3. Apparatus

Usual laboratory equipment including:

1.3.1. analytical balance;

1.3.2. 250 ml conical flask;

1.3.3. 10 ml burette, graduated in 0,05 ml.

1.4. Procedure

1.4.1. Preparation of the specimen for testing

(Carry out the test on the filtered sample. Where moisture and impurities together are less than 1 %, use the specimen without further treatment; where they exceed 1 %, it should be filtered.)

1.4.2. Taking the sample

Take a sample depending on the presumed acid number in accordance with the following table:

Expected acid value	Mass of sample(g)	Weighing accuracy(g)
< 1	20	0,05
1 to 4	10	0,02
4 to 15	2,5	0,01
15 to 75	0,5	0,001
> 75	0,1	0,0002

Weigh the sample in the conical flask (1.3.2).

1.4.3. Determination

Dissolve the sample (1.4.2) in 50 to 150 ml of the previously neutralized mixture of diethyl [^{X1}ether] and ethanol (1.2.1).

Titrate while stirring with the 0,1 mol/l solution of potassium hydroxide (1.2.2) (see Note 2) until the indicator changes (the pink colour of the phenolphthalein persists for at least 10 seconds).

Note 1.

The titrated ethanolic solution of potassium hydroxide (1.2.2) may be replaced by an aqueous solution of potassium or sodium hydroxide provided that the volume of water introduced does not induce phase separation.

Note 2.

If the quantity of 0,1 mol/l potassium hydroxide solution required exceeds 10 ml, use the 0,5 mol/l solution.

Note 3.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

If the solution becomes cloudy during titration, add enough of the solvents (1.2.1) to give a clear solution.

1.5. Acidity: expressed as percentage of oleic acid

Acidity as a percentage by weight is equal to:

$$V \times c \times \frac{M}{1000} \times \frac{100}{m} = \frac{V \times c \times M}{10 \times m}$$

where:

- V = the volume of titrated potassium hydroxide solution used, in millilitres;
 c = the exact concentration in moles per litre of the titrated solution of potassium hydroxide used;
 M = the molar weight in grams per mole of the acid used to express the result (= 282);
 m = the weight in grams of the sample.

[^{X1}Take as the result, the arithmetic mean [^{F11}of two calculations] carried out.]

ANNEX III

DETERMINATION OF PEROXIDE VALUE

1. SCOPE

This Standard describes a method for the determination of the peroxide value of oils and fats.

2. FIELD OF APPLICATION

This Standard is applicable to animal and vegetable oils and fats.

3. DEFINITION

The peroxide value is the quantity of those substances in the sample, expressed in terms of milliequivalents of active oxygen per kilogram, which oxidize potassium iodide under the operating conditions described.

4. PRINCIPLE

Treatment of the test portion, in solution in acetic acid and chloroform, by a solution of potassium iodide. Titration of the liberated iodine with standardized sodium thiosulphate solution.

5. APPARATUS

All the equipment used shall be free from reducing or oxidizing substances.

Note:

Do not grease ground surfaces.

- 5.1. 3 ml glass scoop.
- 5.2. Flasks, with ground necks and stoppers, of about 250 ml capacity, dried beforehand and filled with a pure, dry inert gas (nitrogen or, preferably, carbon dioxide).
- 5.3. 25- or 50-ml burette, graduated in 0,1 ml.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

6. REAGENTS

- 6.1. Chloroform, analytical reagent quality, freed from oxygen by bubbling a current of pure, dry inert gas through it.
- 6.2. Glacial acetic acid, analytical reagent quality, freed from oxygen by bubbling a current of pure, [^{X1}dry inert gas] through it.
- 6.3. Potassium iodide, saturated aqueous solution, recently prepared, free from iodine and iodates.
- 6.4. Sodium thiosulphate, 0,01 or 0,002 [^{X1}Mol/L accurately] standardized aqueous solution, standardized just before use.
- 6.5. Starch solution, 10 g/l aqueous dispersion, recently prepared from natural soluble starch.

7. SAMPLE

Take care that the sample is taken and stored away from the light, kept cold and contained in completely filled glass containers, hermetically sealed with ground-glass or cork stoppers.

8. PROCEDURE

The test shall be carried out in diffuse daylight or in artificial light. Weigh in a glass scoop (5.1) or, failing this, in a flask (5.2), to the nearest 0,001 g, a mass of the sample in accordance with the following table, according to the expected peroxide value:

Expected peroxide value(meq)	Weight of test portion(g)
0 to 12	5,0 to 2,0
12 to 20	2,0 to 1,2
20 to 30	1,2 to 0,8
30 to 50	0,8 to 0,5
50 to 90	0,5 to 0,3

Unstopper a flask (5.2) and introduce the glass scoop containing the test portion. Add 10 ml of chloroform (6.1). Dissolve the test portion rapidly by stirring. Add 15 ml of acetic acid (6.2), then 1 ml of potassium iodide solution (6.3). Insert the stopper quickly, shake for one minute, and leave for exactly five minutes away from the light at a temperature from 15 to 25 °C.

Add about 75 ml of distilled water. Titrate the liberated iodine with the sodium thiosulphate solution [^{X1}(6.4) (0,002 Mol/L solution for expected values less than 12, and 0,01 Mol/L solution] for expected values above 12) shaking vigorously, using starch solution (6.5) as indicator.

Carry out two determinations on the same test sample.

Carry out simultaneously a blank test. If the result of the blank exceeds 0,05 ml of [^{X1}0,01 mol/L sodium] thiosulphate solution (6.4), replace the impure reagents.

9. EXPRESSION OF RESULTS

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The peroxide value (PV), expressed in milliequivalents of active oxygen per kilogram, is given by the formula:

$$PV = \frac{V \times T \times 1000}{m}$$

where:

- V = the number of ml of the standardized sodium thiosulphate solution (6.4) used for the test, corrected to take into account the blank test;
- T = [X¹the exact molarity] of the sodium thiosulphate solution (6.4) used;
- m = the weight in g, of the test portion.

Take as the result the arithmetic mean of the two determinations carried out.

[^{F10}ANNEX IV

DETERMINATION OF WAX CONTENT BY CAPILLARY COLUMN GAS CHROMATOGRAPHY

1. SUBJECT

This method describes a process for determining the wax content of olive oils. Waxes are separated according to the number of their carbon atoms. The method may be used in particular to distinguish between olive oil obtained by pressing and that obtained by extraction (olive-residue oil).

2. PRINCIPLE

Addition of a suitable internal standard to the fat or oil, then fractionation by chromatography on a hydrated silica gel column. Recovery under the test conditions of the fraction eluted first (the polarity of which is less than that of the triglycerides), then direct analysis by capillary column gas chromatography.

3. EQUIPMENT

- 3.1. 25 ml Erlenmeyer flask.
- 3.2. Glass column for gas chromatography, internal diameter 15,0 mm, length 30 to 40 cm, fitted with a stopcock.
- 3.3. Suitable gas chromatograph with a capillary column, equipped with a system for direct introduction into the column comprising the following:
 - 3.3.1. Thermostatic chamber for the columns, equipped with a temperature programmer.
 - 3.3.2. Cold injector for direct introduction into the column.
 - 3.3.3. Flame ionisation detector and converter-amplifier.
 - 3.3.4. Recorder-integrator capable of working with the converter-amplifier (3.3.3), rate of response no slower than 1 second, with variable paper speed. (It is also possible to use computerised systems that allow the acquisition of gas chromatography data via a PC.)
 - 3.3.5. Glass or fused silica capillary column 8 to 12 m long and with an internal diameter of 0,25 to 0,32 mm, with liquid phase, with a uniform film thickness between 0,10 and 0,30 µm. (There are liquid phases suitable for the purpose of type SE-52 or SE-54 available on the market.)

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 3.4. 10 µl microsyringe for on-column injection, equipped with a hardened needle.
- 3.5. Electro vibrator.
- 3.6. Rotary evaporator.
- 3.7. Muffle furnace.
- 3.8. Analytical balance with guaranteed precision of $\pm 0,1$ mg.
- 3.9. Normal laboratory glassware.

4. REAGENTS

- 4.1. Silica gel with a granule size of between 60 and 200 µm.

Place the gel in the furnace at 500 °C for at least four hours. After cooling, add 2 % water in relation to the quantity of sampled silica gel. Shake well to homogenise the slurry. Keep in darkness for at least 12 hours prior to use.

- 4.2. n-hexane, for chromatography.
- 4.3. Ethyl ether, for chromatography.
- 4.4. n-heptane, for chromatography.
- 4.5. Standard solution of lauryl arachidate, at 0,1 % (m/v) in hexane (internal standard). (It is also possible to use palmityl palmitate or myristyl stearate.)
 - 4.5.1. *Sudan 1 (1-phenyl-azo-2-naphthol).*
- 4.6. Carrier gas: hydrogen or helium, gas-chromatographic purity.
- 4.7. Auxiliary gases:
 - pure hydrogen for gas chromatography,
 - pure air for gas chromatography.

5. PROCEDURE

- 5.1. Preparation of the chromatographic column.

Suspend 15 g of silica gel (4.1) in the n-hexane (4.2) and introduce it into the column (3.2). Allow to settle spontaneously. Complete settling with the aid of an electro vibrator (3.5) to make the chromatographic layer more homogeneous. Percolate 30 ml of n-hexane to remove any impurities. Using the balance (3.8) weigh exactly 500 mg of the sample into the 25 ml Erlenmeyer flask (3.1), add the appropriate quantity of the internal standard (4.5) according to the presumed wax content. For example, add 0,1 mg of lauryl arachidate for olive oil, and 0,25 to 0,5 mg for olive-residue oil. Transfer the prepared sample to the chromatography column using two 2 ml portions of n-hexane (4.2).

Allow the solvent to flow away until it reaches 1 mm above the upper level of the absorbant then percolate a further 70 ml of n-hexane in order to eliminate the n-alkanes naturally present. Then start the chromatographic elution by collecting 180 ml of the mixture of n-hexane/ethyl ether (ratio 99:1), keeping a rate of flow of approximately 15 drops every 10 seconds. Elution of the sample must be carried out at a room temperature of 22 ± 4 °C.

NB:

- The n-hexane/ethyl ether mixture (99:1) must be prepared every day.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- For a visual check on the correct elution of the waxes 100 µl of 1 % Sudan in the elution mixture can be added to the sample in solution. Since the colourant has an intermediate retention, between waxes and triglycerides, when the coloration has reached the bottom of the column the elution should be suspended because all the waxes will have been eluted.

Dry the fraction thus obtained in a rotary evaporator (3.6.) until virtually all the solvent has been eliminated. Eliminate the final 2 ml of solvent with the aid of a weak current of nitrogen; then add 2-4 ml n-heptane.

5.2. Analysis by gas chromatography

5.2.1. Preparatory work

Fit the column to the gas chromatograph (3.3) by connecting the inlet port to the on-column system and the outlet port to the detector. Perform a general check on the GC apparatus (operation of gas circuits, detector and recorder efficiency, etc.).

If the column is being used for the first time it should be conditioned first. Pass a little gas through the column, then turn on the GC apparatus. Heat gradually until 350 °C is reached after about four hours. Maintain that temperature for at least two hours then regulate the apparatus to operating conditions (set gas flow, light flame, connect to the electronic recorder (3.3.4), set temperature of column chamber, detector, etc.) and record the signal at a sensitivity at least twice as high as that required for the analysis. The baseline must be linear, with no peaks of any kind, and must not show any deviation.

A negative straight-line drift indicates that the column connections are not tight; a positive drift that the column has not been sufficiently conditioned.

5.2.2. Choice of operating conditions

The operating conditions are generally as follows:

- column temperature:

	20 °C/ minute		5 °C/ minute		20 °C/ minute	
Initially 80 °C (1')	→	240 °C	→	325 °C (6')	→	340 °C (10')

- detector temperature: 350 °C;
- quantity of substance injected: 1 µl of the n-heptane solution (2-4 ml);
- carrier gas: helium or hydrogen at the correct linear velocity for the gas selected (see Appendix);
- instrument sensitivity: suitable for the following conditions:

The conditions may be modified according to the characteristics of the column and the GC apparatus to obtain separation of all the waxes and a satisfactory peak resolution (see figure); the internal standard C₃₂ retention time must be 18 ± 3 minutes. The most representative wax peak must be at least 60 % of the full scale.

The peak integration parameters must be established so as to obtain a correct evaluation of the areas of the peaks in question.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

NB: Given the high final temperature, a positive drift of no more than 10 % of the full scale is permitted.

5.3. Performance of the analysis

Sample 1 µl of the solution using the 10 µl microsyringe; withdraw the syringe plunger so that the needle is empty. Place the needle in the injector and after 1-2 seconds inject quickly; remove the needle slowly after about five seconds.

Record until the waxes are completely eluted.

The base line must always satisfy the required conditions.

5.4. Identification of peaks

Identification of the different peaks should be based on retention time by comparison with wax mixtures of known retention times analysed under the same conditions.

The figure is a chromatogram of the waxes of a virgin olive oil.

5.5. Evaluation of quantity

Calculate the areas of the peaks of the internal standard and the aliphatic esters of C₄₀ to C₄₆ using the integrator.

Calculate the wax content of each of the esters in mg/kg fat using the formula:

$$\text{ester, mg / kg} = \frac{A_x \times m_s \times 1000}{A_s \times m}$$

where:

A _x	=	area of each ester's peak, in square millimetres;
A _s	=	area of the internal standard's peak, in square millimetres;
m _s	=	mass of added internal standard, in milligrams;
m	=	mass of sample for analysis, in grams.

6. EXPRESSION OF RESULTS

Indicate the total of the contents of the various C₄₀ to C₄₆ waxes in mg/kg fat (ppm).


NB: The components to be quantified refer to the peaks with carbon pair numbers between esters C₄₀ and C₄₆, using the example of the olive oil wax chromatogram shown in the figure below. If ester C₄₆ appears twice, it is recommended that to identify it the fraction of the waxes of an olive-residue oil should be analysed where the C₄₆ peak is easy to identify because it is in the clear majority.

The results should be expressed to one decimal place.

Figure Chromatogram of the waxes of an olive oil⁽⁷⁾

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

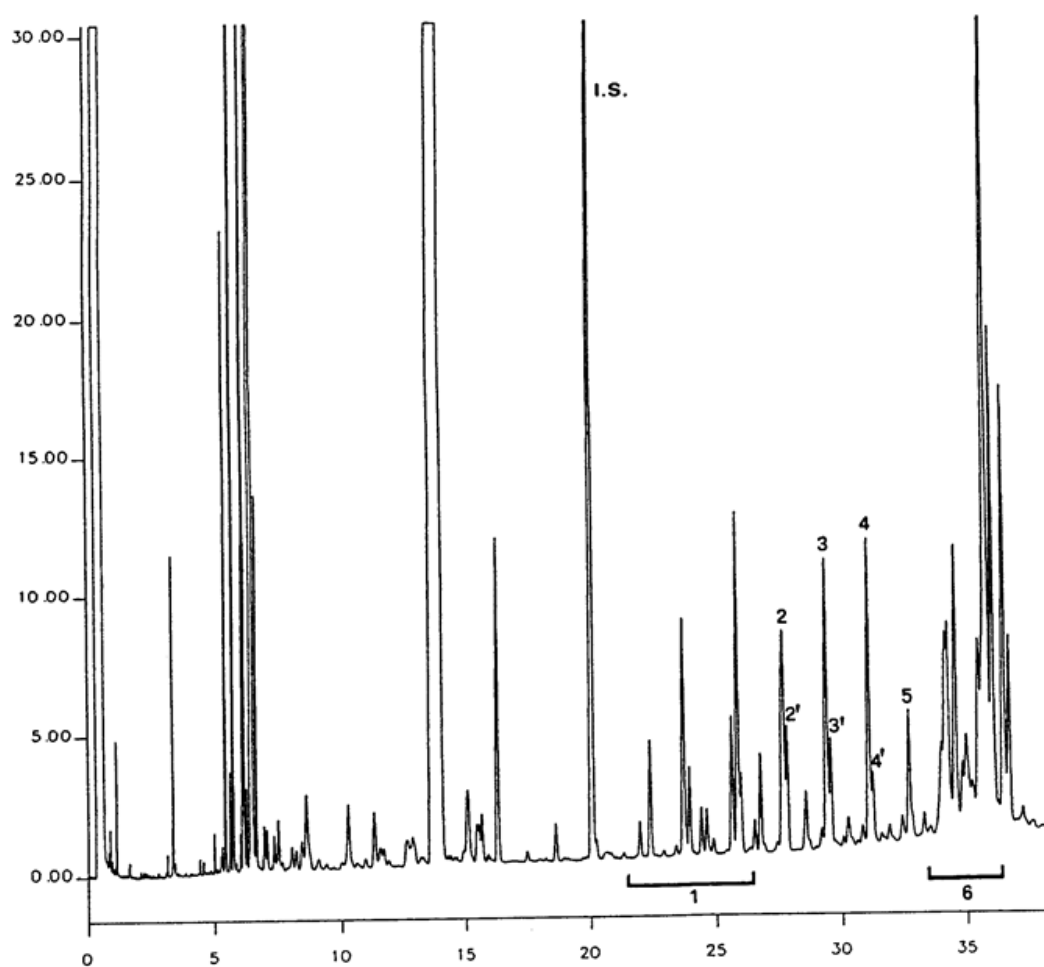

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix

Determination of the linear velocity of the gas

Inject 1-3 µl methane (or propane) into the GC apparatus after it has been regulated to normal operating conditions. Measure the time it takes for the gas to flow through the column from the time it is injected to the time the peak appears (t_M).

The linear velocity in cm/s is given by the formula L/t_M , where L is the length of the column in cm and t_M the time measured in seconds.]

[^{F2}ANNEX V

DETERMINATION OF THE COMPOSITION AND CONTENT OF STEROLS AND TRITERPENES DIALCOHOLS BY CAPILLARY-COLUMN GAS CHROMATOGRAPHY

1. SCOPE

The method describes a procedure for determining the individual and total sterols and triterpene dialcohols content of olive oils and olive-pomace oils.

2. PRINCIPLE

The oil, with added α -cholestanol as an internal standard, is saponified with potassium hydroxide in ethanolic solution and the unsaponifiable matter is then extracted with ethyl ether.

The sterols and triterpene dialcohols fraction is separated from the unsaponifiable matter by thin-layer chromatography on a basic silica gel plate. The fractions recovered from the silica gel are transformed into trimethylsilyl ethers and then analysed by capillary column gas chromatography.

3. APPARATUS

The usual laboratory equipment and in particular the following:

- 3.1. 250 ml flask fitted with a reflux condenser with ground-glass joints.
- 3.2. 500 ml separating funnel.
- 3.3. 250 ml flasks.
- 3.4. Complete apparatus for analysis by thin-layer chromatography using 20 x 20 cm glass plates.
- 3.5. Ultraviolet lamp with a wavelength of 254 or 366 nm.
- 3.6. 100 µl and 500 µl microsyringes.
- 3.7. Cylindrical filter funnel with a G3 porous septum (porosity 15-40 µm) of diameter approximately 2 cm and a depth of 5 cm, suitable for filtration under vacuum with male ground-glass joint.
- 3.8. 50 ml vacuum conical flask with ground-glass female joint, which can be fitted to the filter funnel (point 3.7).
- 3.9. 10 ml test tube with a tapering bottom and a sealing glass stopper.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 3.10. Gas chromatograph suitable for use with a capillary column with split injection system, consisting of:
- 3.10.1. A thermostatic chamber for columns capable of maintaining the desired temperature with an accuracy of $\pm 1^{\circ}\text{C}$;
 - 3.10.2. A temperature-adjustable injection unit with a persilanised glass vaporising element and split system;
 - 3.10.3. A flame ionisation detector (FID);
 - 3.10.4. Data acquisition system suitable for use with the FID detector (point 3.10.3.), capable of manual integration.
- 3.11. Fused-silica capillary column of length 20 to 30 m, internal diameter 0,25 to 0,32 mm, coated with 5 % diphenyl - 95 % dimethylpolysiloxane (SE-52 or SE-54 stationary phase or equivalent), to a uniform thickness between 0,10 and 0,30 μm .
- 3.12. Microsyringe, of 10 ml capacity, for gas chromatography, with cemented needle suitable for split injection.
- 3.13. Calcium dichloride desiccator
4. REAGENTS
- 4.1. Potassium hydroxide minimum titre 85 %.
- 4.2. Potassium hydroxide ethanolic solution, approximately 2 N.
- Dissolve 130 g of potassium hydroxide (point 4.1) with cooling in 200 ml of distilled water and then make up to one litre with ethanol (point 4.10). Keep the solution in well-stoppered dark glass bottles and stored maximum two days.
- 4.3. Ethyl ether, for analysis quality.
- 4.4. Potassium hydroxide ethanolic solution, approximately 0,2 N.
- Dissolve 13 g of potassium hydroxide (point 4.1) in 20 ml of distilled water and make up to one litre with ethanol (point 4.10).
- 4.5. Anhydrous sodium sulphate, for analysis quality.
- 4.6. Glass plates (20x20 cm) coated with silica gel, without fluorescence indicator, thickness 0,25 mm (commercially available ready for use).
- 4.7. Toluene, for chromatography quality.
- 4.8. Acetone, for chromatography quality.
- 4.9. n-Hexane, for chromatography quality.
- 4.10. Ethyl ether, for chromatography quality.
- 4.11. Ethanol of analytical quality.
- 4.12. Ethyl acetate of analytical quality.
- 4.13. Reference solution for thin-layer chromatography: cholesterol or phytosterols, and erythrodiol 5 % solution in ethyl acetate (point 4.11).

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 4.14. 2,7-dichlorofluorescein, 0,2 % in ethanolic solution. Make slightly basic by adding a few drops of 2 N alcoholic potassium hydroxide solution (point 4.2).
- 4.15. Anhydrous pyridine, for chromatography quality (see Note 5).
- 4.16. Hexamethyl disilazane of analytical quality.
- 4.17. Trimethylchlorosilane of analytical quality.
- 4.18. Sample solutions of sterol trimethylsilyl ethers.

To be prepared at the time of use from sterols and erythrodiol obtained from oils containing them.

- 4.19. α -cholestanol, purity more than 99 % (purity must be checked by GC analysis).
- 4.20. α -cholestanol internal standard solution, 0,2 % solution (m/V) in ethyl acetate (point 4.11).
- 4.21. Phenolphthalein solution, 10 g/L in ethanol (point 4.10).
- 4.22. Carrier gas: hydrogen or helium, gas-chromatographic purity.
- 4.23. Auxiliary gases: hydrogen, helium, nitrogen and air, of gas-chromatographic purity.
- 4.24. n-Hexane (point 4.9)/ethyl ether (point 4.10) mixture 65:35 (V/V).
- 4.25. Silylation reagent, consisting of a 9:3:1 (V/V/V) mixture of pyridine/hexamethyl disilazane/trimethylchlorosilane.

5. PROCEDURE

- 5.1. Preparation of the unsaponifiable matter.
 - 5.1.1. Using a 500 μ l microsyringe (point 3.6) introduce into the 250 ml flask (point 3.1) a volume of the α -cholestanol internal standard solution (point 4.20) containing an amount of cholestanol corresponding to approximately 10 % of the sterol content of the sample. For example, for 5 g of olive oil sample add 500 μ l of the α -cholestanol solution (point 4.20) and 1 500 μ l for olive-pomace oil. Evaporate until dryness with a gentle current of nitrogen in a warm water bath, after cooling the flask, weigh $5 \pm 0,01$ g of the dry filtered sample into the same flask.

Note 1: Animal or vegetable oils and fats containing appreciable quantities of cholesterol may show a peak having a retention time close to cholestanol. If this occurs, the sterol fraction will have to be analysed in duplicate with and without internal standard.

- 5.1.2. Add 50 ml of 2 N ethanolic potassium hydroxide solution (point 4.2) and some pumice, fit the reflux condenser and heat to gentle boiling until saponification takes place (the solution becomes clear). Continue heating for a further 20 minutes, then add 50 ml of distilled water from the top of the condenser, detach the condenser and cool the flask to approximately 30 °C.
- 5.1.3. Transfer the contents of the flask quantitatively into a 500 ml separating funnel (point 3.2) using several portions of distilled water (50 ml). Add approximately 80 ml of ethyl ether (point 4.10), shake vigorously for approximately 60 seconds, periodically releasing the pressure by inverting the separating funnel and opening the stopcock. Allow standing until there is complete separation of the two phases (Note 2).

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Then draw off the soap solution as completely as possible into a second separating funnel. Perform two further extractions on the water-alcohol phase in the same way using 60 to 70 ml of ethyl ether (point 4.10).

Note 2: Any emulsion can be destroyed by adding small quantities of ethanol (point 4.11).

5.1.4. Combine the three ether extracts in one separating funnel containing 50 ml of water. Continue to wash with water (50 ml) until the wash water no longer gives a pink colour on the addition of a drop of phenolphthalein solution (point 4.21).

When the wash water has been removed, filter on anhydrous sodium sulphate (point 4.5) into a previously weighed 250 ml flask, washing the funnel and filter with small quantities of ethyl ether (point 4.10).

5.1.5. Evaporate the solvent by distillation on a rotary evaporator at 30 °C under vacuum. Add 5 ml of acetone and remove the volatile solvent completely in a gentle current of air. Dry the residue in the oven at 103±2 °C for 15 min. Cool in desiccators and weigh to the nearest 0,1 mg.

5.2. Separation of the sterol and triterpene dialcohols fraction (erythrodiol + uvaol)

5.2.1. Preparation of the basic thin layer chromatography plates. Immerse the silica gel plates (point 4.6) about 4 cm in the 0,2 N ethanolic potassium hydroxide solution (point 4.5) for 10 seconds, then allow to dry in a fume cupboard for two hours and finally place in an oven at 100 °C for one hour.

Remove from the oven and keep in a calcium chloride desiccator (point 3.13) until required for use (plates treated in this way must be used within 15 days).

Note 3: When basic silica gel plates are used to separate the sterol fraction there is no need to treat the unsaponifiable fraction with alumina. In this way all compounds of an acid nature (fatty acids and others) are retained on the spotting line and the sterols band is clearly separated from the aliphatic and triterpene alcohols band.

5.2.2. Place hexane/ethyl ether mixture (point 4.24) (Note 4) into the development chamber, to a depth of approximately 1 cm. Close the chamber with the appropriate cover and leave thus for at least half an hour, in a cool place, so that liquid-vapour equilibrium is established strips of filter paper dipping into the eluent may be placed on the internal surfaces of the chamber. This reduces developing time by approximately one-third and brings about more uniform and regular elution of the components.

Note 4: The developing mixture needs to be replaced for every test, in order to achieve perfectly reproducible elution conditions, alternative solvent 50:50 (V/V) n-hexane/ethyl ether may be used.

5.2.3. Prepare an approximately 5 % solution of the unsaponifiable (point 5.1.5) in ethyl acetate (point 4.12) and, using the 100 µl microsyringe, depose 0,3 ml of the solution on a narrow and uniform streak on the lower end (2 cm) of the chromatographic plate (point 5.2.1). In line with the streak, place 2 to 3 µl of the material reference solution (point 4.13), so that the sterol and triterpene dialcohols band can be identified after developing.

5.2.4. Place the plate in the developing chamber prepared as specified in point 5.2.2. The ambient temperature should be maintained between 15 and 20 °C (Note 5). Immediately close the chamber with the cover and allow eluting until the solvent front reaches approximately 1 cm from the upper edge of the plate. Remove the plate from

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

the developing chamber and evaporate the solvent in a flow of hot air or by leaving the plate for a short while, under a hood.

Note 5: Higher temperature could worsen the separation.

5.2.5. Spray the plate lightly and uniformly with the 2,7-dichlorofluorescein solution (point 4.14) and then leave to dry. When the plate is observed under ultraviolet light, the sterols and triterpene dialcohols bands can be identified through being aligned with the spots obtained from the reference solution (point 4.13). Mark the limits of the bands along the edges of the fluorescence with a black pencil (see TLC plate figure 3).

5.2.6. By using a metal spatula, scrape off the silica gel of the marked area. Place the finely comminuted material removed into the filter funnel (point 3.7). Add 10 ml of hot ethyl acetate (point 4.12), mix carefully with the metal spatula and filter under vacuum, collecting the filtrate in the conical flask (point 3.8.) attached to the filter funnel.

Wash the residue in the flask three times with ethyl ether (point 4.3) (approximately 10 ml each time), collecting the filtrate in the same flask attached to the funnel, evaporate the filtrate to a volume of 4 to 5 ml, transfer the residual solution to the previously weighed 10 ml test tube (point 3.9), evaporate to dryness by mild heating, in a gentle flow of nitrogen, make up again using a few drops of acetone (point 4.8), evaporate again to dryness,

The residue contained in the test tube must consist of the sterol and triterpene dialcohols fractions.

5.3. Preparation of the trimethylsilyl ethers.

5.3.1. Add the silylation reagent (point 4.25) (Note 6), in the ratio of 50 µl for every milligram of sterols and triterpene dialcohols, in the test tube containing the sterol and triterpene fraction, avoiding any uptake of moisture (Note 7).

Note 6: Ready for use solutions are available commercially. Other silylation reagents, such as, for example, bistrimethylsilyl trifluor acetamide + 1 % trimethylchlorosilane, which has to be diluted with an equal volume of anhydrous pyridine, are also available.

Pyridine can be replaced by the same amount of acetonitrile.

5.3.2. Stopper the test tube, shake carefully (without overturning) until the compounds are completely dissolved. Leave to stand for at least 15 minutes at ambient temperature and then centrifuge for a few minutes. The clear solution is ready for gas chromatographic analysis.

Note 7: The slight opalescence, which may form, is normal and does not cause any anomaly. The formation of a white flock or the appearance of a pink colour is indicative of the presence of moisture or deterioration of the reagent. If this occurs, the test must be repeated (only if hexamethyldisilazane/trimethylchlorosilane is used).

5.4. Gas chromatographic analysis.

5.4.1. Preliminary operations, capillary column conditioning.

5.4.1.1. Fit the column (point 3.11) in the gas chromatograph, by attaching the inlet end to the split injector and the outlet end to the detector.

Carry out general checks on the gas chromatograph unit (leaks from the gas circuits, detector efficiency, efficiency of the splitting system and recording system, etc.).

5.4.1.2. If the column is being used for the first time, it is recommended that it be subjected to conditioning: passing a gentle flow of gas through the column itself, then switch

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

on the gas chromatography unit and begin a gradual heating, up to a temperature of at least 20 °C above the operating temperature (Note 8). Hold this temperature for at least two hours, then place the entire unit in operating mode (adjustment of gas flows and splitting, ignition of the flame, connection with the computing system, adjustment of the column, detector and injector temperature, etc.) and then record the signal with a sensitivity at least two times greater than that one intended for the analysis. The course of the base line must be linear, without peaks of any kind, and must not show drift.

A negative straight-line drift indicates leakage from the column connections; a positive drift indicates inadequate conditioning of the column.

Note 8: The conditioning temperature must always be at least 20 °C less than the maximum temperature specified for the stationary phase used.

5.4.2. Choice of operating conditions.

5.4.2.1. The operating conditions are as follows:

- Column temperature: 260 ± 5 °C;
- Injector temperature: 280-300 °C;
- Detector temperature: 280-300 °C;
- Linear velocity of the carrier gas: helium 20 to 35 cm/s; hydrogen 30 to 50 cm/s;
- Splitting ratio: from 1:50 to 1:100;
- Instrument sensitivity: from 4 to 16 times the minimum attenuation;
- Recording sensitivity: 1 to 2 mV full scale;
- Amount of substance injected: 0,5 to 1 µl of TMSE solution.

These conditions may be changed according to the characteristics of the column and gas chromatograph, so as to obtain chromatograms, which meet the following requirements:

- The retention time for the β-sitosterol peak should be at 20 ± 5 min;
- The campesterol peak should be: for olive oil (mean content 3 %) 20 ± 5 % of full scale; for soybean oil (average content 20 %) 80 ± 10 % of full scale;
- All the present sterols must be separated. In addition to being separated the peaks, they must also be completely resolved, i.e. the peak trace should return to the base line before leaving for the next peak. Incomplete resolution is, however, tolerated, provided that the peak at RRT 1,02 (Sitostanol) can be quantified using the perpendicular.

5.4.3. Analytical procedure

- ##### 5.4.3.1.
- By using the 10 µl microsyringe, take 1 µl of hexane, draw in 0,5 µl of air and then 0,5 to 1 µl of the sample solution. Raise the plunger of the syringe further, so the needle is emptied. Push the needle through the membrane of the injector and after one to two seconds, inject rapidly, and then slowly remove the needle after around five seconds.

An automatic injector can be used as well.

- ##### 5.4.3.2.
- Carry out the recording until the TMSE of the present triterpene dialcohols are completely eluted. The base line must continue to meet the requirements (point 5.4.1.2).

5.4.4. Peak identification

Identify individual peaks on the basis of retention times and by comparison with the mixture of sterol and triterpene dialcohols TMSE, analysed under the same conditions (see Appendix).

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The sterols and triterpene dialcohols are eluted in the following order: cholesterol, brassicasterol, ergosterol, 24-methylen-cholesterol, campesterol, campestanol, stigmasterol, Δ^7 -campesterol, $\Delta^5,23$ -stigmastadienol, clerosterol, β -sistosterol, sitostanol, Δ^5 -avenasterol, $\Delta^5,24$ -stigmastadienol, Δ^7 -stigmastenol, Δ^7 -avenasterol, erythrodiol and uvaol.

The retention times for β -sitosterol, for SE-52 and SE-54 columns, are shown in Table 1.

Figures 1 and 2 show typical chromatograms for some oils.

5.4.5. Quantitative evaluation.

5.4.5.1. Calculate the areas of the α -cholestanol and the sterol and triterpene dialcohols peaks by using the computing system. Ignore peaks for any compound which are not included (ergosterol must not be calculated) among those listed in Table 1. The response factor for α -cholestanol should be considered equal to 1.

5.4.5.2. Calculate the concentration of each individual sterol, in mg/kg of fatty material, as follows:

$$\text{sterol } x = \frac{A_x \times m_s \times 1000}{A_s \times m}$$

where:

- A_x = peak area for sterol x, in computing system counts;
 A_s = area of the α -cholestanol peak, in computing system counts;
 m_s = mass of added α -cholestanol, in milligrams;
 m = mass of the sample used for determination, in grams.

6. EXPRESSION OF THE RESULTS

6.1. Report individual sterol concentrations as mg/kg of fatty material and their sum as 'total sterols'.

The composition of each of the individual sterols and of the erythrodiol and uvaol should be expressed to one decimal point.

Total sterol composition must be expressed without any decimal point.

6.2. Calculate the percentage of each individual sterol from the ratio of the relevant peak area to the total peak area for sterols and erythrodiol and uvaol:

$$\text{sterol } x = \frac{A_x}{\Sigma A} \times 100$$

where:

- A_x = peak area for x;
 ΣA = total peak area for sterols;

6.3. Apparent β -sitosterol: Δ^5 -23-stigmastadienol + clerosterol + β -sitosterol + sitostanol + Δ^5 -avenasterol + Δ^5 -24-stigmastadienol.

6.4. Calculate the percentage of erythrodiol and uvaol:

$$\text{Erythrodiol} + \text{Uvaol} = \frac{E_r + U_v}{E_r + U_v + \Sigma A} \times 100$$

where

- ΣA = sum area for sterol in computing system counts;
 E_r = area of Erythrodiol in computing system counts;

Status: Point in time view as at 01/03/2014.

Changes to legislation: *There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)*

Uv = area of Uvaol in computing system counts;

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix

Determination of the linear speed of the gas

With the gas chromatograph set to normal operating conditions, inject 1 to 3 μl of methane (or propane) and measure the time taken by the gas to pass through the column, from the time of injection to the time at which the peak appears (tM).

The linear speed in cm/s is given by L/tM , where L is the length of the column in centimetres and tM is the measured time, in seconds.

TABLE 1

Relative retention times for sterols

Peak	Identification		Relative retention time	
			SE 54 column	SE 52 column
1	Cholesterol	Δ -5-cholesten-3 β -ol	0,67	0,63
2	Cholestanol	5 α -cholestan-3 β -ol	0,68	0,64
3	Brassicasterol	[24S]-24-methyl- Δ -5,22-cholestadien-3 β -ol	0,73	0,71
*	Ergosterol	[24S] 24-methyl- Δ -5,7,22-cholestatrien-3 β -ol	0,78	0,76
4	24-methylene-cholesterol	24-methylene- Δ -5,24-cholestadien-3 β -ol	0,82	0,8
5	Campesterol	(24R)-24-methyl- Δ -5-cholesten-3 β -ol	0,83	0,81
6	Campestanol	(24R)-24-methyl-cholestan-3 β -ol	0,85	0,82
7	Stigmasterol	(24S)-24-ethyl- Δ -5,22-cholestadien-3 β -ol	0,88	0,87
8	Δ -7-campesterol	(24R)-24-methyl- Δ -7-cholesten-3 β -ol	0,93	0,92
9	Δ -5,23-stigmastadienol	(24R,S)-24-ethyl- Δ -5,23-cholestadien-3 β -ol	0,95	0,95

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

10	Clerosterol	(24S)-24-ethyl- Δ -5,25-cholestadien-3 β -ol	0,96	0,96
11	β -sitosterol	(24R)-24-ethyl- Δ -5-cholesten-3 β -ol	1,0	1,0
12	Sitostanol	24-ethyl-cholestan-3 β -ol	1,02	1,02
13	Δ -5-avenasterol	(24Z)-24-ethylidene- Δ -cholesten-3 β -ol	1,03	1,03
14	Δ -5-24-stigmastadienol	(24R,S)-24-ethyl- Δ -5,24-cholestadien-3 β -ol	1,08	1,08
15	Δ -7-stigmastenol	(24R,S)-24-ethyl- Δ -7-cholesten-3 β -ol	1,12	1,12
16	Δ -7-avenasterol	(24Z)-24-ethylidene- Δ -7-cholesten-3 β -ol	1,16	1,16
17	Erythrodiol	5 α olean-12en-3 β 28diol	1,41	1,41
18	Uvaol	Δ 12-ursen-3 β 28diol	1,52	1,52

Figure 1

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Gas chromatogram of the sterol and triterpene dialchols fraction of a lampante olive oil (spiked with internal standard)

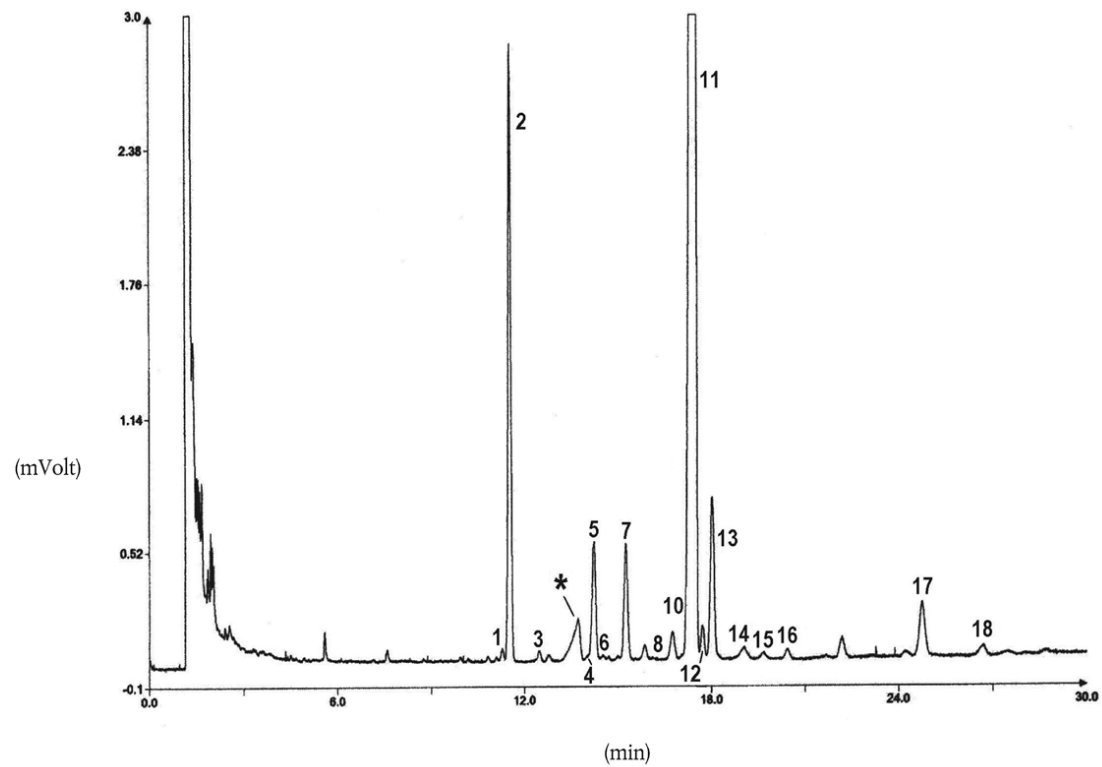


Figure 2

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Gas chromatogram of the sterol and triterpene dialchols fraction of a refined olive oil (spiked with internal standard)


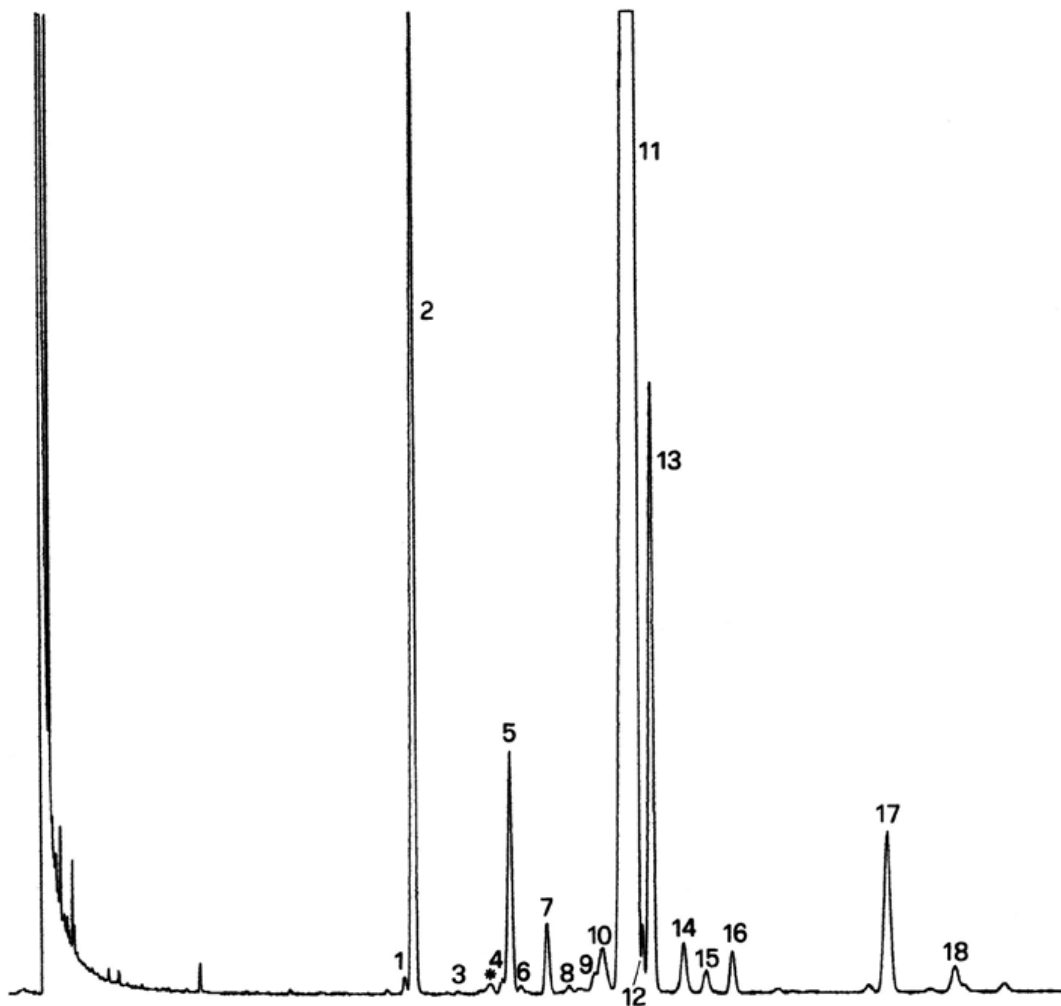


Figure 3

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

TLC plate olive-pomace oil with the zone that must be scraped for sterols and triterpenic dialcohols determination]

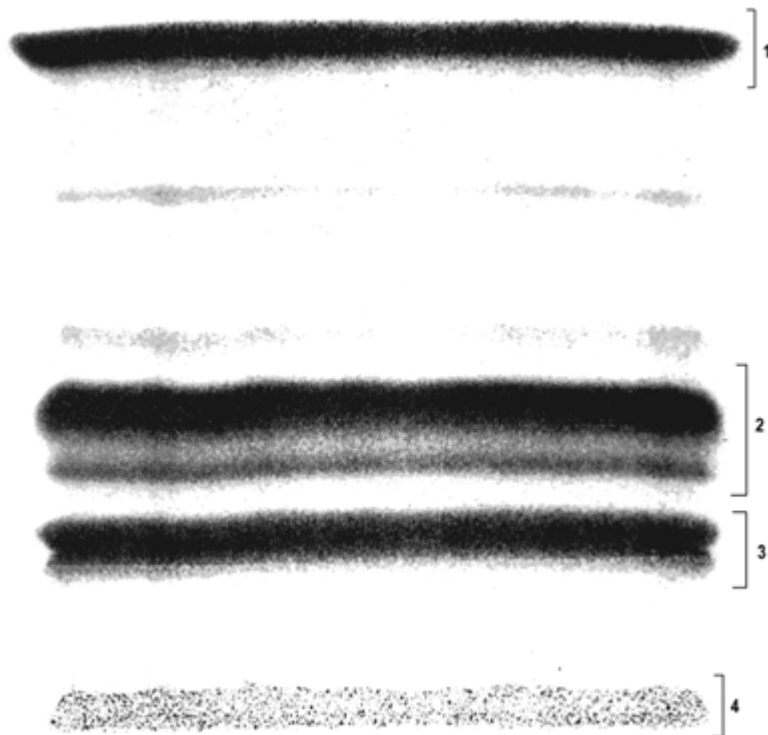

^{F12}ANNEX VI

.....

[^{F10}ANNEX VII

DETERMINATION OF THE PERCENTAGE OF 2-GLYCERYL MONOPALMITATE

1. PURPOSE AND SCOPE

This method describes the analysis procedure for determining the percentage of palmitic acid in position 2 of the triglycerides by evaluating 2-glycerol monopalmitate.

This method can be applied to liquid vegetable oils at ambient temperature (20 °C).

2. PRINCIPLE

After preparation the oil sample is subjected to the action of pancreatic lipase: partial and specific hydrolysis in positions 1 and 3 of the triglyceride molecule causes monoglycerides to appear in position 2. The percentage of 2-glycerol monopalmitate in the monoglyceride fraction is determined after silylation by capillary-column gas chromatography.

3. APPARATUS AND MATERIALS

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 3.1. 25 ml Erlenmeyer flask
- 3.2. 100, 250 and 300 ml beakers
- 3.3. Glass chromatograph column, internal diameter 21-23 mm, length 400 mm, fitted with a sintered glass disc and a stopcock
- 3.4. 10, 50, 100 and 200 ml measuring cylinders
- 3.5. 100 and 250 ml flasks
- 3.6. Rotary evaporator
- 3.7. 10 ml conical-bottomed centrifuge tubes with groundglass stopper
- 3.8. Centrifuge for 10 and 100 ml tubes
- 3.9. Thermostat permitting a stable temperature of $40 \pm 0,5$ °C
- 3.10. 1 and 2 ml graduated pipettes
- 3.11. 1 ml hypodermic syringe
- 3.12. 100 µl microsyringe
- 3.13. 1 000 ml funnel
- 3.14. Capillary gas chromatograph with an on-column cold injector for direct injection of the sample into the column and a furnace able to maintain the selected temperature to approximately 1 °C
- 3.15. On-column cold injector for direct injection of the sample into the column
- 3.16. Flame ionisation detector and electrometer
- 3.17. Recorder-integrator adapted to the electrometer with a response rate no greater than 1 sec and a variable paper roll rate
- 3.18. Capillary column made of glass or fused silica 8-12 metres long, 0,25-0,32 mm internal diameter, covered with methylpolysiloxane or phenyl methylpolysiloxane 5 %, 0,10-0,30 µm thick, useable at 370 °C
- 3.19. 10 µl microsyringe fitted with a hardened needle, at least 7,5 cm long for direct on-column injection.
4. REAGENTS
 - 4.1. Silica gel with a grain size of between 0,063 and 0,200 mm (70/280 mesh) prepared as follows: Place the silica gel in a porcelain capsule, dry in an incubator at 160 °C for four hours, then leave to cool at room temperature in a desiccator. Add water equivalent to 5 % of the mass of the silica gel as follows: Weigh 152 g silica gel into an Erlenmeyer flask then add 8 g of distilled water, stopper and shake gently to distribute the water evenly. Leave to stand for at least 12 hours before use.
 - 4.2. n-hexane (for chromatography)
 - 4.3. Isopropanol
 - 4.4. Isopropanol, 1/1 (v/v) aqueous solution

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 4.5. Pancreatic lipase. It must have an activity of between 2,0 and 10 lipase units per mg. (Pancreatic lipases with an activity of between 2 and 10 units per mg enzyme are commercially available.)
- 4.6. Buffer solution of trishydroxymethylaminomethane: 1 M aqueous solution adjusted to pH 8 (potentiometric control) by conc. HCl (1/1 v/v)
- 4.7. Enzyme-quality sodium cholate, 0,1 % aqueous solution (this solution must be used within two weeks of its preparation)
- 4.8. Calcium chloride, 22 % aqueous solution
- 4.9. Diethyl ether for chromatography
- 4.10. Developer solvent: mixture of n-hexane/diethyl ether (87:13 v:v)
- 4.11. Sodium hydroxide, 12 % by weight solution
- 4.12. Phenolphthalein, 1 % solution in ethanol
- 4.13. Carrier gas: hydrogen or helium, for gas chromatography
- 4.14. Auxiliary gases: hydrogen, 99 % minimum purity, free from moisture and organic substances, and air, for gas chromatography, of the same purity
- 4.15. Silanisation reagent: mixture of pyridine/hexamethyldisilazane, trimethylchlorosilane 9/3/1 (v/v/v). (Ready-to-use solutions are commercially available. Other silylation reagents may be used, particularly bis-trimethylsilyl trifluoroacetamide + 1 % trimethylchlorosilane, diluted with an identical volume of anhydrous pyridine.)
- 4.16. Reference samples: pure monoglycerides or monoglyceride mixtures with a known percentage composition similar to that of the sample.
5. METHOD
 - 5.1. Sample preparation
 - 5.1.1. Oils with a free acidity of less than 3 % do not need to be neutralised before chromatography on a silica gel column. Oils with a free acidity of more than 3 % must be neutralised as per point 5.1.1.1.
 - 5.1.1.1. Pour 50 g of oil and 200 ml n-hexane into the 1 000 ml funnel (3.13). Add 100 ml of isopropanol and a quantity of 12 % sodium hydroxide solution (4.11) equivalent to the free acidity of the oil plus 5 %. Shake vigorously for one minute. Add 100 ml of distilled water, shake again and leave to stand.

After decanting, remove the lower layer containing the soaps. Remove any intermediate layers (mucilage and insoluble substances). Wash the hexane solution of the neutralised oil with successive portions of 50-60 ml of the 1/1 (v/v) isopropanol/water solution (4.4) until the pink colouration of the phenolphthalein disappears.

Remove most of the hexane by vacuum distillation (use a rotary evaporator, for example) and transfer the oil into a 100 ml flask (3.5). Dry the oil in vacuum until the solvent is completely removed.

After that procedure is completed, the acidity of the oil should be less than 0,5 %.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 5.1.2. Put 1,0 g of the oil prepared as above into a 25 ml Erlenmeyer flask (3.1) and dissolve in 10 ml of developer mixture (4.10). Leave the solution to stand for at least 15 minutes before silica gel column chromatography.

If the solution is cloudy centrifuge it to ensure optimum conditions for chromatography. (Ready-to-use 500 mg silica gel SPE cartridges can be used).

- 5.1.3. Preparation of the chromatography column

Pour about 30 ml of the developer solvent (4.10) into the column (3.3), insert a piece of cotton into the bottom part of the column using a glass rod; press to eliminate the air.

In a beaker prepare a suspension of 25 g of silica gel (4.1) in about 80 ml of developer solvent and pour it into the column using a funnel.

Check that all the silica gel is in the column; wash with developer solvent (4.10), open the stopcock and allow the liquid to reach a level about 2 mm above the level of the silica gel.

- 5.1.4. Column chromatography

Weigh accurately 1,0 g of sample prepared as in point 5.1 into a 25 ml Erlenmeyer flask (3.1).

Dissolve the sample in 10 ml of developer solvent (4.10). Pour the solution into the chromatography column prepared as in point 5.1.3. Avoid disturbing the surface of the column.

Open the stopcock and pour the sample solution until it reaches the level of the silica gel. Develop with 150 ml of the developer solvent. Adjust the flow rate to 2 ml/min (so that 150 ml enters the column in about 60-70 minutes).

Recover the eluate in a previously weighed 250 ml flask. Evaporate the solvent under vacuum and remove the final traces of the solvent under a nitrogen current.

Weigh the flask and calculate the recovered extract.

(If ready-to-use silica gel SPE cartridges are used use the following method: Put 1 ml of solution (5.1.2) into the prepared cartridges with 3 ml of n-hexane.

After percolating the solution develop with 4 ml of n-hexane/diethyl ether 9/1 (v/v).

Recover the eluate in a 10 ml tube and evaporate to dry in a nitrogen current.

Expose the dry residue to pancreatic lipase (5.2). (It is essential to check the fatty acid composition before and after crossing the SPE cartridge.)

- 5.2. Hydrolysis by pancreatic lipase

5.2.1. Weigh into the centrifuge tube 0.1 g of the oil prepared as in point 5.1. Add 2 ml of buffer solution (4.6), 0,5 ml of the sodium cholate solution (4.7) and 0,2 ml of the calcium chloride solution, stirring well after each addition. Close the tube with the groundglass stopper and place in the thermostat at $40 \pm 0,5$ °C.

5.2.2. Add 20 mg of lipase, shake carefully (avoid wetting the stopper) and place the tube in the thermostat for exactly two minutes. Then remove it, shake vigorously for exactly 1 minute and leave to cool.

5.2.3. Add 1 ml of diethyl ether, stopper and shake vigorously, then centrifuge and transfer the ether solution into a clean, dry tube using a microsyringe.

- 5.3. Preparation of the silanised derivatives and gas chromatography

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 5.3.1. With a microsyringe insert 100 µl of solution (5.2.3) into a 10 ml conical-bottomed tube.
- 5.3.2. Remove the solvent under a slight nitrogen current, add 200 µl of silanisation reagent (4.15), stopper the tube and leave to stand for 20 minutes.
- 5.3.3. After 20 minutes, add 1 to 5 ml of n-hexane (depending on the chromatography conditions): the resulting solution is ready for gas chromatography.

5.4. Gas chromatography

Operating conditions:

- Injector temperature (on-column injector) lower than solvent boiling point (68 °C);
- Detector temperature: 350 °C;
- Column temperature: programming of furnace temperature: 60 °C for 1 minute, increasing by 15 °C per minute up to 180 °C, then by 5 °C per minute up to 340 °C, then 340 °C for 13 minutes;
- Carrier gas: hydrogen or helium, set at a linear velocity sufficient to obtain the resolution reflected in Figure 1. The retention time of the C₅₄ triglyceride must be 40 ± 5 minutes (see Figure 2). (The operating conditions indicated above are indicative. Operators will have to optimise them to obtain the desired resolution. The peak corresponding to 2-glyceryl monopalmitate must have a minimum height equal to 10 % of the recorder scale.)
- Quantity of substance injected: 0,5-1 µl of the n-hexane solution (5 ml) (5.3.3).

5.4.1. Identification of the peaks

The individual monoglycerides are identified from their retention times and by comparison with those obtained for standard monoglyceride mixtures under the same conditions.

5.4.2. Quantitative evaluation

The area of each peak is calculated using an electronic integrator.

6. EXPRESSION OF RESULTS

The percentage of glyceryl monopalmitate is calculated from the ratio between the area of the corresponding peak and the areas of the peaks of all the monoglycerides (see Figure 2) using the formula:

glyceryl monopalmitate (%):

$$\frac{A_x}{\Sigma A} \times 100$$

where:

A_x = area of the peak corresponding to glyceryl monopalmitate
ΣA = sum of the areas of all the monoglyceride peaks

The result must be to one decimal place.

7. ANALYSIS REPORT

The analysis report must specify:

- reference to this method,
- all the information needed for a full identification of the sample,
- the analysis result,

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- any deviation from the method, whether as the result of a decision by the parties concerned or for another reason,
- details to identify the laboratory, the date of the analysis and the signatures of those responsible for the analysis.

Figure 1

Chromatogram of the products of the silanisation reaction obtained by the action of lipase on a refined olive oil with 20 % esterified oil added (100 %)

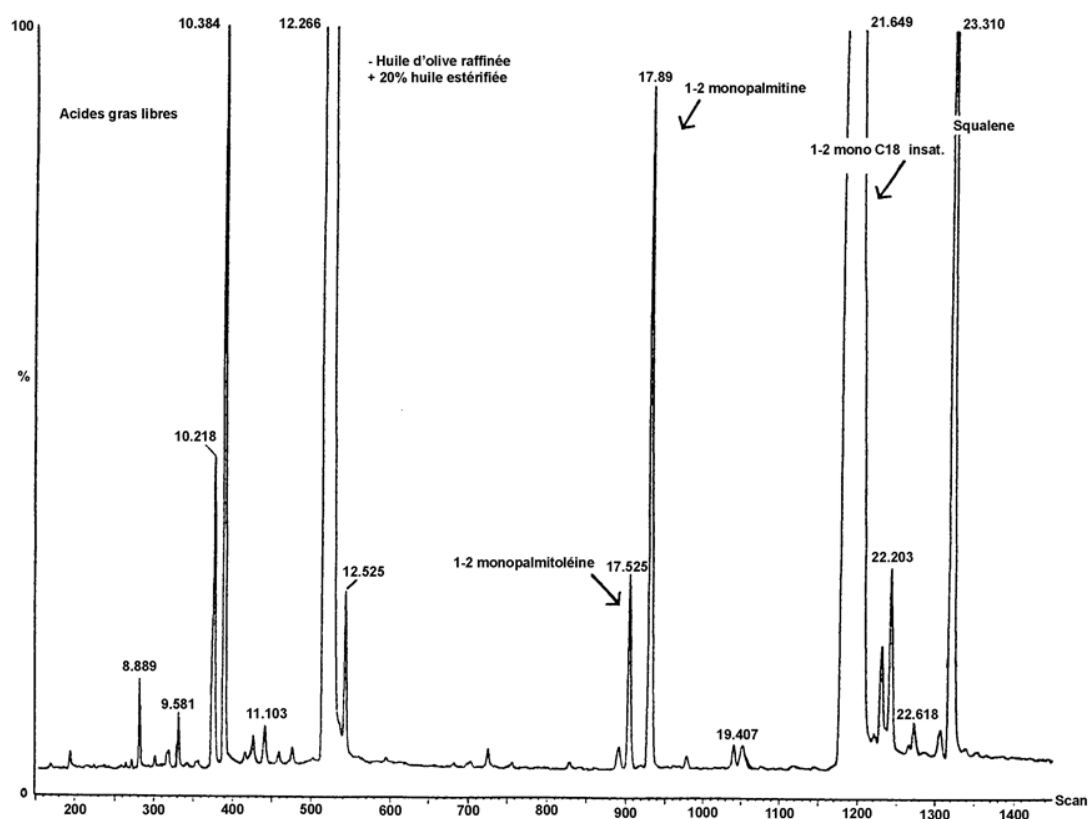


Figure 2


Chromatogram of:

- (A) unesterified olive oil, after lipase; after silanisation; under these conditions (8-12 m capillary column) the wax fraction is eluted at the same time as the diglyceride fraction or slightly afterwards.

After lipase, the triglyceride content should not exceed 15 %

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

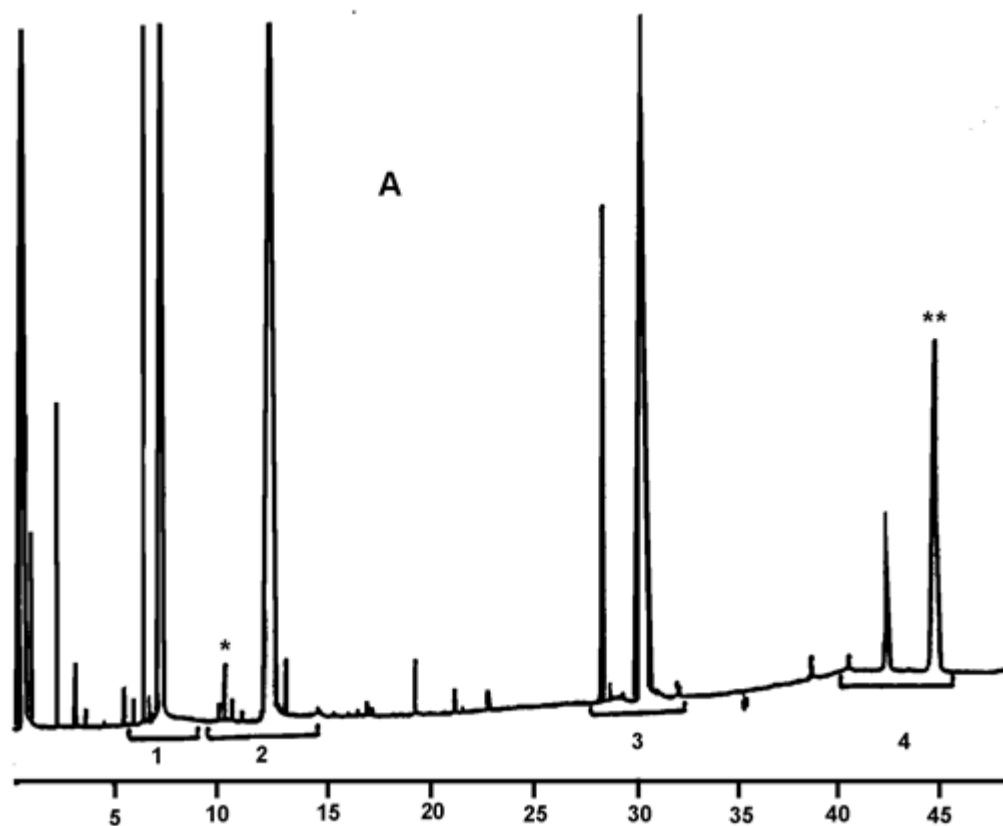


Chromatogram of:

- (B) unesterified oil after lipase; after silanisation; under these conditions (8-12 m capillary column) the wax fraction is eluted at the same time as the diglyceride fraction or slightly afterwards.

After lipase, the triglyceride content should not exceed 15 %.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

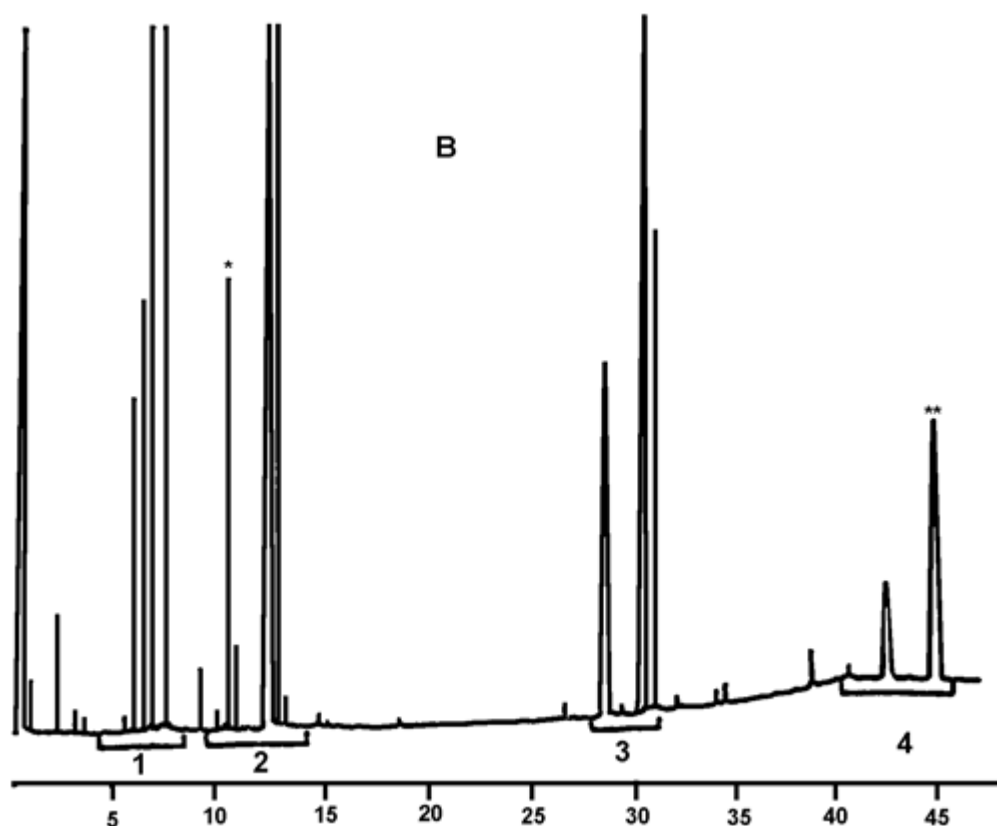

8. NOTES

Note 1. PREPARATION OF THE LIPASE

Lipases with satisfactory activity are commercially available. They can also be prepared in the laboratory in the following manner:

Cool to 0 °C 5 kg of fresh pig's pancreas. Remove the surrounding solid fat and the connective tissue and grind to a liquid paste in a blender. Stir the paste with 2,5 litres of anhydrous acetone for 4-6 hours, then centrifuge. Extract the residue three more times with the same volume of anhydrous acetone, then twice with an acetone/diethyl ether mixture (1/1 v/v) and twice with diethyl ether.

Vacuum-dry the residue for 48 hours to obtain a stable powder which can be stored for a long time in a refrigerator away from moisture.

Note 2. MONITORING LIPASE ACTIVITY

Prepare an olive oil emulsion as follows:

In a mixer stir for 10 minutes a mixture of 165 ml of a 100 g/l gum arabic solution, 15 g of crushed ice and 20 ml of a previously neutralised olive oil.

Pour 10 ml of the emulsion into a 50 ml beaker, then 0,3 ml of a 0,2 g/ml sodium cholate solution and then 20 ml of distilled water.

Put the beaker in a thermostat set at 37 °C; introduce the electrodes of the pH meter and the screw agitator.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Using a burette, add a 0,1 N sodium hydroxide solution drop by drop until a pH of 8,3 is obtained.

Add an aliquot of the lipase powder suspension in water (0,1 g/ml of lipase). As soon as the pH meter reads 8,3, start the chronometer and add the sodium hydroxide solution drop by drop at a rate which maintains the pH at 8,3. Note every minute the volume of solution consumed.

Record the data on an x/y graph with the time on the x-axis and millilitres of 0,1 N alkaline solution consumed to keep a constant pH on the y-axis. A linear graph should be obtained.

Lipase activity, expressed in lipase units per mg, is given by the following formula:

$$A = \frac{V \times N \times 100}{m}$$

where:

A	is activity in lipase units/mg
V	is the number of millilitres of 0,1 N sodium hydroxide solution per minute (calculated on the basis of the graph)
N	is the titre of the sodium hydroxide solution
m	is the mass in mg of the test lipase.

A lipase unit is defined as the quantity of enzyme which releases 10 micro-equivalents of acid per minute.]

^{F13}ANNEX VIII

[^{F13}DETERMINATION OF TRILINOLEIN CONTENT]

1. SCOPE

.....

2. FIELD OF APPLICATION

.....

3. PRINCIPLE

.....

4. APPARATUS

4.1.

4.2.

4.3.

4.4.

4.5.

5. REAGENTS

^{F13}

5.1.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

5.2.

5.3.

5.4.

5.5.

5.6.

6. PREPARATION OF SAMPLES

.....

7. PROCEDURE

7.1.

8. CALCULATION AND EXPRESSION OF RESULTS

F13

F13

F13

Note 1.

.....
Note 2. *Examples:*

.....
Note 3.

.....
Note 4.

.....
Note 5:

.....
F13

.....
F13

.....
F13]

.....

[^{F3} ANNEX IX

**SPECTROPHOTOMETRIC INVESTIGATION IN THE ULTRAVIOLET
FOREWORD**

Spectrophotometric examination in the ultraviolet can provide information on the quality of a fat, its state of preservation and changes brought about in it by technological processes.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The absorption at the wavelengths specified in the method is due to the presence of conjugated diene and triene systems. These absorptions are expressed as specific extinctions $E^{1\%}_{1\text{ cm}}$ (the extinction of 1 % solution of the fat in the specified solvent, in a thickness of 1 cm) conventionally indicated by K (also referred to as 'extinction coefficient').

1. SCOPE

The method describes the procedure for performing a spectrophotometric examination of olive oil (as described in the Appendix) in the ultraviolet.

2. PRINCIPLE OF THE METHOD

The fat in question is dissolved in the required solvent and the extinction of the solution is then determined at the specified wavelengths with reference to pure solvent. Specific extinctions are calculated from the spectrophotometer readings. The specific absorbance at 232 nm and 268 nm in iso-octane or 232 nm and 270 nm in cyclohexane for a concentration of 1 g per 100 ml in a 10 mm cell is calculated.

3. EQUIPMENT

3.1. A spectrophotometer for measuring extinction in the ultraviolet between 220 and 360 nm, with the possibility of reading individual nanometric units. Before use it is recommended that the wavelength and absorbance scales of the spectrometer be checked as follows.

3.1.1. *Wavelength scale:* This may be checked using a reference material consisting of an optical glass filter containing holmium oxide which has distinct absorption bands. The reference material is designed for the verification and calibration of the wavelength scales of visible and ultraviolet spectrophotometers having nominal spectral bandwidths of 5 nm or less. The holmium glass filter is measured in the absorbance mode against an air blank, over the wavelength range of 640 to 240 nm. For each spectral bandwidth (0,10 – 0,25 – 0,50 – 1,00 – 1,50 – 2,00 and 3,00), a baseline correction is performed with an empty cell holder. The wavelengths of the spectral bandwidth are listed in the certificate of the reference material in ISO 3656.

3.1.2. *Absorbance scale:* This may be checked using a reference material consisting of 4 solutions of potassium dichromate in perchloric acid sealed in four UV quartz cells to measure the linearity and photometric accuracy reference in the UV. The potassium dichromate filled cells (40 mg/ml, 60 mg/ml, 80 mg/ml and 100 mg/ml) are measured against a perchloric acid blank. The net absorbance values are listed in the certificate of the reference material in ISO 3656.

3.2. Rectangular quartz cells, with covers, having an optical length of 1 cm. When filled with water or other suitable solvent the cells should not show differences between them of more than 0,01 extinction units.

3.3. 25 ml graduated flasks.

3.4. Analytical balance, capable of being read to the nearest 0,0001 g.

4. REAGENTS

Use only reagents of recognized analytical grade, unless otherwise stated.

Solvent: Iso-octane (2,2,4-trimethylpentane) for the measurement at 232 nm and 268 nm or cyclohexane for the measurement at 232 nm and 270 nm, having an absorbance less than 0,12 at 232 nm and less than 0,05 at 250 nm against distilled water, measured in a 10 mm cell.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

5. PROCEDURE

- 5.1. The sample in question must be perfectly homogeneous and without suspected impurities. Oils which are liquid at ambient temperature are to be filtered through paper at a temperature of approximately 30 °C, hard fats are to be homogenized and filtered at a temperature of not more than 10 °C above the melting point.
- 5.2. Weigh accurately approximately 0,25 g (to the nearest 1 mg) of the sample so prepared into a 25 ml graduated flask, make up to the mark with the solvent specified and homogenize. The resulting solution must be perfectly clear. If opalescence or turbidity is present filter quickly through paper.
- 5.3. Fill a quartz cell with the solution obtained and measure the extinctions at an appropriate wavelength between 232 and 276 nm, using the solvent used as a reference.

The extinction values recorded must lie within the range 0,1 to 0,8. If not the measurements must be repeated using more concentrated or more dilute solutions as appropriate.

NOTE: It may not be necessary to measure the absorbance over the full wavelength range.

6. EXPRESSION OF THE RESULTS

- 6.1. Record the specific extinctions (extinction coefficients) at the various wavelengths calculated as follows:

$$K_{\lambda} = \left(\frac{E_{\lambda}}{c \cdot s} \right)$$

where:

- | | | |
|---------------|---|---|
| K_{λ} | = | specific extinction at wavelength λ , |
| E_{λ} | = | extinction measured at wavelength λ ; |
| c | = | concentration of the solution in g/100 ml; |
| s | = | thickness of the quartz cells in cm. |

The results are to be expressed to two decimal places.

6.2. Variation of the specific extinction (ΔK)

Spectrophotometric analysis of olive oil in accordance with the official method in the Union legislation involves also the determination of the variation of the absolute value of the specific extinction (ΔK), which is given by:

$$\Delta K = \left| K_m - \left(\frac{K_{m-4} + K_{m+4}}{2} \right) \right|$$

where K_m is the specific extinction at wavelength m , the wavelength for maximum absorption depends on the solvent used: 270 for cyclohexane and 268 for iso-octane.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix

OLIVE OIL CHARACTERISTICS

Category	Fatty acid methyl esters (FAMES) and fatty acid ethyl esters (FAEEs)	Acidity (%) ^(*)	Peroxide index mEq/kg ^(*)	Waxes mg/kg ^(**)	Stigmaterpene glycerol monopalmitate (%)	Stigmaterpene diacylglycerols (DPLC) and ECN42 (theoretical calculation)	Stigmaterpene triacylglycerols (TPLC) and ECN42	K ₂₇₀ (*) K ₂₇₀ or K ₂₆₈ ^(a)	Delta-K ^(*) ^(c)	Organoleptic evaluation Median defect (Md) (*)	Organoleptic evaluation Fruity median (Mf) (*)
1.	Extra virgin olive oil Σ FAME + FAEE ≤ 75 mg/kg or 75 mg/kg < Σ FAME + FAEE ≤ 150 mg/kg and (FAEE/FAME) ≤ 1,5	≤ 0,8	≤ 20	≤ 250	≤ 0,9 if total palmitic acid % ≤ 14 % ≤ 1,0 if total palmitic acid % > 14 %	≤ 0,10	≤ 0,2	≤ 2,50	≤ 0,22	≤ 0,01	Md = 0 Mf > 0
2.	Virgin olive oil	—	≤ 2,0	≤ 20	≤ 0,9 if total palmitic acid % ≤ 14 %	≤ 0,10	≤ 0,2	≤ 2,60	≤ 0,25	≤ 0,01	Md ≤ 3,5 Mf > 0

a Total isomers which could (or could not) be separated by capillary column.

b Or where the median defect is less than or equal to 3,5 and the fruity median is equal to 0.

c Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.

d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

e K 270 if solvent is cyclohexane, K 268 if solvent is iso-octane.]

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

						$\leq 1,0$ if total palmitic acid % > 14 %							
3.	Lampante olive oil	> 2,0	—	$\leq 300^c$	$\leq 0,9$ if total palmitic acid % \leq 14 %	$\leq 0,50$	$\leq 0,3$	—	—	—	Md > 3,5 ^b	—	
					$\leq 1,1$ if total palmitic acid % > 14 %								
4.	Refined olive oil	$\leq 0,3$	≤ 5	≤ 350	$\leq 0,9$ if total palmitic acid % \leq 14 %	—	$\leq 0,3$	—	$\leq 1,10$	$\leq 0,16$	—	—	
					$\leq 1,1$ if total palmitic acid % > 14 %								
5.	Olive oil composed of refined and virgin	—	$\leq 1,0$	≤ 15	≤ 350	$\leq 0,9$ if total palmitic acid % \leq 14 %	—	$\leq 0,3$	—	$\leq 0,90$	$\leq 0,15$	—	—

a Total isomers which could (or could not) be separated by capillary column.

b Or where the median defect is less than or equal to 3,5 and the fruity median is equal to 0.

c Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.

d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

e K 270 if solvent is cyclohexane, K 268 if solvent is iso-octane.]

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

	olive oils				$\leq 1,0$ if total palmitic acid % > 14 %								
6.	Crude olive-pomace oil	—	—	—	> 350 ^d	$\leq 1,4$	—	$\leq 0,6$	—	—	—	—	—
7.	Refined olive-pomace oil	—	$\leq 0,3$	≤ 5	> 350	$\leq 1,4$	—	$\leq 0,5$	—	$\leq 2,00$	$\leq 0,20$	—	—
8.	Olive-pomace oil	—	$\leq 1,0$	≤ 15	> 350	$\leq 1,2$	—	$\leq 0,5$	—	$\leq 1,70$	$\leq 0,18$	—	—

a Total isomers which could (or could not) be separated by capillary column.

b Or where the median defect is less than or equal to 3,5 and the fruity median is equal to 0.

c Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be lampante olive oil if the total aliphatic alcohol content is less than or equal to 350 mg/kg or if the erythrodiol and uvaol content is less than or equal to 3,5 %.

d Oils with a wax content of between 300 mg/kg and 350 mg/kg are considered to be crude olive-pomace oil if the total aliphatic alcohol content is above 350 mg/kg and if the erythrodiol and uvaol content is greater than 3,5 %.

e K 270 if solvent is cyclohexane, K 268 if solvent is iso-octane.]

ANNEX X A

ANALYSIS BY GAS CHROMATOGRAPHY OF METHYL ESTERS OF FATTY ACIDS

1. SCOPE

This method gives general guidance for the application of gas chromatography, using packed or capillary columns, to determine the qualitative and quantitative composition of a mixture of fatty acid methyl esters obtained in accordance with the method specified in Annex X B.

The method is not applicable to polymerized fatty acids.

2. REAGENTS

2.1. Carrier gas

Inert gas (nitrogen, helium, argon, hydrogen, etc.), thoroughly dried and with an oxygen content of less than 10 mg/kg.

Note 1.

Hydrogen, which is used as a carrier gas only with capillary columns, can double the speed of analysis but is hazardous. Safety devices are available.

2.2. Auxiliary gases

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

2.2.1. Hydrogen (purity $\geq 99,9\%$), free from organic impurities.

2.2.2. Air or oxygen, free from organic impurities.

2.3. Reference standard

A mixture of methyl esters of pure fatty acids, or the methyl esters of a fat of known composition, preferably similar to that of the fatty matter to be analyzed.

Care shall be taken to prevent the oxidation of polyunsaturated fatty acids.

3. APPARATUS

The instructions given relate to the usual equipment used for gas chromatography, employing packed and/or capillary columns and a flame-ionization detector. Any apparatus giving the efficiency and resolution specified in 4.1.2 is suitable.

3.1. Gas chromatograph

The gas chromatograph shall comprise the following elements.

3.1.1. Injection system

Use an injection system either:

- (a) with packed columns, having the least deadspace possible (in this case the injection system shall be capable of being heated to a temperature 20 to 50 °C higher than that of the column); or
- (b) with capillary columns, in which case the injection system shall be specially designed for use with such columns. It may be of the split type or it may be of the splitless on column injector type.

Note 2.

In the absence of fatty acids with less than 16 carbon atoms, a moving needle injector may be used.

3.1.2. Oven

The oven shall be capable of heating the column to a temperature of at least 260 °C and of maintaining the desired temperature to within 1 °C with a packed column and within 0,1 °C with a capillary column. The last requirement is particularly important when a fused silica tube is used.

The use of temperature-programmed heating is recommended in all cases, and in particular for fatty acids with less than 16 carbon atoms.

3.1.3. Packed column

3.1.3.1. Column, constructed of a material inert to the substances to be analyzed (i.e. glass or stainless steel) having the following dimensions:

- (a) length: 1 to 3 m. A relatively short column should be used when long-chain fatty acids (above C₂₀) are present. When analyzing acids with 4 or 6 carbon atoms, it is recommended that a column 2 m in length is used;
- (b) internal diameter: 2 to 4 mm.

Note 3.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

If polyunsaturated components with more than three double bonds are present, they may be decomposed in a stainless steel column.

Note 4.

A system with packed twin columns may be used.

3.1.3.2. Packing, comprising the following elements:

- (a) *support*: acid-washed and silanized diatomaceous earth, or other suitable inert support with a narrow range of grain size (25 µm range between the limits 125 to 200 µm) the average grain size being related to the internal diameter and length of the column;
- (b) *stationary phase*: polyester type of polar liquid (e.g. diethylene glycol polysuccinate, butanediol polysuccinate, ethyleneglycol polyadipate, etc.) cyanosilicones or any other liquid permitting the chromatographic separation required (see clause 4). The stationary phase should amount to 5 to 20 % (m/m) of the packing. A non-polar stationary phase can be used for certain separations.

3.1.3.3. Conditioning of the column

With the column disconnected, if possible, from the detector, gradually heat the oven to 185 °C and pass a current of inert gas through the freshly prepared column at a rate of 20 to 60 ml/min for at least 16 hours at this temperature, and for a further 2 hours at 195 °C.

3.1.4. Capillary column

3.1.4.1. Tube, made of a material inert to the substances to be analysed (usually glass or fused silica). The internal diameter shall be between 0,2 and 0,8 mm. The internal surface shall undergo an appropriate treatment (e.g. surface preparation, inactivation) before receiving the stationary phase coating. A length of 25 mm is sufficient in most cases.

3.1.4.2. Stationary phase, usually of the type polyglycol (poly(ethylene glycol) 20 000), polyester (butanediol polysuccinate) or polar polysiloxane (cyanosilicones). Bonded (cross-linked) columns are suitable.

Note 5.

There is a risk of polar polysiloxanes giving rise to difficulties in the identification and separation of linolenic acid and C₂₀ acids.

The coatings shall be thin, i.e. 0,1 to 0,2 µm.

3.1.4.3. Assembly and conditioning of the column

Observe the normal precautions for assembling capillary columns, i.e. arrangement of the column in the oven (support), choice and assembly of joints (leak tightness), positioning of the ends of the column in the injector and the detector (reduction of dead-spaces). Place the column under a flow of carrier gas (e.g. 0,3 bar (30 kPa) for a column of length 25 mm and internal diameter 0,3 mm).

Condition the column by temperature programming of the oven at 3 °C/min from ambient temperature to a temperature 10 °C below the decomposure limit of the stationary phase. Maintain the oven at this temperature for one hour until stabilization of the baseline. Return it to 180 °C to work under isothermal conditions.

Note 6.

Suitably pre-conditioned columns are available commercially.

3.1.5. Detector, preferably capable of being heated to a temperature above that of the column.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

3.2. Syringe

The syringe shall have a maximum capacity of 10 µl, and be graduated in 0,1 µl divisions.

3.3. Recorder

If the recorder curve is to be used to calculate the composition of the mixture analysed, an electronic recorder of high precision, compatible with the apparatus used, is required. The recorder shall have the following characteristics:

- (a) rate of response, below 1,5 s, preferably 1 s (the rate of response is the time taken for the recording pen to pass from 0 to 90 % following the sudden introduction of a 100 % signal);
- (b) width of the paper, 20 cm minimum;
- (c) paper speed, adjustable to values between 0,4 and 2,5 cm/min.

3.4. Integrator

Rapid and accurate calculation can be performed with the help of an electronic integrator. This shall give a linear response with adequate sensitivity, and the correction for deviation of the base-line shall be satisfactory.

4. PROCEDURE

The operations described in 4.1 to 4.3 relate to the use of a flame-ionization detector.

As an alternative a gas chromatograph employing a catharometer detector (working on the principle of thermal conductivity changes) may be used. The operating conditions are then modified as described in clause 6.

4.1. Test conditions

4.1.1. Selection of optimum operating conditions

4.1.1.1. Packed column

In the selection of the test conditions, the following variables should be taken into account:

- (a) the length and diameter of the column;
- (b) the nature and amount of the stationary phase;
- (c) the temperature of the column;
- (d) the carrier gas flow;
- (e) the resolution required;
- (f) the size of the test portion, selected in such a way that the assembly of the detector and electrometer gives a linear response;
- (g) the duration of analysis.

In general, the values given in Table 1 and Table 2 will lead to the desired results, i.e. at least 2 000 theoretical plates per metre of column length for methyl stearate and its elution within about 15 minutes.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Where the apparatus allows it, the injector should be at a temperature of about 200 °C and the detector at a temperature equal to or higher than that of the column.

As a rule, the ratio of the flow-rate of the hydrogen supplied to the flame-ionization detector to that of the carrier gas varies from 1:2 to 1:1 depending on the diameter of the column. The flow of oxygen is about 5 to 10 times that of the hydrogen.

Table 1

Internal diameter of columnmm	Carrier gas flowml/min
2	15 to 25
3	20 to 40
4	40 to 60

Table 2

Concentration of stationary phase% (m/m)	Column temperature°C
5	175
10	180
15	185
20	185

4.1.1.2. Capillary column

The properties of efficiency and permeability of capillary columns mean that the separation between constituents and the duration of the analysis are largely dependent on the flow-rate of the carrier-gas in the column. It will therefore be necessary to optimize the operating conditions by acting on this parameter (or more simply on the headloss of the column), according to whether one wishes to improve the separations or to make a rapid analysis.

4.1.2. Determination of the number of theoretical plates (efficiency) and resolution (See Figure 1)

Carry out the analysis of a mixture of methyl stearate and methyl oleate in about equivalent proportions (for example, methyl esters from cocoa butter).

Choose the temperature of the column and the carrier gas flow so that the maximum of the methyl stearate peak is recorded about 15 minutes after the solvent peak. Use a sufficient quantity of the mixture of methyl esters that the methyl stearate peak occupies about three-quarters of the full scale.

Calculate the number of theoretical plates, n (efficiency), using the formula:

$$n = 16 \left[\frac{d_{r1}}{w_1} \right]^2$$

and the resolution, R , using the formula:

$$R = \frac{2\Delta}{w_1 + w_2}$$

where:

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- d_{r1} is the retention distance, in millimetres, from the start of the chromatogram to the maximum of the peak for methyl stearate;
- ω_1 and ω_2 are the widths, in millimetres, of the peaks for methyl stearate and methyl oleate respectively, measured between the points of intersection of the tangents at the points of inflection of the curve with the base-line;
- Δ is the distance, in millimetres, between the peak maxima for methyl stearate and methyl oleate^[F16;]

Textual Amendments

- F16** Substituted by [Commission Regulation \(EEC\) No 1429/92 of 26 May 1992 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)

^[F17]and the resolution index, I_r , using the formula

Textual Amendments

- F17** Inserted by [Commission Regulation \(EEC\) No 1429/92 of 26 May 1992 amending Regulation \(EEC\) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.](#)


$\frac{a}{b}$

where:

- a = the height of the smallest peak, measured from the base line;
- b = the height of the lowest point of the valley between the two adjacent peaks, measured from the base line.]

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

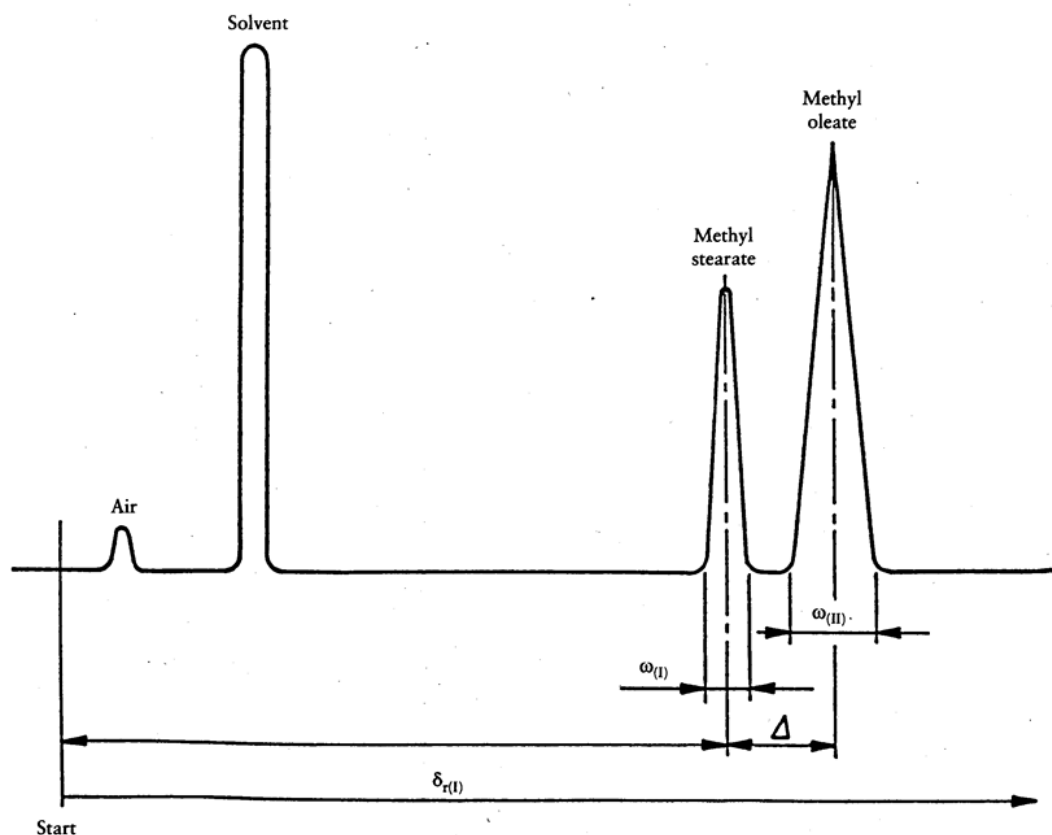

The operating conditions to be selected are those which will afford at least 2 000 theoretical plates per metre of column length for methyl stearate and a resolution of at least 1,25.

4.2. Test portion

Using the syringe (3.2) take 0,1 to 2 μ l of the solution of methyl esters prepared according to Annex X B and inject them into the column.

In the case of esters not in solution, prepare a solution of approximately 100 mg/ml in heptane of chromatographic quality, and inject 0,1 to 1 ml of this solution.

If the analysis is for constituents present only in trace amounts, the size of the test portion may be increased (up to 10-fold).

4.3. Analysis

Generally, the operating conditions shall be those defined in 4.1.1.

Nevertheless, it is possible to work with a lower column temperature when the determination of fatty acids with fewer than 12 carbon atoms is required, or at a higher temperature when determining fatty acids with more than 20 carbon atoms. On occasion, it is possible to employ temperature programming in both these cases. For example, if the sample contains the methyl esters of fatty acids with fewer than 12 carbon atoms, inject the sample at 100 °C (or at 50 to 60 °C if butyric acid is present) and immediately raise the temperature at a rate of 4 to 8 °C/min to the optimum. In certain cases, the two procedures can be combined.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

After the programmed heating, continue the elution at a constant temperature until all the components have been eluted. If the instrument does not have programmed heating, use it at two fixed temperatures between 100 and 195 °C.

If necessary, it is recommended that an analysis be carried out on two fixed phases with different polarities to verify the absence of masked peaks, for example in the case of the simultaneous presence of C_{18:3} and C_{20:0}, or C_{18:3} and C_{18:2} conjugated.

4.4. Preparation of the reference chromatogram and reference graphs

Analyze the reference standard mixture (2.3) using the same operating conditions as those employed for the sample, and measure the retention times or retention distances for the constituent fatty acids. Construct on semi-logarithmic paper, for any degree of unsaturation, the graphs showing the logarithm of retention time or distance as a function of the number of carbon atoms. In isothermal conditions, the graphs for straight-chain acids of the same degree of unsaturation should be straight lines. These lines should be approximately parallel.

It is necessary to avoid conditions such that 'masked peaks' exist, i.e. where the resolution is insufficient to separate two constituents.

5. EXPRESSION OF RESULTS

5.1. Qualitative analysis

Identify the methyl ester peaks for the sample from the graphs prepared in 4.4, if necessary by interpolation.

5.2. Quantitative analysis

5.2.1. Determination of the composition

Apart from exceptional cases, use the internal normalization method, i.e. assume that the whole of the components of the sample are represented on the chromatogram, so that the total of the areas under the peaks represents 100 % of the constituents (total elution).

If the equipment includes an integrator, use the figures obtained therefrom. If not, determine the area under each peak by multiplying the height of the peak by its width at mid-height, and where necessary take into account the various attenuations used during the recording.

5.2.2. Method of calculation

5.2.2.1. General case

Calculate the content of a given component *i*, expressed as a percentage by mass of methyl esters, by determining the percentage represented by the area of the corresponding peak relative to the sum of the areas of all the peaks, using the following formula:

$$\frac{A_i}{A} \times 100$$

where:

A_{*i*} is the area under the peak corresponding to component *i*;

A is the sum of the areas under all the peaks.

Give the result to one decimal place.

Note 7:

In this general case, the result of the calculation based on relative areas is considered to represent a percentage by mass. For the cases in which this assumption is not allowed, see 5.2.2.2.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

5.2.2.2. Use of correction factors

In certain cases, for example in the presence of fatty acids with fewer than eight carbon atoms or of acids with secondary groups, when using thermal conductivity detectors or where the highest degree of accuracy is particularly required, correction factors should be used to convert the percentages of peak areas into mass percentages of the components.

Determine the correction factors with the help of a chromatogram derived from the analysis of a reference mixture of methyl esters of known composition, carried out under operating conditions identical with those used for the sample.

For this reference mixture, the percentage by mass of component i is given by the formula:

$$\frac{m_i}{m} \times 100$$

where:

m_i is the mass of component i in the reference mixture;

m is the total of the masses of the various components of the reference mixture.

From the chromatogram of the reference mixture (4.4) calculate the percentage (area/area) for component i as follows:

$$\frac{A_i}{A} \times 100$$

where:

A_i is the area under the peak corresponding to component i ;

A is the sum of the areas under all the peaks.

The correction factor is then calculated as:

$$K_i = \frac{m_i \times A}{A_i \times m}$$

Commonly, the correction factors are expressed relative to K_{C16} , so that the relative factors become:

$$K'_i = \frac{K_i}{K_{C16}}$$

For the sample, the content of each component i , expressed as a percentage by mass of methyl esters, is:

$$\frac{K'_i \times A_i}{\sum (K'_i \times A_i)} \times 100$$

Give the results to one decimal place.

5.2.2.3. Use of an internal standard

In certain analyses (for example where not all of the fatty acids are quantified, such as when acids with four and six carbons are present alongside acids with 16 and 18 carbons, or when it is necessary to determine the absolute amount of a fatty acid in a sample) it is necessary to use an Internal Standard. Fatty acids with five, 15 or 17 carbons are frequently used. The correction factor (if any) for the Internal Standard should be determined.

The percentage by mass of component i , expressed as methyl esters, is then given by the formula:

$$\frac{m_s \times K'_i \times A_i}{m \times K'_i \times A_i} \times 100$$

where:

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- A_i is the area under the peak corresponding to component *i*;
- A_s is the area under the peak corresponding to the Internal Standard;
- K'_i is the correction factor for component *i* (relative to K_{C1});
- K'_s is the correction factor for the Internal Standard (relative to K_{C16});
- m is the mass, in milligrams, of the test portion;
- m_s is the mass, in milligrams, of the Internal Standard.

Give the results to one decimal place.

[^{F17}6. SPECIAL CASE — DETERMINATION OF TRANS-ISOMERS

It is possible to determine the content of trans-isomers in fatty acids with a number of carbon atoms between 10 and 24 by separating the methyl esters using gas chromatography capillary columns having a specific polarity.

- 6.1. A capillary column made of silica having an internal diameter of between 0,25 mm and 0,32 mm and a length of 50 m, coated with cyanopropylsilicon, the thickness of the coating being between 0,1 and 0,3 μm (type SP 2380, C.P. sil 88, silor 10 and similar types).
- [^{F10}6.2. The methyl esters are prepared using procedure B set out in Annex XB. Fatty substances having a free acidity over 3 % must first be neutralised in accordance with point 5.1.1 of Annex VII.]
- 6.3. The operating conditions for gas chromatography are overall as follows:
- column temperature set between 150 °C and 230 °C (for example 165 °C for 15 minutes then increasing by 5 °C a minute to 200 °C);
 - injector temperature: 250 °C if the splitting system is used or the initial temperature of the column if the on-column system is used;
 - detector temperature: 260 °C;
 - flow rate of the carrier gas (helium and hydrogen): 1,2 ml a minute.


The quantity injected must be such that in the conditions of sensitivity employed the height of the peak corresponding to the methyl ester of the arachidic acid is equal to or greater than 20 % of the bottom of the scale.

- 6.4. Identification of the various methyl esters is effected on the basis of the retention times which are compared with those for the reference mixtures (as indicated at point 2.3).

The esters of trans fatty acids are eluted before the corresponding cis-isomers. An example of a chromatogram is given in figure 2.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

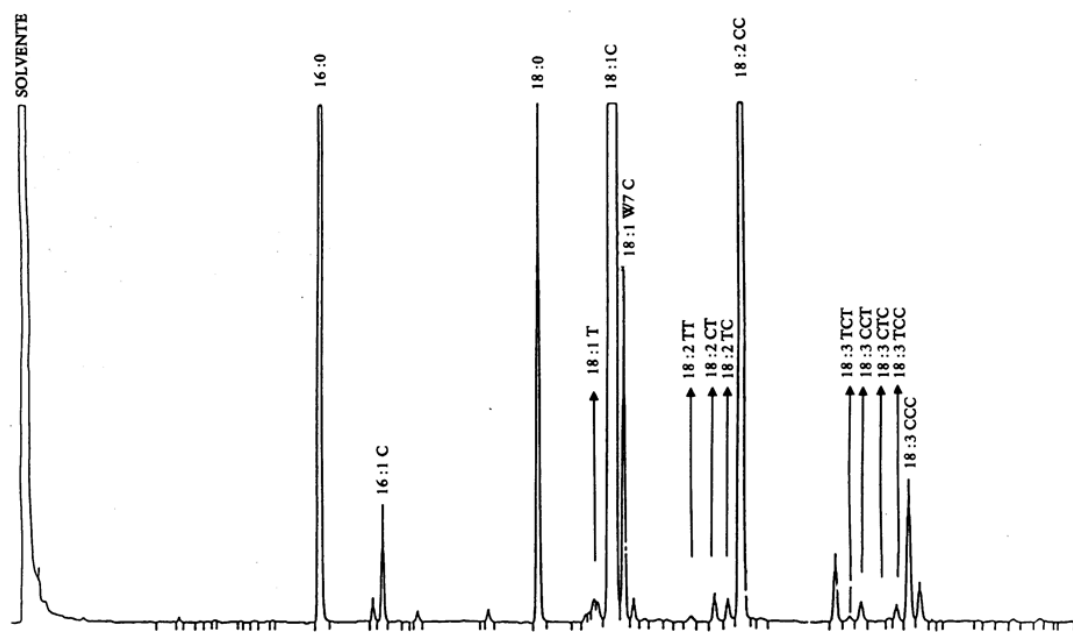

- 6.5. The efficiency of the column determined in accordance with point 4.1.2 must be such as to allow separation of certain critical couples, for example the couple formed by the massif of the transoleic acids and the oleic acid peak, trans C18:1/cis C18:1, with a resolution index greater than 2.
- 6.6. The percentage of the various trans fatty acids is calculated on the basis of the relationship between the surface of the relevant peak and the sum of the surfaces of all the peaks present.

The percentages of:

- the trans octadecenoic acids (T 18: 1) indicated in Annex I to this Regulation as the sum of the transoleic isomers;
- the cis-trans and trans-cis octadecadienoic acids [(CT/TC) 18: 2] indicated in Annex I to this Regulation as the sum of the translinoleic isomers;
- the trans-cis-trans, cis-cis-trans, cis-trans-cis, trans-cis-cis, octadecatrienoic acids [(TCT + CCT + CTC + TCC) 18: 3], indicated in Annex I of this Regulation as the sum of the translinolenic isomers

are taken into account.

Note 8:

Taking into account the particular characteristics of this method, please give the results with 2 decimals.]

[^{F167}.] SPECIAL CASE — USE OF A CATHAROMETER DETECTOR (WORKING ON THE PRINCIPLE OF THERMAL CONDUCTIVITY CHANGES)

A gas chromatograph employing a detector working on the principle of thermal conductivity changes (a catharometer) may also be used for the determination of the qualitative and quantitative composition of a mixture of fatty acid methyl esters. If it is used, the conditions specified in clause 3 and clause 4 should be modified as shown in Table 3.

For quantitative analysis, use the correction factors defined in 5.2.2.2.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Table 3

Variable	Value/condition
Column	Length: 2 to 4 m Internal diameter: 4 mm
Support	Grain size between 160 and 200 µm
Concentration of stationary phase	15 to 25 % (m/m)
Carrier gas	Helium or, failing this, hydrogen, with as low an oxygen content as possible
Auxiliary gases	None
Injector temperature	From 40 to 60 °C above that of the column
Column temperature	180 to 200 °C
Flow of carrier gas	Usually between 60 and 80 ml/min
Size of test portion injected	Usually between 0,5 and 2 µl

[^{F16}8.] TEST REPORT

The test report shall specify the methods used for the preparation of the methyl esters and for the gas chromatographic analysis, and the results obtained. It shall also mention all operating details not specified in this International Standard, or regarded as optional, together with details of any incidents which may have influenced the results.

The test report shall include all information necessary for the complete identification of the sample.

[^{F6}ANNEX X B

PREPARATION OF THE FATTY ACID METHYL ESTERS FROM OLIVE OIL AND OLIVE-POMACE OIL

The following two methods are recommended for preparing the fatty acid methyl esters from olive oils and olive-pomace oils:

- Method A : Trans-esterification with cold methanolic solution of potassium hydroxide
- Method B : Methylation by heating with sodium methylate in methanol followed by esterification in acid medium.

Each method will be applied according to the analytical parameter to be determined and the oil category as indicated below:

- (a) determination of difference between actual and theoretical content of triglycerides with ECN42 (Δ ECN42):
- method A will be applied to samples of all the oil categories after purification of the oil by passing it through a silica gel column;
- (b) determination of the fatty acid composition:

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- method A will be applied directly to samples of the following oil categories:
 - virgin olive oils with an acidity of less than 3,3 %,
 - refined olive oil,
 - olive oil (blend of virgin olive oils and refined olive oil),
 - refined olive-pomace oil,
 - olive-pomace oil (blend of virgin olive oils and refined olive-pomace oil);
 - method B will be applied directly to samples of the following oil categories:
 - virgin olive oil with an acidity of more than 3,3 %,
 - crude olive-pomace oil;
- (c) determination of trans-isomers of fatty acids:
- method A will be applied directly to samples of the following oil categories:
 - virgin olive oils with an acidity of less than 3,3 %,
 - refined olive oil,
 - olive oil (blend of virgin olive oils and refined olive oil),
 - refined olive-pomace oil,
 - olive-pomace oil (blend of virgin olive oils and refined olive-pomace oil);
 - method A will be applied to the following categories of oils after purification of the oil by passing it through a silica gel column:
 - virgin olive oil with an acidity of more than 3,3 %,
 - crude olive-pomace oil.

PURIFICATION OF OIL SAMPLES

When necessary, the samples will be purified by passing the oil through a silica gel column, eluting with hexane/diethyl ether (87:13, v/v) as described in IUPAC method 2.507.

Alternatively, solid-phase extraction on silica gel phase cartridges can be used. A silica gel cartridge (1 g, 6 ml) is placed in a vacuum elution apparatus and washed with 6 ml of hexane. The vacuum is released to prevent the column from becoming dry and then a solution of the oil (0,12 g approximately) in 0,5 ml of hexane is loaded into the column and vacuum is applied. The solution is pulled down and then eluted with 10 ml of hexane/diethyl ether (87:13 v/v) under vacuum. The combined eluates are homogenised and divided in two similar volumes. An aliquot is evaporated to dryness in a rotary evaporator under reduced pressure at room temperature. The pomace is dissolved in 1 ml of heptane and the solution is ready for fatty acid analysis by GC. The second aliquot is evaporated and the pomace is dissolved in 1 ml of acetone for triglyceride analysis by HPLC, if necessary.

METHODS FOR PREPARING THE FATTY ACID METHYL ESTERS

1. Method A: Trans-esterification with cold methanolic solution of potassium hydroxide

1.1. Purpose

This rapid method is applicable to olive oils and olive-pomace oils with a free fatty acid content of less than 3,3 %. Free fatty acids are not esterified by potassium hydroxide. Fatty acid ethyl esters are trans-esterified at a lower rate than glyceridic esters and may be only partially methylated.

1.2. Principle

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Methyl esters are formed by trans-esterification with methanolic potassium hydroxide as an intermediate stage before saponification takes place (title 5 in ISO-5509:2000, title 5 in IUPAC method 2.301).

1.3. Reagents

Methanol containing not more than 0,5 % (m/m) water.

Heptane, chromatographic quality.

Potassium hydroxide, approximately 2 N methanolic solution: dissolve 11,2 g of potassium hydroxide in 100 ml of methanol.

1.4. Apparatus

Screw-top test tubes (5 ml volume) with cap fitted with a PTFE joint.

Graduated or automatic pipettes, 2 ml and 0,2 ml

1.5. Procedure

In a 5 ml screw-top test tube weigh approximately 0,1 g of the oil sample. Add 2 ml of heptane, and shake. Add 0,2 ml of 2 N methanolic potassium hydroxide solution, put on the cap fitted with a PTFE joint, tighten the cap, and shake vigorously for 30 seconds. Leave to stratify until the upper solution becomes clear. Decant the upper layer containing the methyl esters. The heptane solution is suitable for injection into the gas chromatograph. It is advisable to keep the solution in the refrigerator until gas chromatographic analysis. Storage of the solution for more than 12 hours is not recommended.

2. Method B: Methylation by heating with sodium methylate in methanol followed by esterification in acid medium

2.1. Purpose

This method is applicable to olive oils and olive-pomace oils with a free fatty acid content of more than 3,3 %.

2.2. Principle

Neutralisation of the free fatty acids and alkaline methanolysis of the glycerides, followed by esterification of the fatty acids in acid medium (title 4.2. in IUPAC method 2.301).

2.3. Reagents

- heptane, chromatographic quality,
- methanol containing not more than 0,05 % (m/m) water,
- sodium methylate, 0,2 N methanolic solution: dissolve 5 g of sodium in 1 000 ml of methanol (this may be prepared from commercial solutions),
- phenolphthalein, 0,2 % methanolic solution,
- sulphuric acid, 1 N in methanolic solution: add 3 ml of 96 % sulphuric acid to 100 ml of methanol,
- saturated solution of sodium chloride in water.

2.4. Apparatus

- 50 ml flat-bottomed volumetric flask with long, narrow, ground neck,
- reflux condenser: air condenser (1 m long) with ground joint appropriate to the neck of the flask,
- boiling chips,

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

— glass funnel.

2.5. Procedure

Transfer about 0,25 g of the oil sample into a 50 ml ground-necked volumetric flask. With the aid of a funnel, add 10 ml of 0,2 N sodium methylate in methanol and the boiling chips. Fit a reflux condenser, shake, and bring to the boil. The solution should become clear, which usually occurs in about 10 minutes. The reaction is complete after 15 minutes. Remove the flask from the source of heat, wait until the reflux stops, remove the condenser, and add two drops of phenolphthalein solution. Add a few ml of 1 N sulphuric acid in methanol solution until the solution becomes colourless and then add 1 ml in excess. Fit the condenser and boil again for 20 minutes. Withdraw from the source of heat and cool the flask under running water. Remove the condenser, add 20 ml of saturated sodium chloride solution, and shake. Add 5 ml of heptane, plug the flask, and shake vigorously for 15 seconds.

Leave to settle until the two phases have separated. Add saturated sodium chloride solution again until the aqueous layer reaches the lower end of the flask neck. The upper layer containing the methyl esters fills the flask neck. This solution is ready to be injected in the GC.

Caution: Methylation by method B must be done under a hood.

2.6. Alternatives to methylation Method B

2.6.1. Method C

2.6.1.1. Principle

The fatty matter undergoing analysis is treated with methanol-hydrochloric acid, in a sealed vial, at 100 °C.

2.6.1.2. Apparatus

- Strong glass vial of a capacity of about 5 ml (height 40 to 45 mm, diameter 14 to 16 mm).
- 1 and 2 ml graduated pipettes.

2.6.1.3. Reagents

Solution of hydrochloric acid in 2 % methanol. This is prepared from gaseous hydrochloric acid and anhydrous methanol (Note 1).

Hexane, chromatographic quality.

Note 1:

Commercial solutions of hydrogen chloride in methanol can be used. Small amounts of gaseous hydrochloric acid can easily be prepared in the laboratory by simple displacement from the commercial solution ($p = 1,18$) by dripping concentrated sulphuric acid. Since hydrochloric acid is very rapidly absorbed by methanol, it is advisable to take the usual precautions when dissolving it, e.g. introduce the gas through a small inverted funnel with the rim just touching the surface of the liquid. Large quantities of methanolic hydrochloric acid solution can be prepared in advance, as it keeps perfectly in glass-stoppered bottles stored in the dark. Alternatively, this reagent can be prepared by dissolution of acetyl chloride in anhydrous methanol.

2.6.1.4. Procedure

- Place in the glass vial 0,2 g of the fatty matter, which has previously been dried out on sodium sulphate and filtered, and 2 ml of hydrochloric acid-methanol solution. Heat seal the vial.
- Immerse the vial at 100 °C for 40 minutes.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- Cool the vial under running water, open, add 2 ml of distilled water and 1 ml of hexane.
- Centrifuge and remove the hexane phase, which is ready for use.

2.6.2. Method D

2.6.2.1. Principle

The fatty matter undergoing analysis is heated under reflux with methanol-hexane-sulphuric acid. The methyl esters obtained are extracted with petroleum ether.

2.6.2.2. Apparatus

- Test tube of a capacity of about 20 ml, fitted with an air reflux condenser approximately 1 m in length, with ground glass joints.
- 5 ml graduated pipette.
- 50 ml separating funnel.
- 10 ml and 25 ml measuring beakers.
- 15 ml test tube with conical base.

2.6.2.3. Reagents

- Methylation reagent: anhydrous methanol-hexane-concentrated sulphuric acid (p = 1,84) in the ratio 75:25:1 (V/V/V).
- 40 to 60 °C petroleum ether.
- Anhydrous sodium sulphate.

2.6.2.4. Procedure

Place 0,1 g of oil in the 20 ml test tube and add 5 ml of methylation reagent.

Fit the reflux condenser and heat for 30 minutes in a boiling water bath (Note 2).

Transfer quantitatively the mixture into a 50 ml separating funnel, with the aid of 10 ml distilled water and 10 ml petroleum ether. Shake vigorously, and allow the phases to separate, remove the aqueous phase and wash the ether layer twice with 20 ml distilled water. Add to the separating funnel a small quantity of anhydrous sodium sulphate, shake, allow to settle for a few minutes and filter, collecting the filtrate in a 15 ml test tube with a conical base.

Evaporate the solvent over a water bath in a current of nitrogen.

Note 2:

To control boiling, insert a glass rod into the test tube and limit the temperature of the water bath to 90 °C.

3. Precision parameters

The statistical evaluation of the precision of methods A and B was published by the International Olive Oil Council in its method COI/T.20/CO. No 24.

RECOMMENDATIONS FOR GAS CHROMATOGRAPHIC ANALYSIS OF THE FATTY ACID ESTERS FROM OLIVE OIL AND OLIVE-POMACE OIL

1. Procedure

The gas chromatographic analysis of solutions of fatty esters in heptane is to be carried out according to standard ISO-5508 using a capillary column (50 m length × 0,25 or 0,32 mm i.d.) impregnated with cyanopropylsilicone phase as indicated for the determination of fatty acid trans-isomers (COI/T.20/Doc. no. 17).

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Figure 1 gives the typical gas chromatographic profile of an olive-pomace oil containing methyl and ethyl esters of fatty acids, and trans-isomers of methyl esters.

2. Calculations

2.1. For the calculation of the fatty acid composition and Δ ECN42, all the following fatty acids will be taken into account:

Myristic (C14:0).

Palmitic (C16:0). Sum of the areas of the peaks corresponding to the methyl and ethyl esters.

Palmitoleic (C16:1). Sum of the areas of the peaks corresponding to the ω 9 and ω 7 isomers of the methyl ester.

Margaric (C17:0).

Margaroleic (C17:1).

Stearic (C18:0).

Oleic (C18:1). Sum of the areas of the peaks corresponding to the ω 9 and ω 7 isomers of the methyl ester, ethyl ester, and trans-isomers of the methyl ester.

Linoleic (C18:2). Sum of the areas of the peaks corresponding to the methyl and ethyl esters, and the trans-isomers of the methyl ester.

Arachidic (C20:0).

Linolenic (C18:3). Sum of the areas of the methyl ester and the trans-isomers of the methyl ester.

Eicosenoic (C20:1).

Behenic (C22:0).

Lignoceric (C24:0).

Squalene will not be taken into account for the calculation of the total area.

2.2. For the calculation of the percentage of trans-C18:1 the peak corresponding to the methyl esters of this fatty acid is to be used. For the sum [trans-C18:2 + trans-C18:3], all the peaks corresponding to the trans-isomers of these two fatty acids are to be added together. For the calculation of the total area, all the peaks mentioned in 2.1. are to be taken into account (see COI/T.20/Doc. No. 17).


The calculation of the percentage of each fatty acid will be carried out according to the formula:

$$\% X = (\text{Area } X \times 100) / (\text{total area})$$

]

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

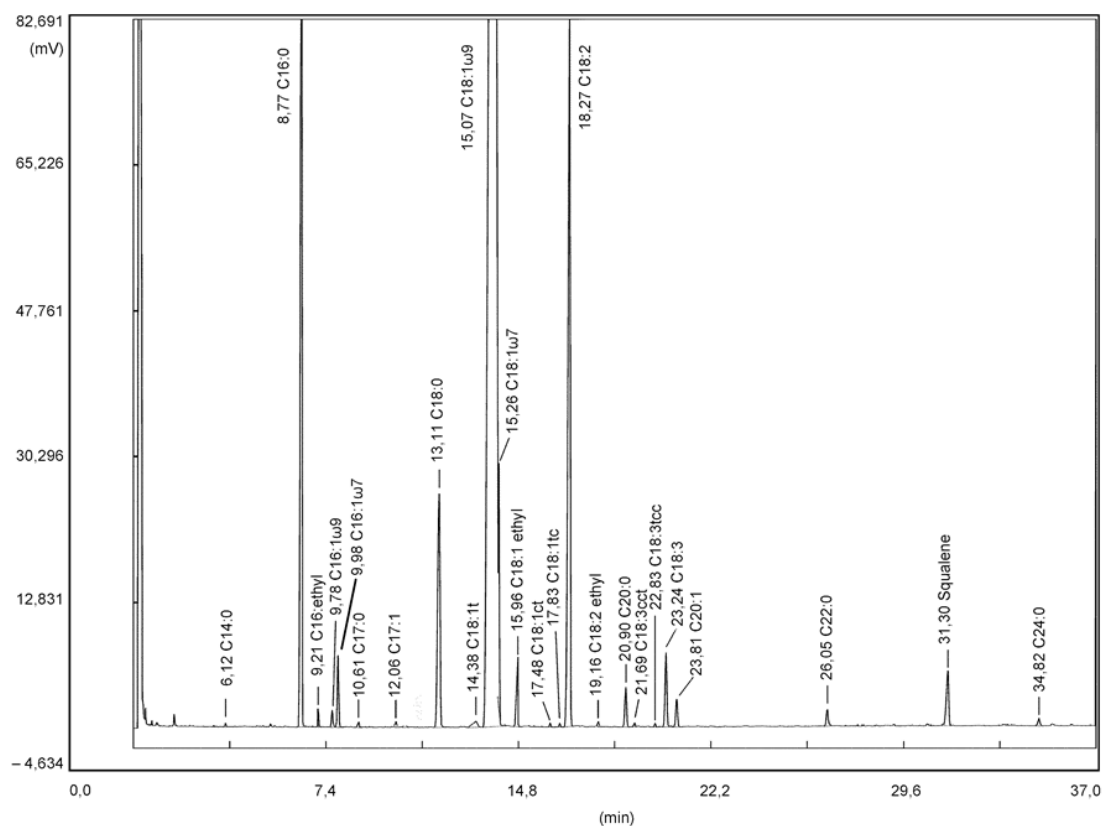

ANNEX XI

DETERMINATION OF VOLATILE HALOGENATED SOLVENTS CONTENT OF OLIVE OIL

1. METHOD

Analysis by gas chromatography using the head space technique.

2. EQUIPMENT

- 2.1. Gas chromatography apparatus fitted with an electron capture detector (ECD).
- 2.2. Head space apparatus.
- 2.3. Gas chromatography column, of glass, 2 m long and 2 mm in diameter, stationary phase. OV101 10 % or equivalent, impregnating a calcined diatomaceous earth, acid washed and silanised and of a particle size of 80 to 100 mesh.
- 2.4. Carrier and auxiliary gas: nitrogen for gas chromatography, suitable for detection by electron capture.
- 2.5. Glass flasks, 10 to 15 ml, with teflon coating and aluminium stopper with fitment for entry of syringe.
- 2.6. Hermetically sealing clamps.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

2.7. Gas syringe 0,5 to 2 ml.

3. REAGENTS

Standard: halogenated solvents of a degree of purity suitable for gas chromatography.

4. PROCEDURE

4.1. Exactly weigh around 3 g of oil in a glass flask (not to be reused); hermetically seal it. Place it in a thermostat at 70 °C for one hour. Using a syringe carefully remove 0,2 to 0,5 ml of the head space. Inject this into the column of the gas chromatography apparatus regulated as follows:

- injector temperature: 150 °C,
- column temperature: 70 to 80 °C,
- detector temperature: 200 to 250 °C.

Other temperatures may also be used provided the results remain equivalent.

4.2. Reference solutions: prepare standard solutions using refined olive oil with no trace of solvents with concentrations ranging from 0,05 to 1 ppm (mg/kg) and corresponding to the presumed content of the sample. The halogenated solvents may be diluted using pentane.

4.3. Quantitative assessment: correlate the surfaces or the elevations of the peaks of the sample and of the standard solution of the concentration presumed closest. If the deviation is greater than 10 % the analysis must be repeated in comparison with another standard solution until the deviation is within 10 %. The content is determined on the basis of the average of the elementary injections.

4.4. Expression of results: in ppm (mg/kg). The detection limit for the method is 0,01 mg/kg.

[F²ANNEX XII

THE INTERNATIONAL OLIVE COUNCIL'S METHOD FOR THE ORGANOLEPTIC ASSESSMENT OF VIRGIN OLIVE OIL

1. PURPOSE AND SCOPE

The purpose of this international method is to determine the procedure for assessing the organoleptic characteristics of virgin olive oil within the meaning of point 1 of Annex XVI to Regulation (EC) No 1234/2007 and to establish the method for its classification on the basis of those characteristics. It also provides indications for optional labelling.

The method described is applicable only to virgin olive oil and to the classification or labelling of such oils according to the intensity of the defects perceived and of the fruitiness, as determined by a group of tasters selected, trained and monitored as a panel.

It also provides for indications for optional labelling.

The IOC standards mentioned in this Annex are used in their last available version.

2. GENERAL BASIC VOCABULARY FOR SENSORY ANALYSIS

Refer to the standard IOC/T.20/Doc. No 4 "Sensory Analysis: General Basic Vocabulary"

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

3. SPECIFIC VOCABULARY

3.1. Negative attributes

Fusty/muddy sediment: Characteristic flavour of oil obtained from olives piled or stored in such conditions as to have undergone an advanced stage of anaerobic fermentation, or of oil which has been left in contact with the sediment that settles in underground tanks and vats and which has also undergone a process of anaerobic fermentation.

Musty-humid-earthy: Characteristic flavour of oils obtained from fruit in which large numbers of fungi and yeasts have developed as a result of its being stored in humid conditions for several days or of oil obtained from olives that have been collected with earth or mud on them and which have not been washed.

Winey-vinegary-acid-sour: Characteristic flavour of certain oils reminiscent of wine or vinegar. This flavour is mainly due to a process of aerobic fermentation in the olives or in olive paste left on pressing mats which have not been properly cleaned and leads to the formation of acetic acid, ethyl acetate and ethanol.

Rancid: Flavour of oils which have undergone an intense process of oxidation.

Frostbitten olives (wet wood): Characteristic flavour of oils extracted from olives which have been injured by frost while on the tree.

3.2. Other negative attributes

Heated or: Characteristic flavour of oils caused by excessive and/or prolonged

Burnt: Heating during processing, particularly when the paste is thermally mixed, if this is done under unsuitable thermal conditions.

Hay-wood: Characteristic flavour of certain oils produced from olives that have dried out.

Rough: Thick, pasty mouth sensation produced by certain old oils.

Greasy: Flavour of oil reminiscent of that of diesel oil, grease or mineral oil.

Vegetable water: Flavour acquired by the oil as a result of prolonged contact with vegetable water which has undergone fermentation processes.

Brine: Flavour of oil extracted from olives which have been preserved in brine.

Metallic: Flavour that is reminiscent of metals. It is characteristic of oil which has been in prolonged contact with metallic surfaces during crushing, mixing, pressing or storage.

Esparto: Characteristic flavour of oil obtained from olives pressed in new esparto mats. The flavour may differ depending on whether the mats are made of green esparto or dried esparto.

Grubby: Flavour of oil obtained from olives which have been heavily attacked by the grubs of the olive fly (*Bactrocera oleae*).

Cucumber: Flavour produced when an oil is hermetically packed for too long, particularly in tin containers, and which is attributed to the formation of 2,6 nonadienal.

3.3. Positive attributes

Fruity: Set of olfactory sensations characteristic of the oil which depends on the variety and comes from sound, fresh olives, either ripe or unripe. It is perceived directly and/or through the back of the nose.

Bitter: Characteristic primary taste of oil obtained from green olives or olives turning colour. It is perceived in the circumvallate papillae on the 'V' region of the tongue.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Pungent: Biting tactile sensation characteristic of oils produced at the start of the crop year, primarily from olives that are still unripe. It can be perceived throughout the whole of the mouth cavity, particularly in the throat.

3.4. **Optional terminology for labelling purposes**

Upon request, the panel leader may certify that the oils which have been assessed comply with the definitions and ranges corresponding to the following adjectives according to the intensity and perception of the attributes.

Positive attributes (fruity, bitter and pungent): According to the intensity of perception:

- *Intense*, when the median of the attribute is more than 6;
- *Medium*, when the median of the attribute is between 3 and 6;
- *Light*, when the median of the attribute is less than 3.

Fruity: Set of olfactory sensations characteristic of the oil which depends on the variety of olive and comes from sound, fresh olives in which neither green nor ripe fruitiness predominates. It is perceived directly and/or through the back of the nose.

Greenly fruity: Set of olfactory sensations characteristic of the oil which is reminiscent of green fruit, depends on the variety of olive and comes from green, sound, fresh olives. It is perceived directly and/or through the back of the nose.

Ripely fruity: Set of olfactory sensations characteristic of the oil which is reminiscent of ripe fruit, depends on the variety of olive and comes from sound, fresh olives. It is perceived directly and/or through the back of the nose.

Well balanced: Oil which does not display a lack of balance, by which is meant the olfactory–gustatory and tactile sensation where the median of the bitter and/or pungent attributes is two points higher than the median of the fruitiness.

Mild oil: Oil for which the median of the bitter and pungent attributes is 2 or less.

4. GLASS FOR OIL TASTING

Refer to the standard IOC/T.20/Doc. No 5, "Glass for Oil Tasting".

5. TEST ROOM

Refer to the standard IOC/T.20/Doc. No 6, "Guide for the Installation of a Test Room".

6. ACCESSORIES

The following accessories, which are required by tasters to perform their task properly, must be supplied in each booth and must be within easy reach:

- glasses (standardised) containing the samples, code numbered, covered with a watch-glass and kept at $28\text{ °C} \pm 2\text{ °C}$;
- profile sheet (see Figure 1) on hard copy, or on soft copy provided that the conditions of the profile sheet are met, together with the instructions for its use if necessary
- pen or indelible ink
- trays with slices of apple and/or water, carbonated water and/or rusks
- glass of water at ambient temperature
- sheet recalling the general rules listed in sections 8.4 and 9.1.1
- spittoons.

7. PANEL LEADER AND TASTERS

7.1. **Panel leader**

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The panel leader must be a suitably trained person with an expert knowledge of the kinds of oils which he or she will come across in the course of their work. They are the key figure in the panel and responsible for its organisation and running.

The work of the panel leader calls for basic training in the tools of sensory analysis, sensory skill, meticulousness in the preparation, organisation and performance of the tests and skill and patience to plan and execute the tests in a scientific manner.

They are the sole person responsible for selecting, training and monitoring the tasters in order to ascertain their level of aptitude. They are thus responsible for the appraisal of the tasters, which must always be objective and for which they must develop specific procedures based on tests and solid acceptance and rejection criteria. See standard IOC/T.20/Doc. No 14, "Guide for the selection, training and monitoring of skilled virgin olive oil tasters".

Panel leaders are responsible for the performance of the panel and hence for its evaluation, of which they must give reliable, objective proof. In any case, they must demonstrate at all times that the method and tasters are under control. Periodic calibration of the panel is recommended (IOC/T.20/Doc. No 14, § 5).

They hold ultimate responsibility for keeping the records of the panel. These records must always be traceable. They must comply with the assurance and quality requirements laid down in international sensory analysis standards and ensure the anonymity of the samples at all times.

They shall be responsible for inventorying and ensuring that the apparatus and equipment needed to comply with the specifications of this method is properly cleaned and maintained and shall keep written proof thereof, as well as of the compliance with the test conditions.

They shall be in charge of the reception and storage of the samples upon their arrival at the laboratory as well as of their storage after being tested. When doing so, they shall ensure at all times that the samples remain anonymous and are properly stored, for which purpose they must develop written procedures in order to ensure that the entire process is traceable and affords guarantees.

In addition, they are responsible for preparing, coding and presenting the samples to the tasters according to an appropriate experimental design in line with pre-established protocols, as well as for assembling and statistically processing the data obtained by the tasters.

They shall be in charge of developing and drafting any other procedures that might be necessary to complement this standard and to ensure that the panel functions properly.

They must seek ways of comparing the results of the panel with those obtained by other panels undertaking the analysis of virgin olive oil in order to ascertain whether the panel is working properly.

It is the duty of the panel leader to motivate the panel members by encouraging interest, curiosity and a competitive spirit among them. To do so, they are strongly recommended to ensure a smooth two-way flow of information with the panel members by keeping them informed about all the tasks they carry out and the results obtained. In addition, they shall ensure that their opinion is not known and shall prevent possible leaders from asserting their criteria over the other tasters.

They shall summon the tasters sufficiently in advance and shall answer any queries regarding the performance of the tests, but shall refrain from suggesting any opinion to them on the sample.

7.2. Tasters

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The people acting as tasters in organoleptic tests carried out on olive oils must do so voluntarily, with all the ensuing consequences of such a voluntary act in terms of obligations and the absence of financial payment. It is therefore advisable for candidates to submit an application in writing. Candidates shall be selected, trained and monitored by the panel leader in accordance with their skills in distinguishing between similar samples; it should be borne in mind that their accuracy will improve with training.

Tasters must act like real sensory observers, setting aside their personal tastes and solely reporting the sensations they perceive. To do so, they must always work in silence, in a relaxed, unhurried manner, paying the fullest possible sensory attention to the sample they are tasting.

Between 8 and 12 tasters are required for each test, although it is wise to keep some extra tasters in reserve to cover possible absences.

8. TEST CONDITIONS

8.1. Presentation of the sample

The oil sample for analysis shall be presented in standardised tasting glasses conforming to the standard IOC/T.20/Doc. No 5 'Glass for oil tasting'.

The glass shall contain 14–16 ml of oil, or between 12,8 and 14,6 g if the samples are to be weighed, and shall be covered with a watch-glass.

Each glass shall be marked with a code made up of digits or a combination of letters and digits chosen at random. The code will be marked by means of an odourfree system.

8.2. Test and sample temperature

The oil samples intended for tasting shall be kept in the glasses at $28\text{ °C} \pm 2\text{ °C}$ throughout the test. This temperature has been chosen because it makes it easier to observe organoleptic differences than at ambient temperature and because at lower temperatures the aromatic compounds peculiar to these oils volatilise poorly while higher temperatures lead to the formation of volatile compounds peculiar to heated oils. See the standard IOC/T.20/Doc. No 5 'Glass for Oil Tasting' for the method which has to be used for heating the samples when in the glass.

The test room must be at a temperature between 20 °C and 25 °C (see IOC/T.20/Doc. No 6).

8.3. Test times

The morning is the best time for tasting oils. It has been proved that there are optimum perception periods as regards taste and smell during the day. Meals are preceded by a period in which olfactory–gustatory sensitivity increases, whereas afterwards this perception decreases.

However, this criterion should not be taken to the extreme where hunger may distract the tasters, thus decreasing their discriminatory capacity; therefore, it is recommended to hold the tasting sessions between 10.00 in the morning and 12 noon.

8.4. Tasters: general rules of conduct

The following recommendations apply to the conduct of the tasters during their work.

When called by the panel leader to participate in an organoleptic test, tasters should be able to attend at the time set beforehand and shall observe the following:

- They shall not smoke or drink coffee at least 30 minutes before the time set for the test.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- They must not have used any fragrance, cosmetic or soap whose smell could linger until the time of the test. They must use an unperfumed soap to wash their hands which they shall then rinse and dry as often as necessary to eliminate any smell.
- They shall fast at least one hour before the tasting is carried out.
- Should they feel physically unwell, and in particular if their sense of smell or taste is affected, or if they are under any psychological effect that prevents them from concentrating on their work, the tasters shall refrain from tasting and shall inform the panel leader accordingly.
- When they have complied with the above, the tasters shall take up their place in the booth allotted to them in an orderly, quiet manner.
- They shall carefully read the instructions given on the profile sheet and shall not begin to examine the sample until fully prepared for the task they have to perform (relaxed and unhurried). If any doubts should arise, they should consult the panel leader in private.
- They must remain silent while performing their tasks.
- They must keep their mobile phone switched off at all times to avoid interfering with the concentration and work of their colleagues.

9. PROCEDURE FOR THE ORGANOLEPTIC ASSESSMENT AND CLASSIFICATION OF VIRGIN OLIVE OIL

9.1. Tasting technique

- 9.1.1. The tasters shall pick up the glass, keeping it covered with the watch-glass, and shall bend it gently; they shall then rotate the glass fully in this position so as to wet the inside as much as possible. Once this stage is completed, they shall remove the watch-glass and smell the sample, taking slow deep breaths to evaluate the oil. Smelling should not exceed 30 seconds. If no conclusion has been reached during this time, they shall take a short rest before trying again.

When the olfactory test has been performed, the tasters shall then evaluate the buccal sensations (overall retronasal olfactory, gustatory and tactile sensations). To do so, they shall take a small sip of approximately 3 ml of oil. It is very important to distribute the oil throughout the whole of the mouth cavity, from the front part of the mouth and tongue along the sides to the back part and to the palate support and throat, since it is a known fact that the perception of tastes and tactile sensations varies in intensity depending on the area of the tongue, palate and throat.

It should be stressed that it is essential for a sufficient amount of the oil to be spread very slowly over the back of the tongue towards the palate support and throat while the taster concentrates on the order in which the bitter and pungent stimuli appear. If this is not done, both of these stimuli may escape notice in some oils or else the bitter stimulus may be obscured by the pungent stimulus.

Taking short, successive breaths, drawing in air through the mouth, enables the taster not only to spread the sample extensively over the whole of the mouth but also to perceive the volatile aromatic compounds via the back of the nose by forcing the use of this channel.

The tactile sensation of pungency should be taken into consideration. For this purpose it is advisable to ingest the oil.

- 9.1.2. When organoleptically assessing a virgin olive oil, it is recommended that FOUR SAMPLES at the most be evaluated in each session with a maximum of three sessions per day, to avoid the contrast effect that could be produced by immediately tasting other samples.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

As successive tastings produce fatigue or loss of sensitivity caused by the preceding samples, it is necessary to use a product that can eliminate the remains of the oil from the preceding tasting from the mouth.

The use of a small slice of apple is recommended which, after being chewed, can be disposed of in the spittoon. Then rinse out the mouth with a little water at ambient temperature. At least 15 minutes shall lapse between the end of one session and the start of the next.

9.2. Use of the profile sheet by tasters

The profile sheet intended for use by tasters is detailed in Figure 1 of this Annex.

Each taster on the panel shall smell and then taste⁽⁸⁾ the oil under consideration. They shall then enter the intensity with which they perceive each of the negative and positive attributes on the 10-cm scale shown in the profile sheet provided.

Should the tasters perceive any negative attributes not listed in section 4, they shall record them under the "others" heading, using the term or terms that most accurately describes the attributes.

9.3. Use of the data by the panel leaders

The panel leader shall collect the profile sheets completed by each taster and shall review the intensities assigned to the different attributes. Should they find any anomaly, they shall invite the taster to revise his or her profile sheet and, if necessary, to repeat the test.

The panel leader shall enter the assessment data of each panel member in a computer program like that provided by the standard IOC/T.20/Doc. No 15) with a view to statistically calculating the results of the analysis, based on the calculation of their median. See sections 9.4 and Appendix to this Annex. The data for a given sample shall be entered with the aid of a matrix comprising 9 columns representing the 9 sensory attributes and n lines representing the n panel members used.

When a defect is perceived and entered under the "others" heading by at least 50 % of the panel, the panel leader shall calculate the median of the defect and shall arrive at the corresponding classification.

The value of the robust coefficient of variation which defines classification (defect with the strongest intensity and fruity attribute) must be no greater than 20 %.

If the opposite is the case, the panel leader must repeat the evaluation of the specific sample in another tasting session.

If this situation arises often, the panel leader is recommended to give the tasters specific additional training (IOC/T.20/Doc. No 14, § 5) and to use the repeatability index and deviation index to check panel performance (IOC/T.20/Doc. No 14, § 6).

9.4. Classification of the oil

The oil is graded as follows in line with the median of the defects and the median for the fruity attribute. The median of the defects is defined as the median of the defect perceived with the greatest intensity. The *median* of the defects and the median of the fruity attribute are expressed to one decimal place.

The oil is graded by comparing the median value of the defects and the median for the fruity attribute with the reference ranges given below. The error of the method has been taken into account when establishing the limits of these ranges, which are therefore considered to be absolute. The software packages allow the grading to be displayed as a table of statistics or a graph.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- (a) Extra virgin olive oil: the median of the defects is 0 and the median of the fruity attribute is above 0;
- (b) Virgin olive oil: the median of the defects is above 0 but not more than 3,5 and the median of the fruity attribute is above 0;
- (c) Lampante olive oil: the median of defect is above 3,5 or the median of the defect is less than or equal to 3,5 and the fruity median is equal to 0.

Note 1:

When the median of the bitter and/or pungent attribute is more than 5,0, the panel leader shall state so on the test certificate.

FIGURE 1

PROFILE SHEET FOR VIRGIN OLIVE OIL

Intensity of perception of defects	
Fusty/muddy sediment ^a	
Musty/humid/earthy ^a	
Winey/vinegary acid/sour ^a	
Frostbitten olives (wet wood)	
Rancid	
Other negative attributes:	
Descriptor:	Metallic # Hay # Grubby # Rough # Brine # Heated or burnt # Vegetable water # Esparto # Cucumber # Greasy #

Intensity of perception of positive attributes	
Fruity	
	Green # Ripe #
Bitter	
Pungent	
Name of taster:	Taster code:
Sample code:	Signature:

^a Delete as appropriate

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix

Method for calculating the median and the confidence intervals

Median

$$Me = \left[p(X < x_m) \leq \frac{1}{2} \wedge p(X \leq x_m) \geq \frac{1}{2} \right]$$

The median is defined as the real number X_m characterised by the fact that the probability (p) that the distribution values (X) are below this number (X_m), is less than and equal to 0,5 and that simultaneously the probability (p) that the distribution values (X) are below or equal to X_m is greater than and equal to 0,5. A more practical definition is that the median is the 50th percentile of a distribution of numbers arranged in increasing order. In simpler terms, it is the midpoint of an ordered set of odd numbers, or the mean of two midpoints of an ordered set of even numbers.

Robust standard deviation

In order to arrive at a reliable estimate of the variability around the mean it is necessary to refer to the robust standard deviation as estimated according to Stuart and Kendall (4). The formula gives the asymptotic robust standard deviation, i.e. the robust estimate of the variability of the data considered where N is the number of observations and IQR is the interquartile range which encompasses exactly 50% of the cases of a given probability distribution:

$$s^* = \frac{1,25 \times \text{IQR}}{1,35 \times \sqrt{N}}$$

The interquartile range is calculated by calculating the magnitude of the difference between the 75th and 25th percentile.

$$\text{IQR} = 75\text{th percentile} - 25\text{th percentile}$$

Where the percentile is the value X_{pc} characterised by the fact that the probability (p) that the distribution values are less than X_{pc} is less than and equal to a specific hundredth and that simultaneously the probability (p) that the distribution values are less than or equal to X_{pc} is greater than and equal to that specific hundredth. The hundredth indicates the distribution fractile chosen. In the case of the median it is equal to 50/100.

$$\text{percentile} = \left[p(X < x_{pc}) \leq \frac{n}{100} \wedge p(X \leq x_{pc}) \geq \frac{n}{100} \right]$$

For practical purposes, the percentile is the distribution value corresponding to a specific area subtended from the distribution or density curve. To give an example, the 25th percentile represents the distribution value corresponding to an area equal to 0,25 or 25/100.

In this method percentiles are computed on the basis of the real values which appear in the data matrix (percentiles computing procedure).

Robust coefficient of variation (%)

The $CV_r\%$ represents a pure number which indicates the percentage variability of the set of numbers analysed. For this reason it is very useful for checking the reliability of the panel assessors.

$$CV_r = \frac{s^*}{Me} \times 100$$

Confidence intervals of the median at 95%

The confidence intervals at 95% (value of the error of the first kind equal to 0,05 or 5%) represent the interval within which the value of the median could vary if it were possible to repeat an experiment an infinite number of times. In practice, it indicates the interval of variability of the

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

test in the operating conditions adopted starting from the assumption that it is possible to repeat it many times. As with the $CVr\%$, the interval helps to assess the reliability of the test.

$$C.I._{upper} = Me + (c \times s^*)$$

$$C.I._{lower} = Me - (c \times s^*)$$

where $C = 1,96$ for the confidence interval at the 95% level.

An example of the calculation sheet is presented in Annex I to the standard IOC/T 20/Doc. No 15.

References

- (1) Wilkinson, L. 1990. Systat: The system for statistics. Evanston, IL.SYSTAT Inc.
- (2) Cicchitelli, G. 1984. Probabilità e Statistica. Maggioli Editore, Rimini.
- (3) Massart, D.L.; Vandeginste, B.G.M.; Deming, Y.; Michotte, L. 1988. Chemometrics. A textbook. Elsevier. Amsterdam.
- (4) Kendall, M.G.; Stuart, A. 1967. The advanced theory of statistics. Vol. 1. Hafner Publishing Co.
- (5) McGill, R.; Tukey, J.W.; Larsen, W.A. 1978. Variation of Box Plots. The American Statistician, 32, (2), 12-16.
- (6) IOC/T.28/Doc. No 1 September 2007, Guidelines for the accreditation of sensory testing laboratories with particular reference to virgin olive oil according to standard ISO/IEC 17025:2005.
- (7) IOC/T.20/Doc. No 14.
- (8) IOC/T.20/Doc. No 15.
- (9) ISO/IEC 17025:05.]

F¹³ ANNEX XIII

.....

F⁴ ANNEX XIV

.....

ANNEX XV

1. OIL CONTENT OF OLIVE RESIDUE
 - 1.1. Apparatus
 - suitable extraction apparatus fitted with a 200 to 250 ml round-bottomed flask,

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- electrically heated bath (e.g., sand bath, water bath) or hotplate,
- analytical balance,
- oven regulated to a maximum of 80° C,
- electrically heated oven fitted with a thermostatic device regulated to 103 ± 2° C and one that can be swept with a stream of air or operated at reduced pressure,
- mechanical mill, easy to clean, and one that allows the olive residues to be ground without a rise in their temperature or any appreciable alteration in their content of moisture, volatile matter or substances extractable with hexane,
- extraction thimble and cotton wool or filter paper from which substances extractable with hexane have already been removed,
- dessicator,
- sieve with 1 mm diameter apertures,
- small particles of previously dried pumice stone.

1.2. Reagent

Normal hexane, technical grade, which must leave a residue of less than 0,002 g per 100 ml, on complete evaporation.

2. PROCEDURE

2.1. Preparation of the test sample

If necessary, use the mechanical mill, which has previously been properly cleaned, to grind the laboratory sample in order to reduce it to particles that can pass completely through the sieve.

Use about one twentieth of the sample to complete the process of cleaning the mill, discard the ground material, grind the remainder and collect, mix carefully and analyze without delay.

2.2. Test portion

As soon as the grinding operation has been completed, weigh out about 10 g of the sample to the nearest 0,01 g for testing.

2.3. Preparation of the extraction thimble

Place the test portion in the thimble and plug with cotton wool. If a filter paper is used, envelope the test portion in it.

2.4. Preliminary drying

If the olive residues are very moist (i.e., moisture and volatile matter content more than 10 %), carry out preliminary drying by placing the loaded thimble (or filter paper) in the oven heated for an appropriate time at not more than 80° C in order to reduce the moisture and volatile matter content to less than 10 %.

2.5. Preparation of the round-bottomed flask

Weigh to the nearest 1 mg the flask containing one or two particles of pumice stone, previously dried in the stove at 103 ± 2° C and then cooled in a dessicator for not less than one hour.

2.6. Initial extraction

Into the extraction apparatus insert the thimble (or filter paper) containing the test portion. Pour into the flask the requisite quantity of hexane. Fit the flask to the extraction apparatus and place the whole on the electrically heated bath. Adjust the rate of heating in such a way that the reflux rate is not less than three drops per second (moderate, not violent boiling). After four hours

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

extraction, allow to cool. Remove the thimble from the extraction apparatus and place it in a stream of air in order to drive off most of the impregnating solvent.

2.7. Second extraction

Tip the contents of the thimble into the micro-grinder and grind as finely as possible. Return the ground mixture to the thimble without loss and place it back in the extraction apparatus.

Continue the extraction for a further two hours using the same round-bottomed flask containing the initial extract.

The resultant solution in the extraction flask must be clear. If not, filter it through a filter paper and wash the original flask and the filter paper several times with hexane. Collect the filtrate and the washing solvent in a second round-bottomed flask which has been dried and tared to the nearest 1 mg.

2.8. Removal of solvent and weighing of extract

Remove the greater part of the solvent by distillation on an electrically heated bath. Remove the last traces of solvent by heating the flask in the oven at $103 \pm 2^\circ \text{C}$ for 20 minutes. Assist the elimination process either by blowing in air, or preferably an inert gas, at intervals or by using reduced pressure.

Leave the flask in a dessicator to cool for at least one hour and weigh to the nearest 1 mg.

Heat again for 10 minutes under the same conditions, cool in a dessicator and reweigh.

The difference between the two weighings shall not exceed 10 mg. If it does, heat again for periods of 10 minutes followed by cooling and weighing until the weight difference is 10 mg or less. Note the last weight of the flask.

Carry out duplicate determinations on the test sample.

3. EXPRESSION OF RESULTS

3.1. Method of calculation and formula

(a) The extract expressed as a percentage by mass of the product as received is equal to:

$$S = m_1 \times \frac{100}{m_0}$$

where:	S	= is the percentage by mass of extract of the product as received,
	m0	= is the mass, in grams, of the test portion,
	m1	= is the mass, in grams, of the extract after drying.

Take as the result the arithmetic mean of the duplicate determinations, providing the repeatability conditions are satisfied.

Express the result to the first decimal place.

(b) The extract is expressed on a dry matter basis by using the formula:

$$S \times \frac{100}{100-U} = \text{oil percentage of extract on a dry basis}$$

where:

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- S = is the percentage of extract by means of the product as received (see (a)),
- U = is its moisture and volatile matter content.

3.2. Repeatability

The difference between the duplicate determinations carried out simultaneously or in rapid succession by the same analyst shall not exceed 0,2 g of hexane extract per 100 g of sample.

If this condition is not satisfied, repeat the analysis on two other test portions. If, in this case too, the difference exceeds 0,2 g, take as the result the arithmetic mean of the four determinations.

ANNEX XVI

DETERMINATION OF IODINE VALUE

1. SCOPE

This International Standard specifies a method for the determination of the iodine value of animal and vegetable fats and oils, referred to hereafter as fats.

2. DEFINITION

For the purposes of this International Standard, the following definition applies:

- 2.1. *iodine value*. The mass of iodine absorbed by the sample under the operating conditions specified in this International Standard.

The iodine value is expressed as grams of iodine per 100 g of sample.

3. PRINCIPLE

Dissolution of a test portion in solvent and addition of Wijs reagent. After a specified time, addition of potassium iodide solution and water, and titration of the liberated iodine with sodium thiosulfate solution.

4. REAGENTS

All reagents shall be of recognized analytical grade:

- 4.1. *water*, complying with the requirements of ISO 3696, Grade 3.
- 4.2. *potassium iodide*, 100 g/l solution, not containing iodate or free iodine.
- 4.3. *starch*, solution.

Mix 5 g of soluble starch in 30 ml of water, add this mixture to 1 000 ml of boiling water, boil for three minutes and allow to cool.

- 4.4. *sodium thiosulfate*, standard volumetric solution $c(\text{Na}_2\text{S}_2\text{O}_3 \cdot 5\text{H}_2\text{O}) = 0,1 \text{ mol/l}$, standardized not more than seven days before use.
- 4.5. *solvent*, prepared by mixing equal volumes of cyclohexane and acetic acid.
- 4.6. *Wijs reagent*, containing iodine monochloride in acetic acid. Commercially available Wijs reagent shall be used.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

5. APPARATUS

Usual laboratory apparatus and, in particular, the following:

- 5.1. *glass weighing scoops*, suitable for the test portion and for inserting into the flasks (6.2).
- 5.2. *conical flasks*, of 500 ml capacity, fitted with ground glass stoppers and completely dry.

6. PREPARATION OF THE TEST SAMPLE

The homogenized sample is dried over sodium sulphate and filtered.

7. PROCEDURE

7.1. Test portion

The mass of the test portion varies according to its expected iodine value as shown in Table 1.

Expected iodine value	Mass of test portion(g)
less than 5	3,0
5 to 20	1,0
21 to 50	0,4
51 to 100	0,2
101 to 150	0,13
151 to 200	0,1

Weigh the test portion to the nearest 0,1 mg in a glass weighing scoop (5.1).

7.2. Determination

Place the test portion in a 500 ml flask (6.2). Add 20 ml of the solvent (4.5) to dissolve the fat. Add exactly 25 ml of the Wijs reagent (4.6), insert the stopper, swirl the contents and place the flask in the dark. Do not use a mouth pipette for the Wijs reagent.

Similarly, prepare a blank with the solvent and the reagent but omitting the test portion.

For samples having an iodine value below 150, leave the flasks in the dark for one hour; for those with an iodine value above 150 and for polymerized products or products oxidized to a considerable extent, leave for two hours.

At the end of the time, add 20 ml of the potassium iodide solution (4.2) and 150 ml of water (4.1) to each of the flasks.

Titrate with the standard volumetric sodium thiosulfate solution (4.4) until the yellow colour due to iodine has almost disappeared. Add a few drops of the starch solution (4.3) and continue the titration until the blue colour just disappears after very vigorous shaking.

Note:

Potentiometric determination of the end point is permissible.

7.3. Number of determinations

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Carry out two determinations on the same test sample.

8. EXPRESSION OF RESULTS

The iodine value is given by the expression

$$\frac{12,69 \cdot c(V_1 - V_2)}{m}$$

where:

- c = is the numerical value of the exact concentration, in moles per litre, of the standard volumetric sodium thiosulfate solution (4.4) used;
- V₁ = is the numerical value of the volume, in millilitres, of the standard volumetric sodium thiosulfate solution (4.4) used for the blank test;
- V₂ = is the numerical value of the volume, in millilitres, of the standard volumetric sodium thiosulfate solution (4.4) used for the determination;
- m = is the numerical value of the mass, in grams, of the test portion (7.1).

Take as the result the arithmetic mean of the two determinations, provided that the requirement for repeatability (9.2) is satisfied.

[^{F18}ANNEX XVII

METHOD FOR THE DETERMINATION OF STIGMASTADIENES IN VEGETABLE OILS

Textual Amendments

F18 Inserted by Commission Regulation (EC) No 656/95 of 28 March 1995 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis and Council Regulation (EEC) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.

1. PURPOSE

Determination of stigmastadienes in vegetable oils containing low concentrations of these hydrocarbons, particularly in virgin olive oil and crude olive-residue oil.

2. SCOPE

The standard may be applied to all vegetable oils although measurements are reliable only where the content of these hydrocarbons lies between 0,01 and 4,0 mg/kg. The method is particularly suited to detecting the presence of refined vegetable oils (olive, olive residue, sunflower, palm, etc.) in virgin olive oil since refined oils contained stigmastadienes and virgin oils do not.

3. PRINCIPLE

Isolation of unsaponifiable matter. Separation of steroidal hydrocarbon fraction by column chromatography on silica gel and analysis by capillary gas chromatography.

4. APPARATUS

- 4.1. 250 ml flasks suitable for use with a reflux condenser.
- 4.2. Separating funnels of 500 ml capacity.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 4.3. 100 ml round-bottom flasks.
- 4.4. Rotary evaporator.
- 4.5. Glass chromatography column (1,5 to 2,0 cm internal diameter by 50 cm length) with Teflon tap and a plug of glass wool fibre or sintered glass disc at the bottom. To prepare silica gel column, pour hexane into the chromatography column to a depth of approximately 5 cm and then fill with a slurry of silica gel in hexane (15 g in 40 ml) with the help of hexane portions. Allow to settle and finish settling by applying slight vibration. Add anhydrous sodium sulphate to a height of approximately 0,5 cm, finally elute the excess hexane.
- 4.6. Gas chromatograph with flame ionization detector, split or cold on-column injector and oven programmable to within ± 1 °C.
- 4.7. Fused silica capillary column for gas chromatography (0,25 or 0,32 mm internal diameter by 25 m length) coated with 5 %-phenylmethylsilicone phase, 0,25 mm film thickness.

Note 1:

Other columns of similar or lower polarity can be used.

- 4.8. Integrator-recorder with possibility of valley-valley integration mode.
- 4.9. 5 to 10 ml microsyringe for gas chromatography with cemented needle.
- 4.10. Electrical heating mantle or hot place.

5. REAGENTS

All reagents should be of analytical grade unless otherwise specified. The water used should be distilled water, or water of at least equivalent purity.

- 5.1. Hexane or mixture of alkanes of bp interval 65 to 70 °C, distilled with rectifying column.

Note 2:

The solvent must be distilled to remove impurities.

- 5.2. 96 v/v ethanol.
- 5.3. Anhydrous sodium sulphate.
- 5.4. Alcoholic potassium hydroxide solution at 10 %. Add 10 ml of water to 50 g potassium hydroxide, stir, and then dissolve the mixture in ethanol to 500 ml.

Note 3:

Alcoholic potash turns brown on standing. It should be prepared freshly each day and kept in well stoppered dark glass bottles.

- 5.5. Silica gel 60 for column chromatography, 70 to 230 mesh, (Merck, reference 7734 or similar).

Note 4:

Usually, silica gel can be used directly from the container without any treatment. However, some batches of silica gel may show low activity resulting in bad chromatographic separations. Under this circumstance, the silica gel should be treated in the following way: Activate the silica gel by heating for a minimum of four hours at 550 °C. After heating, place the silica gel in a desiccator

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

while the gel is cooling and then transfer the silica gel to a stoppered flask. Add 2 % of water and shake until no lumps can be seen and the powder flows freely.

If batches of silica gel result in chromatograms with interfering peaks, the silica gel should be treated as above. An alternative could be the use of extra pure silica gel 60 (Merck, reference 7754).

5.6. Stock solution (200 ppm) of cholesta-3,5-diene (Sigma, 99 % purity) in hexane (10 mg in 50 ml).

5.7. Standard solution of cholesta-3,5-diene hexane at concentration of 20 ppm, obtained by dilution of above solution.

Note 5:

The solutions 5.6 and 5.7 are stable for a period of at least four months if kept at less than 4 °C.

5.8. Solution of n-nonacosane in hexane at concentration of approximately 100 ppm.

5.9. Carrier gas for chromatography: helium or hydrogen of 99,9990 % purity.

5.10. Auxiliary gases for flame ionization detector: hydrogen of 99,9990 % purity and purified air.

6. PROCEDURE

6.1. Preparation of unsaponifiable matter

6.1.1. Weigh $20 \pm 0,1$ g of oil into a 250-ml flask (4.1), add 1 ml of the standard solution of cholesta-3,5-diene (20 μ g) and 75 ml of alcoholic potash at 10 %, fit reflux condenser, and heat to slight boiling for 30 minutes, Remove the flask containing the sample from the heat and allow the solution to cool slightly (do not allow to cool completely as the sample will set). Add 100 ml of water and transfer the solution to a separating funnel (4.2) with the aid of 100 ml of hexane. Shake the mixture vigorously for 30 seconds and allow the separate.

Note 6:

If an emulsion is produced which does not rapidly disappear, add small quantities of ethanol.

6.1.2. Transfer the aqueous phase beneath to a second separating funnel and extract again with 100 ml of hexane. Once more run off the lower phase and wash the hexane extracts (combined in another separating funnel) three times with 100 ml each time of a mixture of ethanol-water (1: 1) until neutral pH is reached.

6.1.3. Pass the hexane solution through anhydrous sodium sulphate (50 g), wash with 20 ml hexane and evaporate in a rotary evaporator at 30 °C under reduced pressure until dryness.

6.2. Separation of steroidal hydrocarbon fraction

6.2.1. Take the residue to the fractioning column with the aid of two 1-ml portions of hexane, run the sample onto the column by allowing the solution level to drop to the top of the sodium sulphate and start the chromatographic elution with hexane at a flow rate of 1 ml/min approximately. Discard the first 25 to 30 ml of eluate and then collect the following 40 ml fraction. After collection, transfer this fraction to a 100-ml round bottomed flask (4.3).

Note 7:

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The first fraction contains saturated hydrocarbons (Figure 1 a) and the second fraction the steroidal ones. Further elution provides squalene and related compounds. To achieve a good separation between saturated and steroidal hydrocarbons, the optimization of fraction volumes is required. For this, the volume of the first fraction should be adjusted so that when the second fraction is analysed the peaks representing the saturated hydrocarbons are low (see Figure 1 c); if they do not appear but the intensity of the standard peak is low, the volume should be reduced. Anyway, a complete separation between the components of the first and second fractions is unnecessary; as there is no overlapping of peaks during GC analysis if GC conditions are adjusted as indicated in 6.3.1. The optimization of the volume of the second fraction is generally not needed as a good separation exists with the further components. Nevertheless, the presence of a large peak at approximately 1,5 minutes lower retention time than the standard is due to squalene, and it is indicative of a bad separation.

6.2.2. Evaporate the second fraction in a rotary evaporator at 30 °C under reduced pressure until dryness, and immediately dissolve the residue in 0,2 ml of hexane. Keep the solution in the refrigerator until analysis.

Note 8:

Residues 6.1.3 and 6.2.2 should not be kept dry and at room temperature. As soon as they are obtained, the solvent should be added and the solutions should be kept in the refrigerator.

6.3. Gas chromatography

6.3.1. Working conditions for split injection:

- injector temperature: 300 °C,
- detector temperature: 320 °C,
- integrator-recorder: the parameters for integration should be fixed so as to give a correct assessment of the areas. Valley-valley integration mode is recommended,
- sensitivity: about 16 times the minimum attenuation,
- amount of solution injected: 1 µl,
- oven programming temperatures: initial 235 °C for six minutes and then rising at 2 °C/minute up to 285 °C,
- injector with 1: 15 flow divider,
- carrier: helium or hydrogen at about 120 kPa pressure.

These conditions may be adjusted in accordance with the characteristics of the chromatograph and the column to give chromatograms meeting the following requirements: internal standard peak within approximately five minutes of the time given in 6.3.2; the internal standard peak should be at least 80 % of the full scale.

The gas chromatographic system must be checked injecting a mixture of the stock solution of cholestadiene (5.6) and n-nonacosane solution (5.8). The cholesta-3,5-diene peak must appear before the n-nonacosane (Figure 1c); if it does not occur two actions can be undertaken: reduce the oven temperature and/or use a less polar column.

6.3.2. Peak identification

The internal standard peak appears at approximately 19 minutes and the 3,5-stigmastadiene at a relative retention time of approximately 1,29 (see Figure 1b). The 3,5-stigmastadiene occurs with small quantities of an isomer, and usually, both elute together as a single chromatographic peak. Nevertheless, if the column is too polar or shows a high resolving power, the isomer can appear as a small peak before and close to that of stigmasta-3,5-diene (Figure 2). In order to ensure that the stigmastadienes are eluted as one peak, it is advisable to replace the column by one which is either less polar or has a wider internal diameter.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Note 9:

Stigmastadienes for reference can be obtained from the analysis of a refined vegetable oil by using less amount of sample (1 to 2 g). Stigmastadienes originate a prominent and easily identifiable peak.

6.3.3. Quantitative analysis

The stigmastadienes content is determined according to the formula:

mg/kg of stigmastadienes =

$$\frac{A_s \times M_c}{A_c \times M_o}$$

where:	A_s = area of stigmastadienes peak (if the peak is resolved into two isomers, sum of areas of the two peaks), A_c = area of internal standard (cholestadiene), M_c = mass of standard added, in micrograms, M_o = mass of oil taken, in grams.
--------	---

Detection limit: about 0,01 mg/kg.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

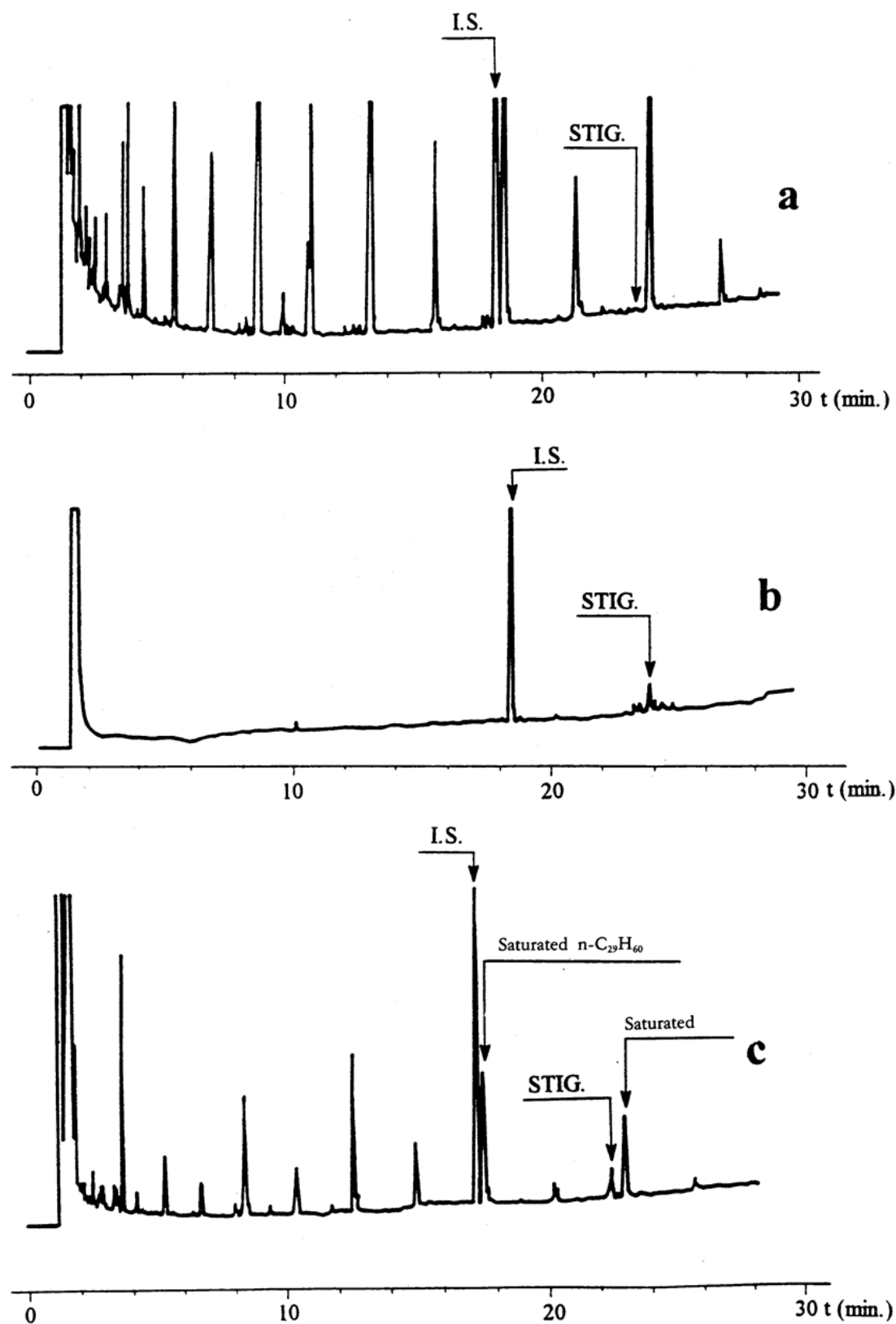



Figure 1

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Gas chromatograms obtained from olive oil samples analysed on a fused silica capillary column (0,25 mm internal diameter by 25 m) coated with 5 %-phenylmethylsilicone, 0,25 µm film thickness.

- (a) First fraction (30 ml) from a virgin oil, spiked with standard.
- (b) Second fraction (40 ml) from an olive oil containing 0,10 mg/kg of stigmastadienes.
- (c) Second fraction (40 ml) containing a small proportion of the first fraction.

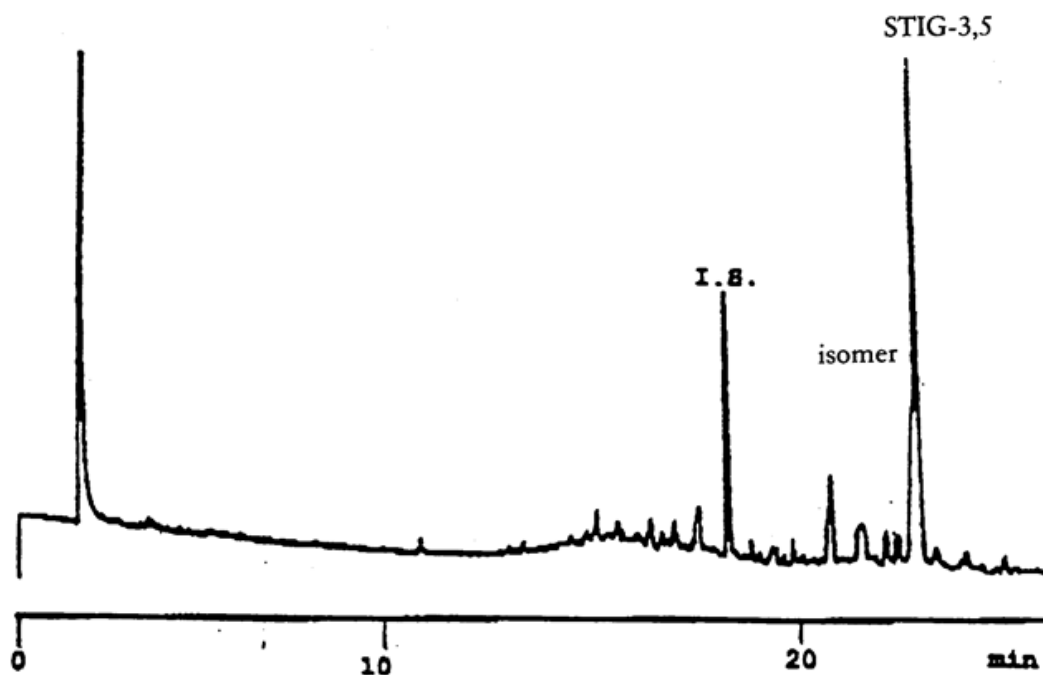

Gas chromatogram obtained from a refined olive oil sample analysed on DB-5 column showing the isomer of 3,5-stigmastadiene.]

[F³ ANNEX XVIII

DETERMINATION OF THE DIFFERENCE BETWEEN ACTUAL AND THEORETICAL CONTENT OF TRIACYLGLYCEROLS WITH ECN 42

1. SCOPE

Determination of the absolute difference between the experimental values of triacylglycerols (TAGs) with equivalent carbon number 42 (ECN₄₂^{HPLC}) obtained by determination in the oil by high performance liquid chromatography and the theoretical value of TAGs with an equivalent carbon number of 42 (ECN 42_{theoretical}) calculated from the fatty acid composition.

2. FIELD OF APPLICATION

The standard is applicable to olive oils. The method is applicable to the detection of the presence of small amounts of seed oils (rich in linoleic acid) in every class of olive oils.

3. PRINCIPLE

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The content of triacylglycerols with ECN 42 determined by HPLC analysis and the theoretical content of triacylglycerols with ECN 42 (calculated on the basis of GLC determination of fatty acid composition) correspond within a certain limit for genuine olive oils. A difference larger than the values adopted for each type of oil points out that the oil contains seed oils.

4. METHOD

The method for the calculation of the theoretical content of triacylglycerols with ECN 42 and of the difference with respect to the HPLC data is essentially made by the coordination of analytical data obtained by means of other methods. It is possible to distinguish three phases: determination of fatty acid composition by capillary gas chromatography, calculation of theoretical composition of triacylglycerols with ECN 42, HPLC determination of ECN 42 triacylglycerols.

4.1. Apparatus

- 4.1.1. Round-bottomed flasks, 250 and 500 ml.
- 4.1.2. Beakers 100 ml.
- 4.1.3. Glass chromatographic column, 21 mm internal diameter, 450 mm length, with cock and normalised cone (female) at the top.
- 4.1.4. Separating funnels, 250 ml, with normalised cone (male) at the bottom, suitable for connection to the top of the column.
- 4.1.5. Glass rod, 600 mm length.
- 4.1.6. Glass funnel, 80 mm diameter.
- 4.1.7. Volumetric flasks, 50 ml.
- 4.1.8. Volumetric flasks, 20 ml.
- 4.1.9. Rotary evaporator.
- 4.1.10. High performance liquid chromatograph, allowing thermostatic control of column temperature.
- 4.1.11. Injection units for 10 µl delivery.
- 4.1.12. Detector: differential refractometer. The full scale sensitivity should be at least 10^{-4} units of refractive index.
- 4.1.13. Column: stainless steel tube 250 mm length x 4,5 mm internal diameter packed with 5 µm diameter particles of silica with 22 to 23 % carbon in the form of octadecylsilane.
- 4.1.14. Data processing software.
- 4.1.15. Vials, of about 2 ml volumes, with Teflon-layered septa and screw caps.

4.2. Reagents

The reagents should be of analytical purity. Elution solvents should be de-gassed, and may be recycled several times without effect on the separations.

- 4.2.1. Petroleum ether 40– 60 °C chromatographic grade or hexane.
- 4.2.2. Ethyl ether, peroxide-free, freshly distilled.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 4.2.3. Elution solvent for purifying the oil by column chromatography mixture petroleum ether/ethyl ether 87/13 (v/v).
- 4.2.4. Silica gel, 70-230 mesh, type Merck 7734, with water content standardised at 5 % (w/w/).
- 4.2.5. Glass wool.
- 4.2.6. Acetone for HPLC.
- 4.2.7. Acetonitrile or propionitrile for HPLC.
- 4.2.8. HPLC elution solvent: acetonitrile + acetone (proportions to be adjusted to obtain the desired separation; begin with 50:50 mixture) or propionitrile.
- 4.2.9. Solubilisation solvent: acetone.
- 4.2.10. Reference triglycerides: commercial triglycerides (tripalmitin, triolein, etc.) may be used and the retention times then plotted in accordance with the equivalent carbon number, or alternatively reference chromatograms obtained from soya oil, mixture 30:70 soya oil — olive oil and pure olive oil (see notes 1 and 2 and figures 1 to 4).
- 4.2.11. Solid phase extraction column with silica phase 1 g, 6 ml.

4.3. **Sample preparation**

As a number of interfering substances can give rise to false positive results, the sample must always be purified according to IUPAC method 2.507, used for the determination of polar compounds in frying fats.

4.3.1. *Chromatographic column preparation*

Fill the column (4.1.3) with about 30 ml of elution solvent (4.2.3), then introduce inside the column some glass wool (4.2.5) pushing it to the bottom of the column by means of the glass rod (4.1.5).

In a 100 ml beaker, suspend 25 g of silica gel (4.2.4) in 80 ml of elution mixture (4.2.3), then transfer it to the column by means of a glass funnel (4.1.6).

To ensure the complete transfer of the silica gel to the column, wash the beaker with the elution mixture and transfer the washing portions to the column too.

Open the cock and let the solvent elute from the column until its level is about 1 cm over the silica gel.

4.3.2. *Column chromatography*

Weigh with the accuracy of 0,001 g, $2,5 \pm 0,1$ g of oil, previously filtered, homogenised and anhydriified, if necessary, in a 50 ml volumetric flask (4.1.7).

Dissolve it in about 20 ml of elution solvent (4.2.3). If necessary, slightly heat it to make the dissolution easily. Cool at room temperature and adjust the volume with elution solvent.

By means of a volumetric pipette, introduce 20 ml of solution inside the column prepared according to 4.3.1, open the cock and let the solvent elute to the silica gel layer level.

Then elute with 150 ml of elution solvent (4.2.3), adjusting the solvent rate at about 2 ml/min (150 ml will take about 60-70 minutes to pass through the column).

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The eluate is recovered in a 250 ml round-bottomed flask (4.1.1) previously tared in an oven and exactly weighed. Eliminate the solvent at reduced pressure in a rotary evaporator (4.1.9) and weigh the residue that will be used to prepare the solution for HPLC analysis and for methyl ester preparation.

The sample recovery from the column must be 90 % at least for the extra virgin, virgin, ordinary, refined and olive oil categories, and a minimum of 80 % for lampante and olive-pomace oils.

4.3.3. *SPE purification*

Silica SPE column is activated by passing 6 ml of hexane (4.2.3) under vacuum, avoiding dryness.

Weigh to an accuracy of 0,001 g, 0,12 g in a 2 ml vial (4.1.15) and dissolve with 0,5 ml of hexane (4.2.3).

Load the SPE column with the solution and elute with 10 ml of hexane-diethyl ether (87:13 v/v) (4.2.3) under vacuum.

The collected fraction is evaporated to dryness in a rotary evaporator (4.1.9) under reduced pressure at room temperature. The residue is dissolved in 2 ml of acetone (4.2.6) for triacylglycerol (TAG) analysis.

4.4. **HPLC analysis**

4.4.1. *Preparation of the samples for chromatographic analysis*

A 5 % solution of the sample to be analysed is prepared by weighing $0,5 \pm 0,001$ g of the sample into a 10 ml graduated flask and making up to 10 ml with the solubilisation solvent (4.2.9).

4.4.2. *Procedure*

Set up the chromatographic system. Pump elution solvent (4.2.8) at a rate of 1,5 ml/min to purge the entire system. Wait until a stable base line is obtained.

Inject 10 μ l of the sample prepared as in point 4.3.

4.4.3. *Calculation and expression of results*

Use the area normalisation method, i.e. assume that the sum of the areas of the peaks corresponding to TAGs from ECN 42 up to ECN 52 is equal to 100 %.

Calculate the relative percentage of each triglyceride using the formula:

$$\% \text{ triglyceride} = \text{area of peak} \times 100 / \text{sum of peak areas}$$

.

The results should be given to at least two decimal places.

See notes 1 to 4.

4.5. **Calculation of triacylglycerols composition (moles %) from fatty acid composition data (area %)**

4.5.1. *Determination of fatty acid composition*

Fatty acid composition is determined by ISO 5508 by means of a capillary column. The methyl esters are prepared according to COI/T.20/Doc. No 24.

4.5.2. *Fatty acids for calculation*

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Glycerides are grouped by their Equivalent Carbon Number (ECN), taking into account the following equivalencies between ECN and fatty acids. Only fatty acids with 16 and 18 carbon atoms were taken into consideration, because only these are important for olive oil. The fatty acids should be normalised to 100 %.

Fatty acid (FA)	Abbreviation	Molecular weight(MW)	ECN
Palmitic acid	P	256,4	16
Palmitoleic acid	Po	254,4	14
Stearic acid	S	284,5	18
Oleic acid	O	282,5	16
Linoleic acid	L	280,4	14
Linolenic acid	Ln	278,4	12

4.5.3. Conversion of area % into moles for all fatty acids (1)

$\text{moles P} = \frac{\text{area \% P}}{\text{MW P}}$	$\text{moles S} = \frac{\text{area \% S}}{\text{MW S}}$	$\text{moles Po} = \frac{\text{area \% Po}}{\text{MW Po}}$
$\text{moles O} = \frac{\text{area \% O}}{\text{MW O}}$	$\text{moles L} = \frac{\text{area \% L}}{\text{MW L}}$	$\text{moles Ln} = \frac{\text{area \% Ln}}{\text{MW Ln}}$

4.5.4. Normalisation of fatty acid moles to 100 % (2)

$$\text{moles \% P(1,2,3)} = \frac{\text{moles P} * 100}{\text{moles(P+S+Po+O+L+Ln)}}$$

$$\text{moles \% S(1,2,3)} = \frac{\text{moles S} * 100}{\text{moles(P+S+Po+O+L+Ln)}}$$

$$\text{moles \% Po(1,2,3)} = \frac{\text{moles Po} * 100}{\text{moles(P+S+Po+O+L+Ln)}}$$

$$\text{moles \% O(1,2,3)} = \frac{\text{moles O} * 100}{\text{moles(P+S+Po+O+L+Ln)}}$$

$$\text{moles \% L(1,2,3)} = \frac{\text{moles L} * 100}{\text{moles(P+S+Po+O+L+Ln)}}$$

$$\text{moles \% Ln(1,2,3)} = \frac{\text{moles Ln} * 100}{\text{moles(P+S+Po+O+L+Ln)}}$$

The result gives the percentage of each fatty acid in moles % in the overall (1, 2, 3-) position of the TAGs.

Then the sum of the saturated fatty acids P and S (SFA) and the unsaturated fatty acids Po, O, L and Ln (UFA) are calculated (3):

$$\text{moles \% SFA} = \text{moles \% P} + \text{moles \% S}$$

$$\text{moles \% UFA} = 100 - \text{moles \% SFA}$$

4.5.5. Calculation of the fatty acid composition in 2- and 1, 3- positions of TAGs

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The fatty acids are distributed to three pools as follows: one for 2- position and two identical for 1- and 3- positions, with different coefficients for the saturated (P and S) and unsaturated acids (Po, O, L and Ln).

4.5.5.1. Saturated fatty acids in 2-position [P(2) and S(2)] (4):

$$\text{moles \% P(2)} = \text{moles \% P(1,2,3)} * 0,06$$

$$\text{moles \% S(2)} = \text{moles \% S(1,2,3)} * 0,06$$

4.5.5.2. Unsaturated fatty acids in 2-position [Po(2), O(2), L(2) and Ln(2)] (5):

$$\text{moles \% Po(2)} = \frac{\text{moles \% Po(1,2,3)}}{\text{moles \% UFA}} * (100 - \text{moles \% P(2)} - \text{moles \% S(2)})$$

$$\text{moles \% O(2)} = \frac{\text{moles \% O(1,2,3)}}{\text{moles \% UFA}} * (100 - \text{moles \% P(2)} - \text{moles \% S(2)})$$

$$\text{moles \% L(2)} = \frac{\text{moles \% L(1,2,3)}}{\text{moles \% UFA}} * (100 - \text{moles \% P(2)} - \text{moles \% S(2)})$$

$$\text{moles \% Ln(2)} = \frac{\text{moles \% Ln(1,2,3)}}{\text{moles \% UFA}} * (100 - \text{moles \% P(2)} - \text{moles \% S(2)})$$

4.5.5.3. Fatty acids in 1,3-positions [P(1,3), S(1,3), Po(1,3), O(1,3), L(1,3) and Ln(1,3)] (6):

$$\text{moles \% P(1,3)} = \frac{\text{moles \% P(1,2,3)} - \text{moles \% P(2)}}{2} + \text{moles \% P(1,2,3)}$$

$$\text{moles \% S(1,3)} = \frac{\text{moles \% S(1,2,3)} - \text{moles \% S(2)}}{2} + \text{moles \% S(1,2,3)}$$

$$\text{moles \% Po(1,3)} = \frac{\text{moles \% Po(1,2,3)} - \text{moles \% Po(2)}}{2} + \text{moles \% Po(1,2,3)}$$

$$\text{moles \% O(1,3)} = \frac{\text{moles \% O(1,2,3)} - \text{moles \% O(2)}}{2} + \text{moles \% O(1,2,3)}$$

$$\text{moles \% L(1,3)} = \frac{\text{moles \% L(1,2,3)} - \text{moles \% L(2)}}{2} + \text{moles \% L(1,2,3)}$$

$$\text{moles \% Ln(1,3)} = \frac{\text{moles \% Ln(1,2,3)} - \text{moles \% Ln(2)}}{2} + \text{moles \% Ln(1,2,3)}$$

4.5.6. Calculation of triacylglycerols

4.5.6.1. TAGs with one fatty acid (AAA, here LLL, PoPoPo) (7)

$$\text{moles \% AAA} = \frac{\text{moles \% A(1,3)} * \text{moles \% A(2)} * \text{moles \% A(1,3)}}{10000}$$

4.5.6.2. TAGs with two fatty acids (AAB, here PoPoL, PoLL) (8)

$$\text{moles \% AAB} = \frac{\text{moles \% A(1,3)} * \text{moles \% A(2)} * \text{moles \% B(1,3)} * 2}{10000}$$

$$\text{moles \% ABA} = \frac{\text{moles \% A(1,3)} * \text{moles \% B(2)} * \text{moles \% A(1,3)}}{10000}$$

4.5.6.3. TAGs with three different fatty acids (ABC, here OLLn, PLLn, PoOLn, PPOLn) (9)

$$\text{moles \% ABC} = \frac{\text{moles \% A(1,3)} * \text{moles \% B(2)} * \text{moles \% C(1,3)} * 2}{10000}$$

$$\text{moles \% BCA} = \frac{\text{moles \% B(1,3)} * \text{moles \% C(2)} * \text{moles \% A(1,3)} * 2}{10000}$$

$$\text{moles \% CAB} = \frac{\text{moles \% C(1,3)} * \text{moles \% A(2)} * \text{moles \% B(1,3)} * 2}{10000}$$

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

4.5.6.4. Triacylglycerols with ECN42

The triacylglycerols with ECN42 are calculated according to equations 7, 8 and 9 and are then given in order of expected elution in HPLC (normally only three peaks).

LLL

PoLL and the positional isomer LPoL

OLLn and the positional isomers OLnL and LnOL

PoPoL and the positional isomer PoLPo

PoOLn and the positional isomers OPoLn and OLnPo

PLLn and the positional isomers LLnP and LnPL

PoPoPo

SLnLn and the positional isomer LnSLn

PPoLn and the positional isomers PLnPo and PoPLn

The triacylglycerols with ECN42 are given by the sum of the nine triacylglycerols including their positional isomers. The results should be given to at least two decimal places.

5. EVALUATION OF THE RESULTS

The calculated theoretical content and the content determined by the HPLC analysis are compared. If the difference in the absolute value of the HPLC data minus the theoretical data is greater than the values stated for the appropriate oil category in the standard, the sample contains seed oil.

Results are given to two decimal figures.

6. EXAMPLE (THE NUMBERS REFER TO THE SECTIONS IN THE TEXT OF THE METHOD)

— 4.5.1. Calculation of moles % fatty acids from GLC data (normalised area %)

The following data are obtained for the fatty acid composition by GLC:

FA	P	S	Po	O	L	Ln
MW	256,4	284,5	254,4	282,5	280,4	278,4
Area %	10,0	3,0	1,0	75,0	10,0	1,0

— 4.5.3 Conversion of area % into moles for all fatty acids (see formula (1))

$$\begin{aligned} \text{moles P} &= \frac{10}{256,4} = 0,03900 \text{ moles P} \\ \text{moles S} &= \frac{3}{284,5} = 0,01054 \text{ moles S} \\ \text{moles Po} &= \frac{1}{254,4} = 0,00393 \text{ moles Po} \\ \text{moles O} &= \frac{75}{282,5} = 0,26549 \text{ moles O} \\ \text{moles L} &= \frac{10}{280,4} = 0,03566 \text{ moles L} \\ \text{moles Ln} &= \end{aligned}$$

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

$$\frac{1}{278,4} = 0,00359 \text{ moles Ln}$$
$$\text{Total} = 0,35821 \text{ moles TAGs}$$

— 4.5.4 *Normalisation of fatty acid moles to 100 % (see formula (2))*

$$\text{moles \% P(1,2,3)} = \frac{0,03900 \text{ moles P} * 100}{0,35821 \text{ moles}} = 10,887 \%$$

$$\text{moles \% S(1,2,3)} = \frac{0,01054 \text{ moles S} * 100}{0,35821 \text{ moles}} = 2,942 \%$$

$$\text{moles \% Po(1,2,3)} = \frac{0,00393 \text{ moles Po} * 100}{0,35821 \text{ moles}} = 1,097 \%$$

$$\text{moles \% O(1,2,3)} = \frac{0,26549 \text{ moles O} * 100}{0,35821 \text{ moles}} = 74,116 \%$$

$$\text{moles \% L(1,2,3)} = \frac{0,03566 \text{ moles L} * 100}{0,35821 \text{ moles}} = 9,955 \%$$

$$\text{moles \% Ln(1,2,3)} = \frac{0,00359 \text{ moles Ln} * 100}{0,35821 \text{ moles}} = 1,002 \%$$

$$\text{Total moles \%} = 100 \%$$

Sum of the saturated and unsaturated fatty acids in the 1,2,3-position of TAGs (see formula (3)):

$$\text{moles \% SFA} = 10,887 \% + 2,942 \% = 13,829 \%$$

$$\text{moles \% UFA} = 100,000 \% - 13,829 \% = 86,171 \%$$

— 4.5.5 *Calculation of the fatty acid composition in 2- and 1,3-positions of the TAGs*

— Saturated fatty acids in 2-position [P(2) and S(2)] (see formula (4))

4.5.5.1

$$\text{moles \% P(2)} = 10,887 \% * 0,06 = 0,653 \text{ moles \%}$$

$$\text{moles \% S(2)} = 2,942 \% * 0,06 = 0,177 \text{ moles \%}$$

— Unsaturated fatty acids in 2-position [Po(1,3), O(1,3), L(1,3) and Ln(1,3)] (see formula 4.5.5.2 (5))

$$\text{moles \% Po(2)} = \frac{1,097 \%}{86,171 \%} * (100 - 0,653 - 0,177) = 1,262 \text{ moles \%}$$

$$\text{moles \% O(2)} = \frac{74,116 \%}{86,171 \%} * (100 - 0,653 - 0,177) = 85,296 \text{ moles \%}$$

$$\text{moles \% L(2)} = \frac{9,955 \%}{86,171 \%} * (100 - 0,653 - 0,177) = 11,457 \text{ moles \%}$$

$$\text{moles \% Ln(2)} = \frac{1,002 \%}{86,171 \%} * (100 - 0,653 - 0,177) = 1,153 \text{ moles \%}$$

— Fatty acids in 1,3-positions [P(1,3), S(1,3), Po(1,3), O(1,3), L(1,3) and Ln(1,3)] (see 4.5.5.3 formula (6))

$$\text{moles \% P(1,3)} = \frac{10,887 - 0,653}{2} + 10,887 = 16,004 \text{ moles \%}$$

$$\text{moles \% S(1,3)} = \frac{2,942 - 0,177}{2} + 2,942 = 4,325 \text{ moles \%}$$

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

$$\text{moles \% Po(1,3)} = \frac{1,097 - 1,262}{2} + 1,097 = 1,015 \text{ moles \%}$$

$$\text{moles \% O(1,3)} = \frac{74,116 - 85,296}{2} + 74,116 = 68,526 \text{ moles \%}$$

$$\text{moles \% L(1,3)} = \frac{9,955 - 11,457}{2} + 9,955 = 9,204 \text{ moles \%}$$

$$\text{moles \% Ln(1,3)} = \frac{1,002 - 1,153}{2} + 1,002 = 0,927 \text{ moles \%}$$

— 4.5.6. Calculation of triacylglycerols

From the calculated fatty acid composition in sn-2- and sn-1,3-positions:

FA in	1,3-pos	2-pos
P	16,004 %	0,653 %
S	4,325 %	0,177 %
Po	1,015 %	1,262 %
O	68,526 %	85,296 %
L	9,204 %	11,457 %
Ln	0,927 %	1,153 %
Sum	100,0 %	100,0 %

the following triacylglycerols are calculated:

LLL

PoPoPo

PoLL with 1 positional isomer

SLLn with 1 positional isomer

PoPoL with 1 positional isomer

PPoLn with 2 positional isomers

OLLn with 2 positional isomers

PLLn with 2 positional isomers

PoOLn with 2 positional isomers

— TAGs with one fatty acid (LLL, PoPoPo) (see formula (7))

4.5.6.1.

$$\text{mol \% LLL} = \frac{9,204 \% * 11,457 \% * 9,204 \%}{10000}$$

$$, = 0,09706 \text{ mol LLL}$$

$$\text{mol \% PoPoPo} = \frac{1,015 \% * 1,262 \% * 1,015 \%}{10000} = 0,00013 \text{ mol PoPoPo}$$

— TAGs with two fatty acids (PoLL, SLLn, PoPoL) (see formula (8))

4.5.6.2

$$\text{mol \% PoLL} + \text{LLPo} = \frac{1,015 \% * 11,457 \% * 9,204 \% * 2}{10000} = 0,02141$$

$$\text{mol \% LPoL} = \frac{9,204 \% * 1,262 \% * 9,204 \%}{10000} = 0,01069$$

$$0,03210 \text{ mol PoLL}$$

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

$$\text{mol \% SLnLn} + \text{LnLnS} = \frac{4,325 \% * 1,153 \% * 0,927 \% * 2}{10000} = 0,00092$$

$$\text{mol \% LnSLn} = \frac{0,927 \% * 0,177 \% * 0,927 \%}{10000} = 0,00002$$

0,00094 mol SLnLn

$$\text{mol \% PoPoL} + \text{LPoPo} = \frac{1,015 \% * 1,262 \% * 9,204 \% * 2}{10000} = 0,00236$$

$$\text{mol \% PoLPo} = \frac{1,015 \% * 11,457 \% * 1,015 \%}{10000} = 0,00118$$

0,00354 mol PoPoL

— TAGs with three different fatty acids (PoPLn, OLLn, PLLn, PoOLn) See formula (9)
4.5.6.3

$$\text{mol \% PPOln} = \frac{16,004 \% * 1,262 \% * 0,927 \% * 2}{10000} = 0,00374$$

$$\text{mol \% LnPPo} = \frac{0,927 \% * 0,653 \% * 1,015 \% * 2}{10000} = 0,00012$$

$$\text{mol \% PoLnP} = \frac{1,015 \% * 1,153 \% * 16,004 \% * 2}{10000} = 0,00375$$

0,00761 mol PPOln

$$\text{mol \% OLLn} = \frac{68,526 \% * 11,457 \% * 0,927 \% * 2}{10000} = 0,14556$$

$$\text{mol \% LnOL} = \frac{0,927 \% * 85,296 \% * 9,204 \% * 2}{10000} = 0,14555$$

$$\text{mol \% LLnO} = \frac{9,204 \% * 1,153 \% * 68,526 \% * 2}{10000} = 0,14544$$

0,43655 mol OLLn

$$\text{mol \% PLLn} = \frac{16,004 \% * 11,457 \% * 0,927 \% * 2}{10000} = 0,03399$$

$$\text{mol \% LnPL} = \frac{0,927 \% * 0,653 \% * 9,204 \% * 2}{10000} = 0,00111$$

$$\text{mol \% LLnP} = \frac{9,204 \% * 1,153 \% * 16,004 \% * 2}{10000} = 0,03397$$

0,06907 mol PLLn

$$\text{mol \% PoOLn} = \frac{1,015 \% * 85,296 \% * 0,927 \% * 2}{10000} = 0,01605$$

$$\text{mol \% LnPoO} = \frac{0,927 \% * 1,262 \% * 68,526 \% * 2}{10000} = 0,01603$$

$$\text{mol \% OLnPo} = \frac{68,526 \% * 1,153 \% * 1,015 \% * 2}{10000} = 0,01604$$

0,04812 mol PoOLn

ECN42 = 0,69512 mol TAGs

Note 1: The elution order can be determined by calculating the equivalent carbon numbers, often defined by the relation

$$\text{ECN} = \text{CN} - 2n$$

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

, where CN is the carbon number and n is the number of double bonds; it can be calculated more precisely by taking into account the origin of the double bond. If n_o , n_l and n_{ln} are the numbers of double bonds attributed to oleic, linoleic and linolenic acids respectively, the equivalent carbon number can be calculated by means of the relation of the formula:

$$EN = CN - d_o n_o - d_l n_l - d_{ln} n_{ln}$$

where the coefficient d_o , d_l and d_{ln} can be calculated by means of the reference triglycerides. Under the conditions specified in this method, the relation obtained will be close to:

$$ECN = CN - (2,60 n_o) - (2,35 n_l) - (2,17 n_{ln})$$

Note 2: With several reference triglycerides, it is also possible to calculate the resolution with respect to triolein:

$$\alpha = RT^1 / RT$$

triolein

by use of the reduced retention time

$$RT^1 = RT - RT \text{ solvent}$$

The graph of $\log \alpha$ against f (number of double bonds) enables the retention values to be determined for all the triglycerides of fatty acids contained in the reference triglycerides — see Figure 1.

Note 3: The efficiency of the column should permit clear separation of the peak of trilinolein from the peaks of the triglycerides with an adjacent RT. The elution is carried out up to ECN 52 peak.

Note 4: A correct measure of the areas of all peaks of interest for the present determination is ensured if the second peak corresponding to ECN 50 is 50 % of full scale of the recorder.

Figure 1

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Graph of $\log \alpha$ against f (number of double bonds)

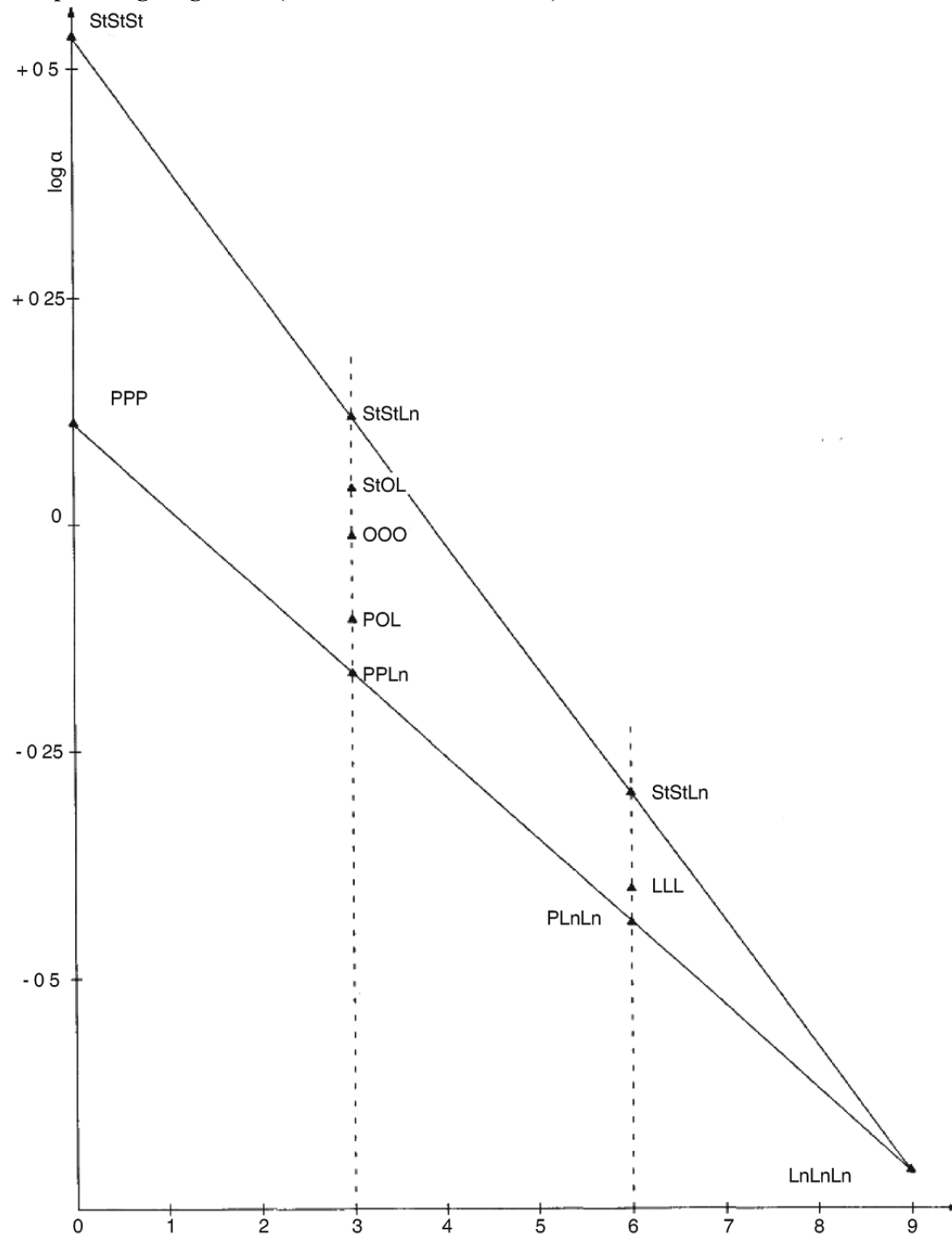



Figure 2

Low linoleic olive oil

(a)

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

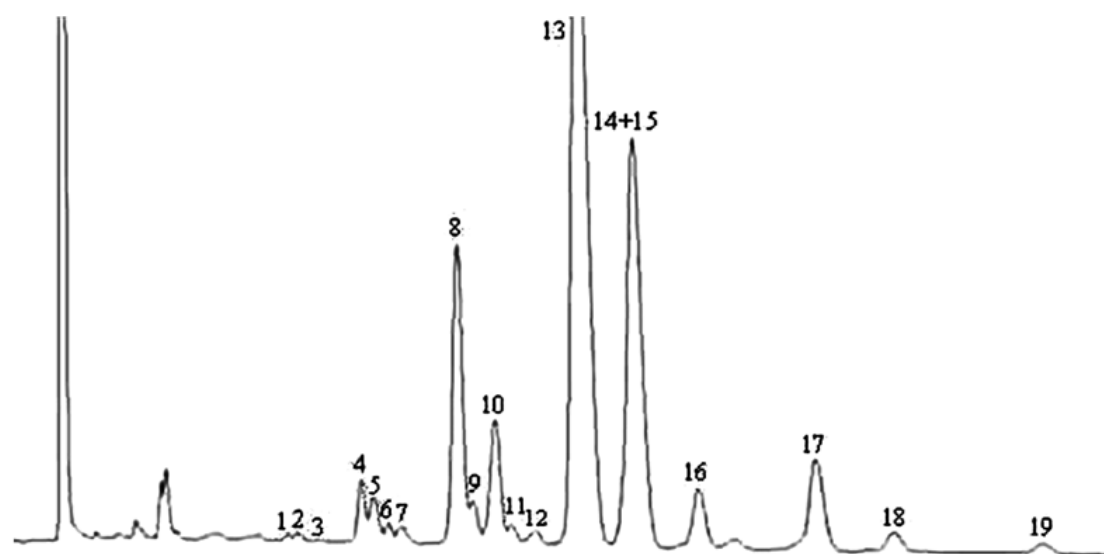

(b)

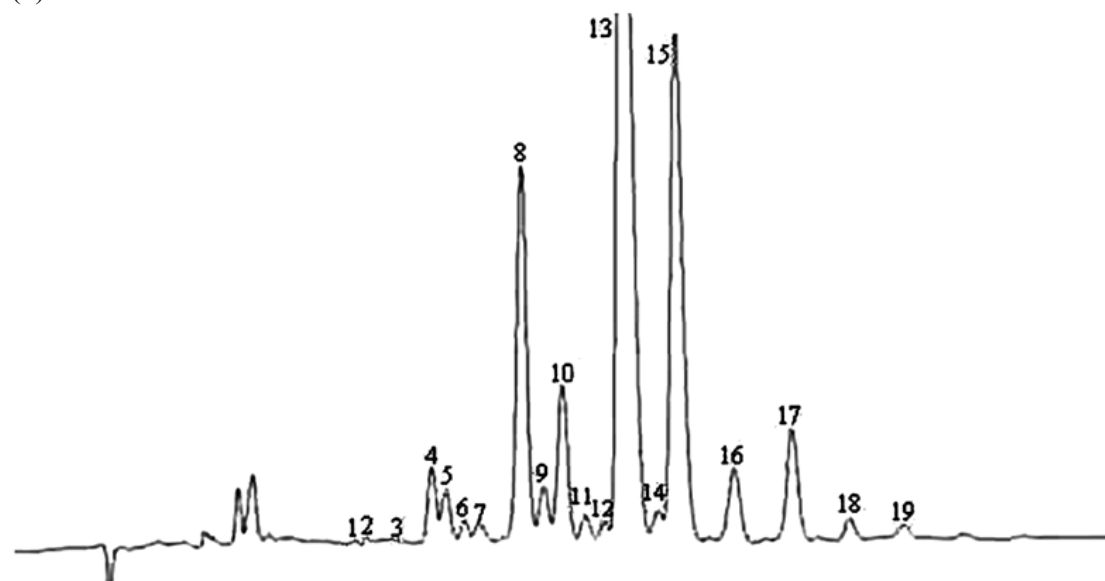



Figure 3

High linoleic olive oil


(a)

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

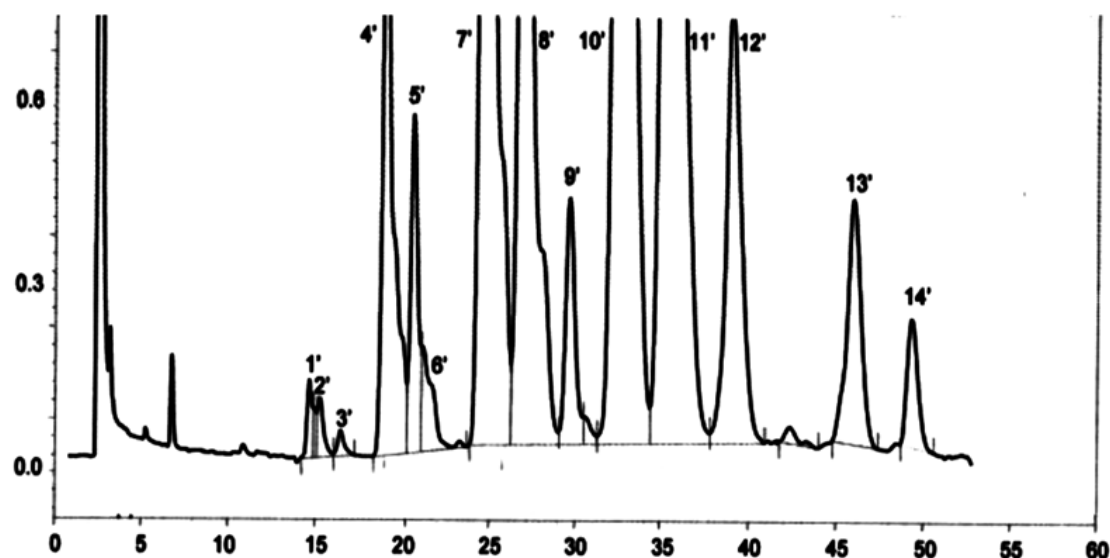

(b)

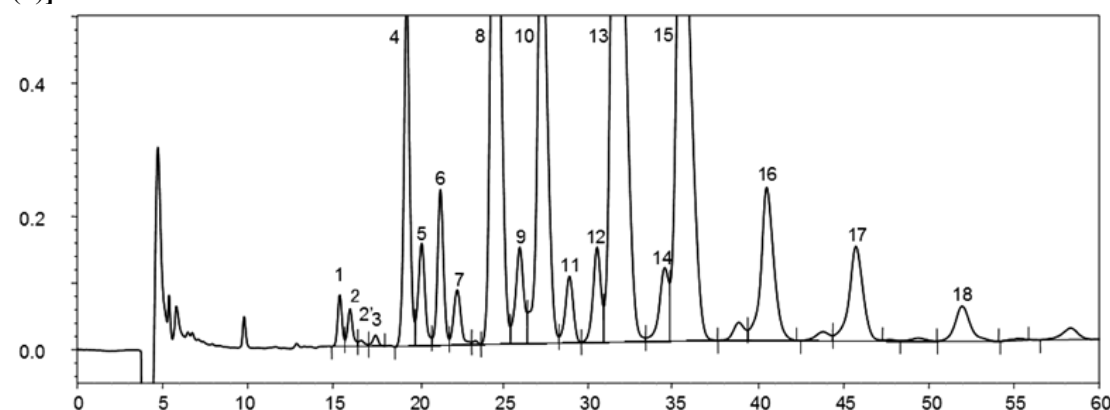

[^{F19}ANNEX XIX

DETERMINATION OF ALIPHATIC ALCOHOLS CONTENT BY CAPILLARY GAS CHROMATOGRAPHY

Textual Amendments

F19 Inserted by Commission Regulation (EC) No 796/2002 of 6 May 2002 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis and the additional notes in the Annex to Council Regulation (EEC) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.

1. OBJECT

The procedure describes a method for the determination of aliphatic alcohols content in oils and fats.

2. PRINCIPLE OF THE METHOD

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The fatty substance, with 1-eicosanol added as internal standard, is saponified with ethanolic potassium hydroxide and then the unsaponifiable matter extracted with ethyl ether. The alcoholic fraction is separated from the unsaponifiable matter by chromatography on a basic silica gel plate; the alcohols recovered from the silica gel are transformed into trimethylsilyl ethers and analysed by capillary gas chromatography.

3. EQUIPMENT

- 3.1. 250 ml round-bottomed flask fitted with a reflux condenser having ground-glass joints.
- 3.2. 500 ml separating funnel.
- 3.3. 250 ml round-bottomed flasks.
- 3.4. Chromatographic tank for thin-layer chromatographic analysis, for glass plates of dimensions 20 x 20 cm.
- 3.5. Ultraviolet lamp having a wavelength of 366 or 254 nm.
- 3.6. 100 µl and 500 µl microsyringes.
- 3.7. A cylindrical filter funnel with a G3 porous septum (porosity 15 to 40 µm) of diameter approximately 2 cm and a depth of some 5 cm, with an attachment suitable for filtration under vacuum and a 12/21 male ground glass joint.
- 3.8. 50 ml vacuum conical flask with a 12/21 ground-glass female joint which can be fitted to the filter funnel (3.7).
- 3.9. A 10 ml test tube with a tapering bottom and a sealing stopper.
- 3.10. Gas chromatograph for use with a capillary column, and provided with a splitting system composed of:
 - 3.10.1. Thermostatic chamber for columns (column oven) to hold the temperature desired with a precision of ± 1 °C.
 - 3.10.2. A temperature-adjustable injection unit with a persilanised glass vaporising element.
 - 3.10.3. A flame ionisation detector and converter-amplifier.
 - 3.10.4. Recorder-integrator for operation with the converter-amplifier (3.10.3), with response time not exceeding one second and with variable paper-speed.
- 3.11. Glass or fused silica capillary column, of length 20 to 30 m, internal diameter 0,25 to 0,32 mm, with SE-52 or SE-54 liquid phase or equivalent, with a film thickness between 0,10 and 0,30 µm.
- 3.12. Microsyringe for gas chromatography, of 10 µl capacity with hardened needle.
- 3.13. Analytical balance sensitive to 1 mg (with 0,1 mg display).

4. REAGENTS

- 4.1. Potassium hydroxide, approximately 2 N ethanolic solution: 130 g potassium hydroxide (minimum concentration 85 %) is dissolved, with cooling, in 200 ml distilled water and then made up to one litre with ethanol. The solution should be stored in a well-stoppered opaque glass bottle.
- 4.2. Ethyl ether, pure for analysis.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 4.3. Anhydrous sodium sulphate, analytical purity.
- 4.4. Glass plates coated with silica gel, without fluorescence indicator, thickness 0,25 mm (commercially available ready for use).
- 4.5. Potassium hydroxide, approximately 0,2 N ethanolic solution; 13 g of potassium hydroxide are dissolved in 20 ml of distilled water and made up to one litre with ethanol.
- 4.6. Benzene, for chromatography (see 5.2.2).
- 4.7. Acetone, for chromatography (See 5.2.2).
- 4.8. Hexane, for chromatography (see 5.2.2).
- 4.9. Ethyl ether, for chromatography (see 5.2.2).
- 4.10. Chloroform, for chromatography.
- 4.11. Reference solution for thin-layer chromatography: cholesterol or phytosterols, 5 % solution in chloroform.
- 4.12. 0,2 % solution of 2', 7'-dichlorofluorescein in ethanol. Make slightly basic by adding a few drops of 2 N alcoholic potassium hydroxide solution.
- 4.13. Anhydrous pyridine, for chromatography.
- 4.14. Hexamethyl disilazane.
- 4.15. Trimethylchlorosilane.
- 4.16. Standard solutions of trimethylsilyl ethers of aliphatic alcohols from C₂₀ to C₂₈. They may be prepared from mixtures of pure alcohols at the time they are required for use.
- 4.17. A 0,1 % (m/v) solution of 1-eicosanol in chloroform (internal standard).
- 4.18. Carrier gas: hydrogen or helium, gas-chromatographic purity.
- 4.19. Auxiliary gas: nitrogen, gas-chromatographic purity.

5. PROCEDURE

5.1. Preparation of the unsaponifiables

- 5.1.1. Using a 500 µl microsyringe place, into a 250 ml round-bottom flask, a volume of 0,1 % 1-eicosanol solution in chloroform (4.17) containing a quantity of 1-eicosanol approximately equal to 10 % of the aliphatic alcohols content in that portion of sample to be taken for analysis. For example, to 5 g of sample add 250 µl of the 0,1 % 1-eicosanol solution if olive oil and 1 500 µl if olive pomace oil.

Evaporate to dryness in current of nitrogen and then weigh accurately 5 g of the dry filtered sample into the same flask.

- 5.1.2. Add 50 ml of 2 N potassium hydroxide ethanolic solution, fit the reflux condenser and heat the apparatus to slight boiling on a steam bath, stirring continuously throughout the heating process until saponification has taken place (the solution becomes clear). Continue heating for a further 20 minutes and then add 50 ml of distilled water through the condenser. The condenser is then disconnected and the flask cooled to approximately 30 °C.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- 5.1.3. The contents of the flask are quantitatively transferred to a separating funnel of 500 ml capacity by adding distilled water several times, using a total of around 50 ml distilled water. Add approximately 80 ml of ethyl ether, shake vigorously for approximately 30 seconds and allow to settle (Note 1).

Separate off the lower aqueous phase collecting it in a second separating funnel. Two further extractions are effected on the aqueous phase, in the same manner, using each time 60 to 70 ml ethyl ether.

Note 1:

Emulsions may be eliminated by adding, using as a spray, small quantities of ethyl alcohol or methyl alcohol.

- 5.1.4. The ethyl ether extracts are combined in a separating funnel and washed with distilled water (50 ml at a time) until the washing water gives a neutral reaction.

Discard the washing water, dry with anhydrous sodium sulphate and filter, into a flask of 250 ml capacity which has been weighed beforehand, the funnel and filter being washed with small quantities of ethyl ether which are added to the total.

- 5.1.5. Distil the ether down to a few ml, then bring to dryness under a slight vacuum or in a current of nitrogen, completing drying in an oven at 100 °C for approximately a quarter of an hour, and then weigh after cooling in a desiccator.

5.2. Separation of alcoholic fractions

- 5.2.1. Preparation of basic TLC plates: the silica gel plates (4.4) are immersed completely, in 0,2 N potassium hydroxide solution (4.5) for 10 seconds, and then left to dry for two hours under an extractor hood and finally placed in an oven at 100 °C for one hour.

Remove from the oven and keep in a calcium chloride desiccator until required for use (plates treated in this way must be used within 15 days).

Note 2:

When basic silica gel plates are used to separate the alcoholic fraction there is no need to treat the unsaponifiables with alumina. It follows that all acid compounds (fatty acids and others) are retained at the origin thereby obtaining both aliphatic alcohol and terpenic alcohol bands which are both separated distinctly from the sterol band.

- 5.2.2. Place a 65/35 by volume hexane/ethyl ether mixture in the plate-developing chamber to a depth of approximately 1 cm⁽⁹⁾.

Close the chamber using an appropriate cover and leave for half an hour to allow equilibration between vapour and liquid. Strips of filter paper dipping into the eluent may be affixed to the inside surfaces of the tank to reduce the development time by approximately one third and obtain more uniform, regular elution of the components.

Note 3:

The developing solution must be replaced for each analysis in order to obtain reproducible developing conditions.

- 5.2.3. An approximately 5 % solution of unsaponifiable matter (5.1.5) in chloroform is prepared and 0,3 ml of the solution is streaked as a uniform strip of minimum thickness, using the 100 µl microsyringe, on a TLC plate at approximately 2 cm from the bottom of the TLC plate. Aligned with the origin, 2 to 3 µl of the aliphatic alcohols

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

reference solution (4.11) are spotted for the identification of the aliphatic alcohols band after development has been completed.

- 5.2.4. Place the plate inside the development tank as stated in 5.2.2. The ambient temperature should be maintained between 15 and 20 °C. Immediately close the chamber with the cover and allow to elute until the solvent front reaches approximately 1 cm from the upper edge of the plate.

The plate is then removed from the development chamber and the solvent evaporated under a hot air current or the plate is left for a while under the extractor hood.

- 5.2.5. The plate is sprayed lightly and evenly with the solution of 2', 7'-dichlorofluorescein when the plate is observed under ultra violet light. The aliphatic alcohols band can be identified through being aligned with the stain obtained from the reference solution: mark the limits of the band with a black pencil; outlining the band of aliphatic alcohols and the band immediately above that, which is the terpenic alcohols band, together.

Note 4:

The aliphatic alcohols band and the terpenic alcohols band are to be grouped together in view of the possible migration of some aliphatic alcohols into the triterpenic alcohols band.

- 5.2.6. Using a metal spatula scrape off the silica gel in the marked area. Place the finely comminuted material removed into the filter funnel (3.7). Add 10 ml of hot chloroform, mix carefully with the metal spatula and filter under vacuum, collecting the filtrate in the conical flask (3.8) attached to the filter funnel.

Wash the pomace in the flask three times with ethyl ether (approximately 10 ml each time) collecting the filtrate in the same flask attached to the funnel. Evaporate the filtrate to a volume of 4 to 5 ml, transfer the residual solution to the previously weighed 10 ml test tube (3.9), evaporate to dryness by mild heating in a gentle flow of nitrogen, make up again using a few drops of acetone, evaporate again to dryness, place in an oven at 105 °C for approximately 10 minutes and then allow to cool in a desiccator and weigh.

The pomace inside the test tube is composed of the alcoholic fraction.

- 5.3. Preparation of the trimethylsilyl ethers

- 5.3.1. The reagent for silylation, consisting of a mixture of 9:3:1 by volume (Note 5) of pyridine-hexamethyldisilazane-trimethylchlorosilane in the proportion of 50 µl for each milligram of aliphatic alcohols, is added to the test tube containing the alcoholic fraction, avoiding all absorption of moisture (Note 6).

Note 5:

Solutions which are ready for use are available commercially. Other silanising reagents such as, for example, bis-trimethylsilyl, trifluor acetamide + 1 % trimethyl chlorosilane, which has to be diluted with an equal volume of anhydrous pyridine, are also available.

Note 6:

The slight opalescence which may form is normal and does not cause any interference. The formation of a white floc or the appearance of a pink colour are indicative of the presence of moisture or deterioration of the reagent. If these occur the test must be repeated.

- 5.3.2. Stopper the test tube, shake carefully (without overturning) until the aliphatic alcohols are completely dissolved. Stand for at least 15 minutes at ambient temperature and then centrifuge for a few minutes. The clear solution is ready for gas chromatographic analysis.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

5.4. Gas chromatography analysis

5.4.1. Preliminary operations, column packing

5.4.1.1. Fit the column in the gas chromatograph, attaching the inlet end to the injector connected to the splitting system and the outlet end to the detector. Carry out a general check of the gas chromatography assembly (tightness of gas fittings, efficiency of the detector, efficiency of the splitting system and of the recording system, etc.).

5.4.1.2. If the column is being used for the first time it is recommended that it should be subjected to conditioning. A little carrier gas is passed through the capillary column and then the gas chromatography assembly is switched on and gradually heated until a temperature not less than 20 °C above the operating temperature (see Note 7) is attained. That temperature is held for not less than two hours and then the assembly is brought to the operating conditions (regulation of gas flow, split flame ignition, connection to the electronic recorder, adjustment of the temperature of the capillary column oven, the detector and the injector, etc.) and the signal is adjusted to a sensitivity not less than twice the highest level contemplated for the execution of the analysis. The course of the base line must be linear, without peaks of any kind, and must not drift. A negative straight-line drift indicates leakage from the column connections; a positive drift indicates inadequate conditioning of the column.

Note 7:

The conditioning temperature shall be at least 20 °C less than the maximum temperature contemplated for the liquid phase employed.

5.4.2. Choice of operating conditions

5.4.2.1. The guideline operating conditions are as follows:

- column temperature: the initial isotherm is set at 180 °C for eight minutes and then programmed at 5 °C/minute to 260 °C and a further 15 minutes at 260 °C,
- temperature of evaporator: 280 °C,
- temperature of detector: 290 °C,
- linear velocity of carrier gas: helium 20 to 35 cm/s, hydrogen 30 to 50 cm/s,
- splitting ratio: 1:50 to 1:100,
- sensitivity of instrument: 4 to 16 times the minimum attenuation,
- sensitivity of recording: 1 to 2 mV fs,
- paper speed: 30 to 60 cm/h,
- quantity of substance injected: 0,5 to 1 µl of TMSE solution.

The above conditions may be modified according to the characteristics of the column and of the gas chromatograph to obtain chromatograms satisfying the following conditions:

- alcohol C₂₆ retention time shall be 18 ± 5 minutes,
- the alcohol C₂₂ peak shall be 80 ± 20 % of the full scale value for olive oil and 40 ± 20 % of the full scale value for seed oil.

5.4.2.2. The above requirements are checked by repeated injection of the standard TMSE mixture of alcohols and the operating conditions are adjusted to yield the best possible results.

5.4.2.3. The parameters for the integration of peaks shall be set so that a correct appraisal of the areas of the peaks considered is obtained.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

5.4.3. Analytical procedure

5.4.3.1. Using the microsyringe of 10 µl capacity draw in 1 µl of hexane followed by 0,5 µl of air and subsequently 0,5 to 1 µl of the sample solution; raise the plunger of the syringe further so the needle is emptied. Push the needle through the membrane of the injection unit and after one to two seconds inject rapidly, then slowly remove the needle after some five seconds.

5.4.3.2. Recording is effected until the TMSE of the aliphatic alcohols present have been eluted completely. The base line shall always correspond to the requirements of 5.4.1.2.

5.4.4. Peak identification

The identification of individual peaks is effected according to the retention times and by comparison with the standard TMSE mixture, analysed under the same conditions.

A chromatogram of the alcoholic fraction of a virgin olive oil is shown in Figure 1.

5.4.5. Quantitative evaluation

5.4.5.1. The peak areas of 1-eicosanol and of the aliphatic alcohols C₂₂, C₂₄, C₂₆ and C₂₈ are calculated by electronic integration.

5.4.5.2. The contents of each aliphatic alcohol, expressed in mg/

1 000

g fatty substance, are calculated as follows:

$$\text{alcohol x} = \frac{A_x \cdot m_s \cdot 1\,000}{A_s \cdot m}$$

where:

A_x = area of the alcohol peak x

A_s = area of 1-eicosanol

m_s = mass of 1-eicosanol in milligrams

m = mass of sample drawn for determination, in grams.

6. EXPRESSION OF THE RESULTS

The contents of the individual aliphatic alcohols in mg/1 000 g of fatty substance and the sum of the 'total aliphatic alcohols' are reported.

Status: Point in time view as at 01/03/2014.


Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

APPENDIX

Determination of the linear velocity of the gas

1 to 3 μl of methane or propane are injected into the gas chromatograph set at normal operating conditions and the time taken for the methane or propane to flow through the column from the instant of injection to the instant the peak elutes (t_M) is measured using a stop clock.

The linear velocity in cm/s is given by L/t_M , where L is the length of the column in centimetres and t_M is the measured time in seconds.]

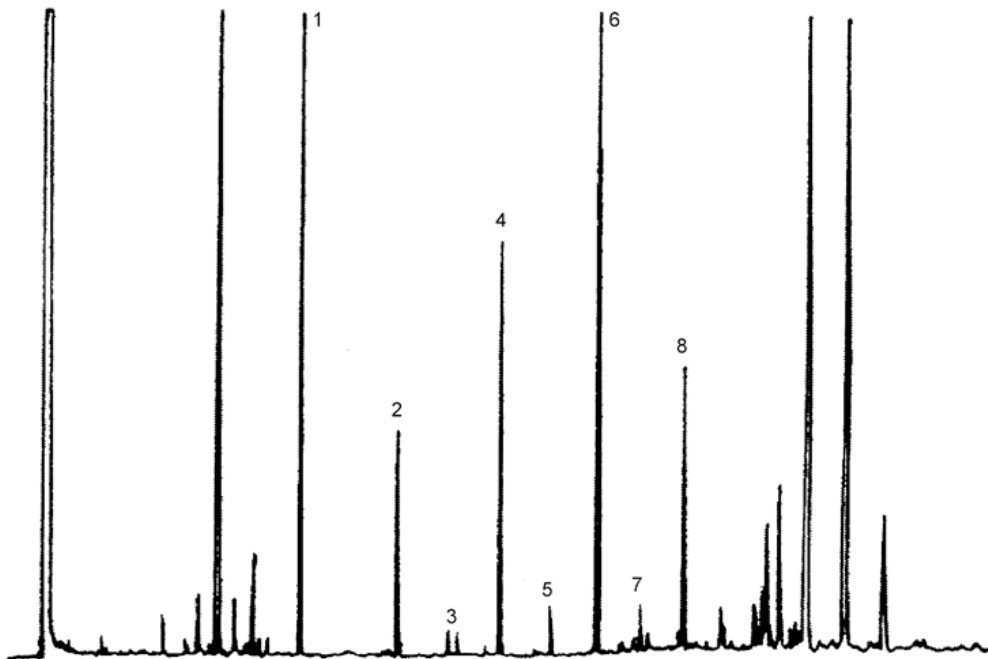

[^{F14}ANNEX XX

Method for the determination of the content of waxes, fatty acid methyl esters and fatty acid ethyl esters by capillary gas chromatography

1. PURPOSE

This method is for the determination of the content of waxes, fatty acid methyl and ethyl esters in olive oils. The individual waxes and alkyl esters are separated according to the number of carbon atoms. The method is recommended as a tool for distinguishing between olive oil and olive-pomace oil and as a quality parameter for extra virgin olive oils enabling the detection of fraudulent mixtures of extra virgin olive oils with lower quality oils whether they are virgin, lampante or some deodorised oils.

2. PRINCIPLE

Addition of suitable internal standards to the oil and fractionation by chromatography on a hydrated silica gel column. Recovery of the fraction eluted under the test conditions

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

(with a lower polarity than that of the triacylglycerols) and direct analysis by capillary gas chromatography.

3. APPARATUS

3.1. **Erlenmeyer flask, 25 ml.**

3.2. **Glass column** for liquid chromatography, internal diameter 15 mm, length 30-40 cm, fitted with a suitable stopcock.

3.3. **Gas chromatograph** suitable for use with a capillary column, equipped with a system for direct, on-column injection comprising:

3.3.1. **Thermostat-controlled oven with temperature programming.**

3.3.2. **Cold injector** for direct on-column injection

3.3.3. **Flame ionisation detector and converter-amplifier.**

3.3.4. **Recorder-integrator** (Note 1) for use with the converter-amplifier (point 3.3.3), with a response time of not more than 1 s and a variable paper speed.

Note 1: Computerised systems may also be used where the gas chromatography data are entered through a PC.

3.3.5. **Capillary column, fused silica (for analysis of the waxes and methyl and ethyl esters)**, length 8-12 m, internal diameter 0,25-0,32 mm, internally coated with liquid phase (Note 2) to a uniform thickness of 0,10-0,30 μm .

Note 2: Suitable commercial liquid phases are available for this purpose such as SE52, SE54, etc.

3.4. **Microsyringe**, 10 μl , with hardened needle, for direct on-column injection.

3.5. **Electric shaker.**

3.6. **Rotary evaporator.**

3.7. **Muffle oven.**

3.8. **Analytical balance** for weighing to an accuracy of $\pm 0,1$ mg.

3.9. Usual laboratory glassware.

4. REAGENTS

4.1. **Silica gel**, 60-200 μm mesh. Place the silica gel in the muffle oven at 500 $^{\circ}\text{C}$ for at least 4 h. Allow to cool and then add 2 % water in relation to the quantity of silica gel used. Shake well to homogenise slurry and keep in the desiccator for at least 12 h prior to use.

4.2. **n-hexane**, chromatography grade or residue grade (the purity must be checked).

WARNING – Fumes may ignite. Keep away from sources of heat, sparks or naked flames. Make sure the bottles are always properly closed. Ensure proper ventilation during usage. Avoid build-up of fumes and remove any possible fire risk, such as heaters or electric apparatus not manufactured from non-inflammable material. Pernicious if inhaled, because it may cause nerve cell damage. Avoid breathing in the fumes. Use a suitable respiratory apparatus if necessary. Avoid contact with eyes and skin.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

4.3. Ethyl ether, chromatography grade

WARNING – Highly inflammable and moderately toxic. Irritates the skin. Pernicious if inhaled. May cause damage to eyes. Effects may be delayed. It can form explosive peroxides. Fumes may ignite. Keep away from sources of heat, sparks or naked flames. Make sure the bottles are always properly closed. Ensure proper ventilation during usage. Avoid build-up of fumes and remove any possible fire risk, such as heaters or electric apparatus not manufactured from non-inflammable material. Do not evaporate to dryness or near dryness. The addition of water or an appropriate reducing agent can reduce peroxide formation. Do not drink. Avoid breathing in the fumes. Avoid prolonged or repeated contact with skin.

4.4. n-heptane, chromatography grade, or iso-octane

WARNING – Inflammable. Pernicious if inhaled. Keep away from sources of heat, sparks or naked flames. Make sure the bottles are always properly closed. Ensure proper ventilation during usage. Avoid breathing in the fumes. Avoid prolonged or repeated contact with skin.

4.5. Standard solution of lauryl arachidate (Note 3), at 0,05 % (m/V) in heptane (internal standard for waxes).

Note 3: Palmityl palmitate, myristyl stearate or arachidyl laureate may also be used.

4.6. Standard solution of methyl heptadecanoate, at 0,02 % (m/V) in heptane (internal standard for methyl and ethyl esters).

4.7. Sudan 1 (1-phenylazo-2-naphthol).

4.8. Carrier gas: hydrogen or helium, pure, gas chromatography grade.

WARNING

Hydrogen. Highly inflammable, under pressure. Keep away from sources of heat, sparks, naked flames or electric apparatus not manufactured from non-inflammable material. Make sure the bottle valve is shut when not in use. Always use with a pressure reducer. Release the tension of the reducer spring before opening the bottle valve. Do not stand in front of the bottle outlet when opening the valve. Ensure proper ventilation during usage. Do not transfer hydrogen from one bottle to another. Do not mix gas in the bottle. Make sure the bottles cannot be knocked over. Keep them away from sunlight and sources of heat. Store in a corrosive-free environment. Do not use damaged or unlabelled bottles.

Helium. Compressed gas at high pressure. It reduces the amount of oxygen available for breathing. Keep the bottle shut. Ensure proper ventilation during usage. Do not enter storage areas unless they are properly ventilated. Always use with a pressure reducer. Release the tension of the reducer spring before opening the bottle valve. Do not transfer gas from one bottle to another. Make sure the bottles cannot be knocked over. Do not stand in front of the bottle outlet when opening the valve. Keep them away from sunlight and sources of heat. Store in a corrosive-free environment. Do not use damaged or unlabelled bottles. Do not inhale. Use solely for technical purposes.

4.9. Auxiliary gases:

— Hydrogen, pure, gas chromatography grade.

— Air, pure, gas chromatography grade.

WARNING

Air. Compressed gas at high pressure. Use with caution in the presence of combustible substances as the self-ignition temperature of most of the organic compounds in the air is considerably lower under high pressure. Make sure the bottle valve is shut when not in use.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Always use a pressure reducer. Release the tension of the reducer spring before opening the bottle valve. Do not stand in front of the bottle outlet when opening the valve. Do not transfer gas from one bottle to another. Do not mix gas in the bottle. Make sure the bottles cannot be knocked over. Keep them away from sunlight and sources of heat. Store in a corrosive-free environment. Do not use damaged or unlabelled bottles. Air intended for technical purposes must not be used for inhaling or respiratory apparatus.

5. PROCEDURE

5.1. Preparation of the chromatography column

Suspend 15 g of silica gel (point 4.1) in n-hexane (point 4.2) and introduce into the column (point 3.2). Allow to settle spontaneously. Complete settling with the aid of an electric shaker to make the chromatographic bed more homogeneous. Percolate 30 ml of n-hexane to remove any impurities. Weigh exactly about 500 mg of the sample into the 25-ml flask (point 3.1), using the analytical balance (point 3.8), and add a suitable amount of internal standard (point 4.5), depending on the assumed wax content, e.g. add 0,1 mg of lauryl arachidate in the case of olive oil, 0,25-0,50 mg in the case of olive-pomace oil and 0,05 mg of methyl heptadecanoate for olive oils (point 4.6).

Transfer the prepared sample to the chromatography column with the aid of two 2-ml portions of n-hexane (point 4.2).

Allow the solvent to flow to 1 mm above the upper level of the absorbent. Percolate a further of n-hexane/ethyl ether (99:1) and collect 220 ml at a flow of about 15 drops every 10 seconds. **(This fraction contains the methyl and ethyl esters and waxes).** (Note 4) (Note 5).

Note 4: The n-hexane/ethyl ether (99:1) mixture should be freshly prepared every day

Note 5: 100 µl of Sudan I dye at 1 % in the elution mixture can be added to the sample solution to check visually that the waxes are eluted properly.

The retention time of the dye lies in between that of the waxes and triacylglycerols. Hence, when the dye reaches the bottom of the chromatography column, elution has to be suspended because all the waxes have been eluted.

Evaporate the resultant fractions in a rotary evaporator until the solvent is almost removed. Remove the last 2 ml under a weak current of nitrogen. Collect the fraction containing the methyl and ethyl esters is diluted with 2-4 ml of n-heptane or iso-octane.

5.2. Gas chromatography analysis

5.2.1. Preliminary procedure

Fit the column to the gas chromatograph (point 3.3), connecting the inlet port to the on-column system and the outlet port to the detector. Check the gas chromatography apparatus (operation of gas loops, efficiency of detector and recorder system, etc.).

If the column is being used for the first time, it is advisable to condition it. Run a light flow of gas through the column, then switch on the gas chromatography apparatus. Gradually heat until a temperature of 350 °C is reached after approximately 4 h.

Maintain this temperature for at least 2 h, then regulate the apparatus to the operating conditions (regulate gas flow, light flame, connect to electronic recorder (point 3.3.4), regulate oven temperature for column, regulate detector, etc.). Record the signal at a sensitivity at least twice as high as that required for the analysis. The base line should be linear, with no peaks of any kind, and must not have any drift.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Negative straight-line drift indicates that the column connections are not correct while positive drift indicates that the column has not been properly conditioned.

5.2.2. Choice of operating conditions for waxes and methyl and ethyl esters (Note 6).

The operating conditions are generally as follows:

—Column temperature	: 20 °C/min 5 °C/min
—Detector temperature	: 80 °C at first (1') ————— 140 °C ————— 335 °C (20)
—Amount injected	: 350 °C.
—Carrier gas	: 1 µl of n-heptane solution (2-4 ml).
—Instrument sensitivity	: helium or hydrogen at the optimal linear speed for the gas chosen (see Appendix A).
	: suitable for fulfilling the above conditions.

Note 6: Due to the high final temperature, positive drift is allowed but may not exceed more than 10 % of the full-scale value.

These conditions may be modified to suit the characteristics of the column and the gas chromatograph in order to separate all the waxes and fatty acid methyl and ethyl esters and to obtain satisfactory peak separation (see Figures 2, 3 and 4) and a retention time of 18 ± 3 minutes for the lauryl arachidate internal standard. The most representative peak of the waxes must be over 60 % of the full-scale value while the methyl heptadecanoate internal standard for the methyl and ethyl esters must reach the full-scale value.

The peak integration parameters should be determined in such a way as to obtain a correct evaluation of the peak areas considered.

5.3. Performance of the analysis

Take up 10 µl of the solution with the aid of the 10 µl micro-syringe, drawing back the plunger until the needle is empty. Introduce the needle into the injection system and inject quickly after 1–2 s. After about 5 s, gently extract the needle.

Perform the recording until the waxes or stigmastadienes are completely eluted, depending on the fraction being analysed.

The base line must always meet the required conditions.

5.4. Peak identification

Identify the peaks from the retention times by comparing them with mixtures of waxes with known retention times, analysed under the same conditions. The alkyl esters are identified from mixtures of methyl and ethyl esters of the chief fatty acids in olive oils (palmitic and oleic).

Figure 1 provides a chromatogram of the waxes in a virgin olive oil. Figures 2 and 3 show the chromatograms of two retail extra virgin olive oils, one with methyl and ethyl esters and the other without them. Figure 4 gives the chromatograms for a top-quality extra virgin olive oil and the same oil spiked with 20 % deodorised oil.

5.5. Quantitative analysis of the waxes

Determine the area of the peaks corresponding to the lauryl arachidate internal standard and the aliphatic esters from C₄₀ to C₄₆ with the aid of the integrator.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Determine the total waxes content by adding each individual wax, in mg/kg of fat, as follows:

$$\text{Waxes, mg / kg} = \frac{(\sum A_x) \times m_s \times 1000}{A_s \times m}$$

where:

- A_x = area corresponding to the peak for the individual ester, in computer counts
 A_s = area corresponding to the peak for the lauryl arachidate internal standard, in computer counts
 m_s = mass of the lauryl arachidate internal standard added, in milligrams;
 m = mass of the sample taken for determination, in grams.

5.5.1. Quantitative analysis of the methyl and ethyl esters

With the aid of the integrator, determine the areas of the peaks corresponding to the methyl heptadecanoate internal standard, the methyl esters of the C_{16} and C_{18} fatty acids and the ethyl esters of the C_{16} and C_{18} fatty acids.

Determine the content of each alkyl ester, in mg/kg of fat, as follows:

$$\text{Ester, mg / kg} = \frac{A_x \times m_s \times 1000}{A_s \times m}$$

where:

- A_x = area corresponding to the peak for the individual C_{16} and C_{18} ester, in computer counts
 A_s = area corresponding to the peak for the methyl heptadecanoate internal standard, in computer counts
 m_s = mass of the methyl heptadecanoate internal standard added, in milligrams;
 m = mass of the sample taken for determination, in grams.

6. EXPRESSION OF RESULTS

Report the sum of the contents of the different waxes from C_{40} to C_{46} (*Note 7*) in milligrams per kilograms of fat.

Report the sum of the contents of the methyl esters and ethyl esters from C_{16} to C_{18} and the total of the two.

Results should be expressed to the nearest mg/kg.

Note 7: The components for quantification refer to the peaks with even carbon numbers amongst the C_{40} - C_{46} esters, according to the specimen chromatogram of the waxes in olive oil provided in the attached figure. For identification purposes, if the C_{46} ester is split, it is recommended to analyse the wax fraction of an olive-pomace oil where the C_{46} peak is distinguishable because it is clearly predominant.

Report the ratio between ethyl esters and methyl esters

Figure 1

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Example of a gas chromatogram of the wax fraction of an olive oil⁽¹⁰⁾

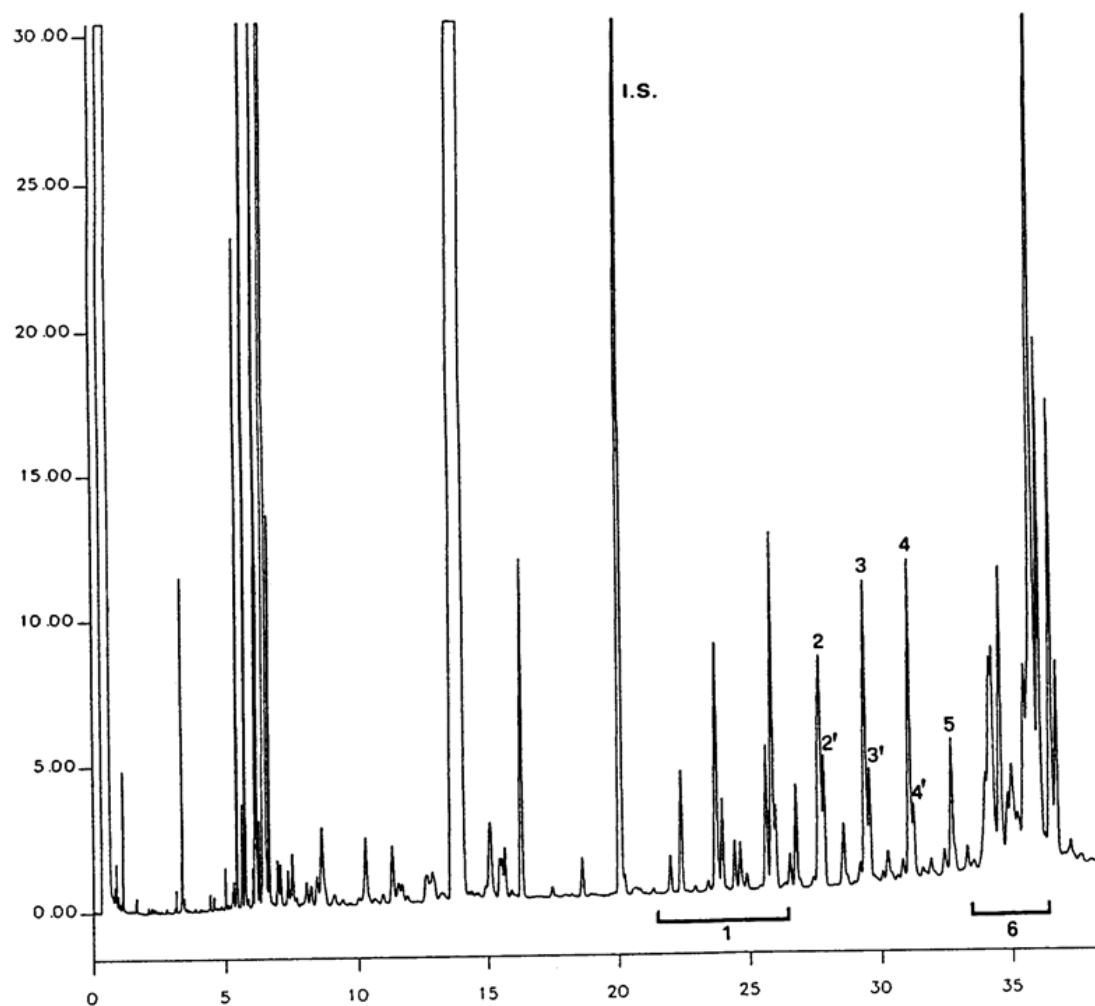


Figure 2

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Methyl esters, ethyl esters and waxes in a virgin olive oil

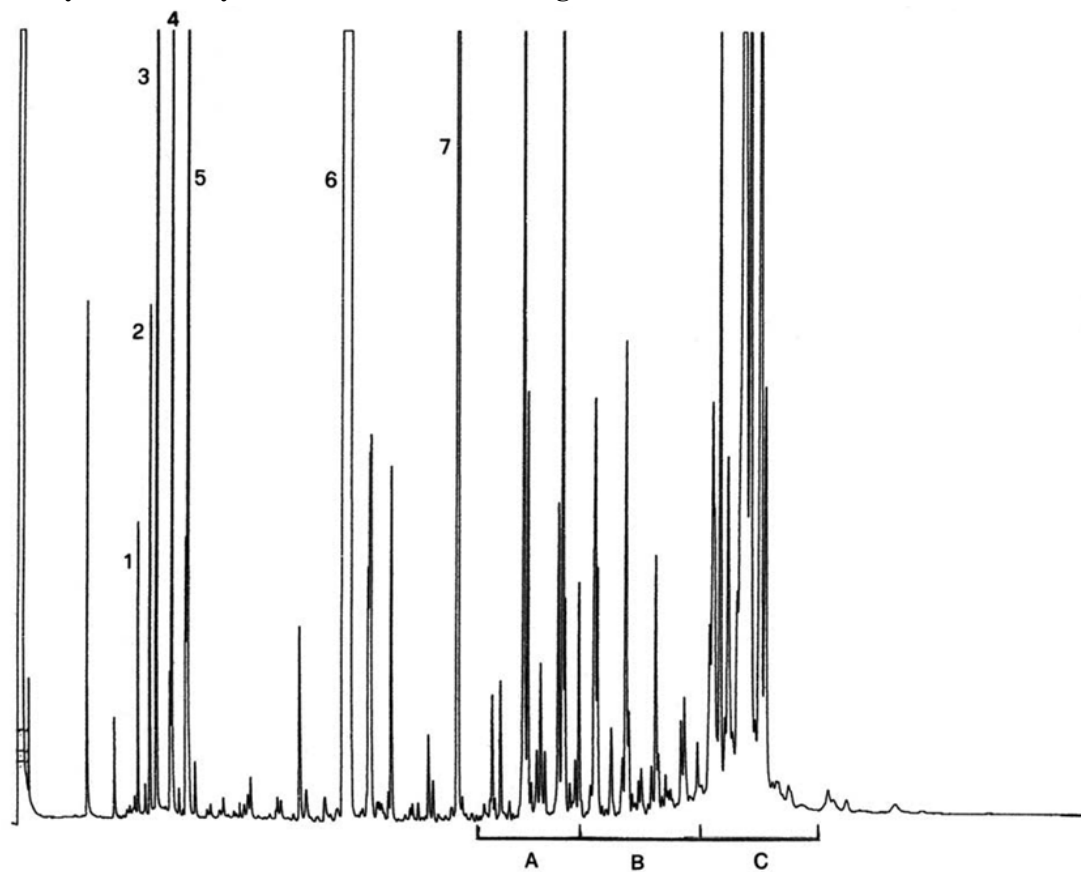


Figure 3

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Methyl esters, ethyl esters and waxes in an extra virgin olive oil


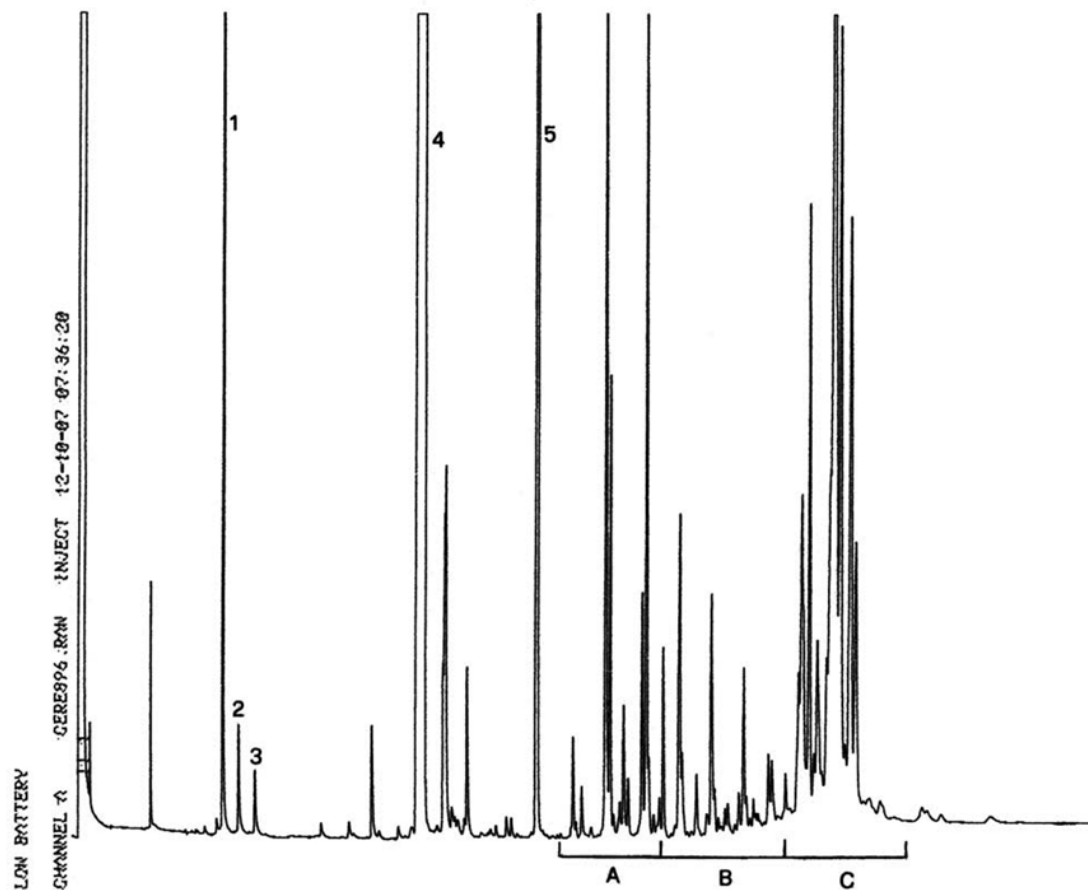


Figure 4

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Part of a chromatogram of an extra virgin olive oil and the same oil spiked with deodorised oil

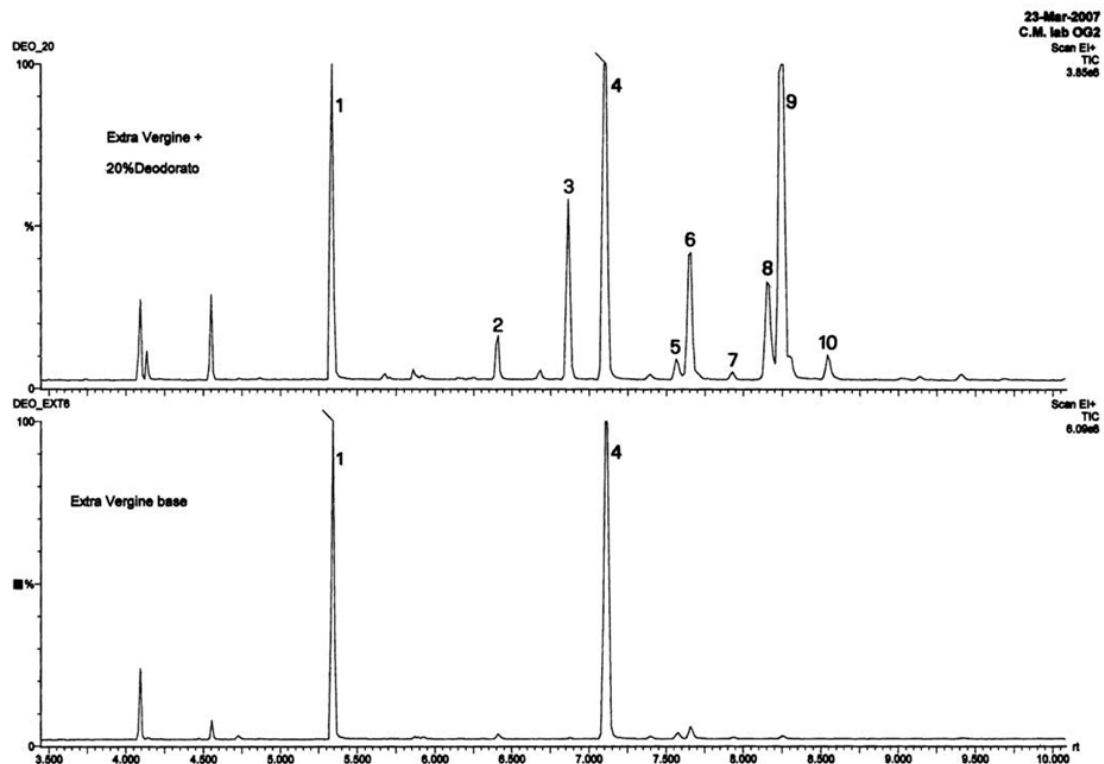

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Appendix A

Determination of linear gas speed

Inject 1:3 µl of methane (or propane) into the gas chromatograph after adjusting it to the normal operating conditions. Measure the time the gas takes to run through the column from the moment it is injected until the peak emerges (tM).

The linear speed in cm/s is given by L/tM where L is the length of the column, in cm, and tM is the time measured in s.]

[^{F15}ANNEX XXa

METHOD FOR THE DETECTION OF EXTRANEEOUS OILS IN OLIVE OILS

1. SCOPE

This method is used to detect the presence of extraneous vegetable oils in olive oils. High linoleic vegetable oils (soybean, rapeseed, sunflower, etc.), and some high oleic vegetable oils - such as hazelnut, high oleic sunflower and olive-pomace oils - can be detected in olive oils. The level detected depends on the type of extraneous oil and the variety of olive. For hazelnut oil, a detection level between 5 and 15 % is common. The method is unable to identify the type of extraneous oil detected, and only indicates if the olive oil is genuine or non-genuine.

2. PRINCIPLE

The oil is purified by solid phase extraction (SPE) on silica gel cartridges. The triacylglycerol (TAG) composition is determined by reverse phase high resolution liquid chromatography using a refractive index detector and propionitrile as the mobile phase. Fatty acid methyl esters (FAMES) are prepared from purified oil by methylation with a cold solution of KOH in methanol (Annex X B) and then the esters are analysed by capillary gas chromatography using high polar columns (Annex X A). The theoretical triacylglycerol composition is calculated from the fatty acid composition by a computer program assuming a 1,3-random, 2-random distribution of fatty acids in the triacylglycerol, with restrictions for saturated fatty acids in the 2-position. The calculation method is a modification of the procedure described in Annex XVIII. Several mathematical algorithms are calculated from theoretical and experimental (HPLC) triacylglycerol compositions, and the resulting values are compared with those contained in a database built from genuine olive oils.

3. MATERIAL AND REAGENTS

3.1. Oil purification

- 3.1.1. 25-ml conical flasks.
- 3.1.2. 5-ml screw top glass tubes and caps fitted with PTFE joint.
- 3.1.3. Silica gel cartridges, 1 g (6 ml), for solid phase extraction (for example, Waters, Massachusetts, USA).
- 3.1.4. *n*-hexane, analytical grade.
- 3.1.5. Solvent mixture of hexane/diethyl ether (87:13, v/v).
- 3.1.6. *N*-heptane, analytical grade.
- 3.1.7. Acetone, analytical grade.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

3.2. HPLC analysis of triacylglycerols

- 3.2.1. Micro syringes (50 µL) and needles for HPLC injection.
- 3.2.2. Propionitrile, super purity or HPLC grade (for example, ROMIL, Cambridge, United Kingdom), used as mobile phase.
- 3.2.3. HPLC column (25 cm × 4 mm internal diameter), packed with RP-18 phase (4 µm particle size).

3.3. Preparation of fatty acid methyl esters

(See Annex X B)

- 3.3.1. Methanol containing not more than 0,5 % water.
- 3.3.2. Heptane, analytical grade.
- 3.3.3. A 2N solution of potassium hydroxide in methanol. Dissolve 1,1 g of potassium hydroxide in 10 ml of methanol.
- 3.3.4. 5-ml screw top glass tubes and caps provided with PTFE joint.

3.4. GC analysis of FAMES

(See method for the determination of *trans*-unsaturated fatty acids by capillary column gas chromatography set out in Annex X A).

- 3.4.1. Micro syringes (5 µL) and needles for GC injection.
- 3.4.2. Hydrogen or helium as carrier gas.
- 3.4.3. Hydrogen and oxygen for FID detector.
- 3.4.4. Nitrogen or helium as auxiliary carrier gas.
- 3.4.5. Fused silica capillary column (50-60 m × 0,25 – 0,30 mm internal diameter) coated with cyanopropylpolysiloxane or cyanopropylphenylsiloxane phases (SP-2380 or similar) with 0,20-0,25 µm of film thickness.

4. APPARATUS

- 4.1. Vacuum apparatus for solid phase extraction.
- 4.2. Rotary evaporator.
- 4.3. HPLC equipment composed of:
 - 4.3.1. Degasser for the mobile phase.
 - 4.3.2. Rheodyne injector valve with a 10 µL loop.
 - 4.3.3. High pressure pump unit.
 - 4.3.4. Thermostatic oven for the HPLC column capable of maintaining sub-ambient temperatures (15-20 °C), (for example, Peltier type).
 - 4.3.5. Refractive index detector.
 - 4.3.6. Computerised data acquisition system provided with an integration program.

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

4.4 Capillary gas chromatography equipment described in Annex X A, provided with:

4.4.1. Split injector.

4.4.2. Flame ionisation detector (FID).

4.4.3. Oven with programmable temperature.

4.4.4. Computerised data acquisition system provided with an integration program.

4.5. Computer with Microsoft EXCEL program.

5. ANALYTICAL PROCEDURE

5.1. Oil purification

An SPE silica gel cartridge is placed in a vacuum elution apparatus and washed under vacuum with 6 ml of hexane. The vacuum is released to prevent the column from drying and a conical flask is placed under the cartridge. A solution of the oil (0,12 g, approximately) in 0,5 ml of hexane is loaded into the column and the solution is pulled through and then eluted with 10 ml of the solvent mixture (3.1.5) of hexane-diethyl ether (87:13 v/v) under vacuum. The eluted solvent is homogenised and approximately half of the volume is poured into another conical flask. Both solutions are separately evaporated to dryness in a rotary evaporator under reduced pressure at room temperature. For triacylglycerol analysis, one of the residues is dissolved in 1 ml of acetone (See first paragraph of point 5.2) and poured into a 5-ml screw top glass tube. The other residue is dissolved in 1 ml of *n*-heptane and poured into a second 5-ml screw top glass tube for preparing the fatty acid methyl esters.

Note: Oil purification may be done using a silica gel column, as described in IUPAC method 2.507.

5.2. HPLC analysis of triacylglycerols

Set up the HPLC system, maintaining the column temperature at 20 °C and using propionitrile as the mobile phase at a flow rate of 0,6 ml/min. When the baseline is stable run a solvent injection; if the base line appears disturbed in the region from 12 to 25 min, use another type of acetone or a mixture of propionitrile/acetone (25:75) to dissolve the sample.

Note: Some types of acetone produce disturbances of the baseline in the above-mentioned region.

Inject a 10 µl aliquot of the solution of purified oil in acetone (5 %). The run takes approximately 60 min. Oven temperature and/or flow rate must be adjusted to achieve a chromatogram similar to that depicted in Figure 1 where trilinolein (peak 1) elutes at 15,5 min and the resolutions between the pairs LLL/OLLn (peaks 1 and 2) and OLL/OOLn (peaks 4 and 5) are good.

The height of peak 2 (OLLn+PoLL) must reach at least 3 % of the full scale.

5.3. Preparation of fatty acid methyl esters

Add 0,1 mL of a 2N solution of potassium hydroxide in methanol to the solution of purified oil in 1 mL of *n*-heptane. Cap the tube and screw tight. Shake the tube vigorously for 15 seconds and leave to stratify until the upper layer becomes clear (5 minutes). The *n*-heptane solution is ready to be injected into the gas chromatograph. The solution may be left at room temperature for a maximum of 12 hours.

5.4. GC analysis of fatty acid methyl esters

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

The procedure described in the method for the determination of *trans*-unsaturated fatty acids must be used (see Annex X A).

The GC system is set up at an oven temperature of 165 °C. The recommended oven temperature is isothermal at 165 °C for 10 min, then raising it to 200 °C at 1,5 °C/min. An injector temperature between 220 °C and 250 °C is recommended to minimise the formation of *trans*-fatty acids (see Annex X A). Detector temperature 250 °C. Hydrogen or helium must be used as the carrier gas at a column head pressure of 130 kPa, approximately. Injection volume 1 µL in split injection mode.

A GC profile similar to that shown in Figure 2 must be obtained. Special attention must be paid to the resolution between C18:3 and C20:1 (the C18:3 peak must appear before the C20:1). To achieve these conditions, the initial temperature and/or the column head pressure must be optimised. Adjust the injector conditions (temperature, split ratio and volume injection) to minimise the discrimination of palmitic and palmitoleic acid.

The height of the C20:0 peak must be about 20 % of full scale to quantify the *trans* isomers. If the C18:0 peak appears distorted, reduce the sample amount.

6. INTEGRATION OF CHROMATOGRAPHIC PEAKS

6.1. HPLC chromatogram

Figure 1 shows a typical HPLC chromatogram of the triacylglycerols of a purified olive oil. For peak integration, three baselines must be traced: the first between the start of peak 1 and the end of peak 3; the second between the start of peak 4 and the valley before peak 8; the third between the valley preceding peak 8 and the end of peak 18.

The total area is the sum of the areas of all the peaks (identified and not identified) from peak 1 to peak 18. The percentage of each peak is given by

$$TAG_x(\%) = 100(A_x + A_T)$$

The percentages have to be given to two decimal figures.

6.2. GC chromatogram

Figure 2 shows a GC chromatogram of fatty acid alkyl esters obtained from a purified olive oil. Percentages of the following fatty acids must be calculated:

Palmitic;	P (C16:0)	=	methyl ester + ethyl ester
Stearic;	S (C18:0)	=	methyl ester
Palmitoleic;	Po (C16:1)	=	sum of methyl esters of the two <i>cis</i> -isomers
Oleic;	O (C18:1)	=	sum of methyl esters of the two <i>cis</i> -isomers + ethyl ester + <i>trans</i> -isomers
Linoleic;	L (C18:2)	=	methyl ester+ ethyl ester + <i>trans</i> -isomers
Linolenic;	Ln (C18:3)	=	methyl ester + <i>trans</i> -isomers

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

Arachidic;	A (C20:0)	=	methyl ester
Eicosenoic (gondoic);	G (C20:1)	=	methyl ester

Ethyl and *trans*-isomers esters may be absent in the GC chromatogram.

Total area (AT) is the sum of all the peaks appearing in the chromatogram from C14:0 to C24:0, except that corresponding to squalene. The percentage of each peak is calculated as follows:

$$FA_x(\%) = 100(A_x + A_T)$$

The results have to be expressed to two decimal places.

For the calculations of the computer programs, it is not necessary to normalise to 100 because this is done automatically.

Figure 1

HPLC chromatogram of TAGs of a 'Chamlali' virgin olive oil. Main components of chromatographic peaks

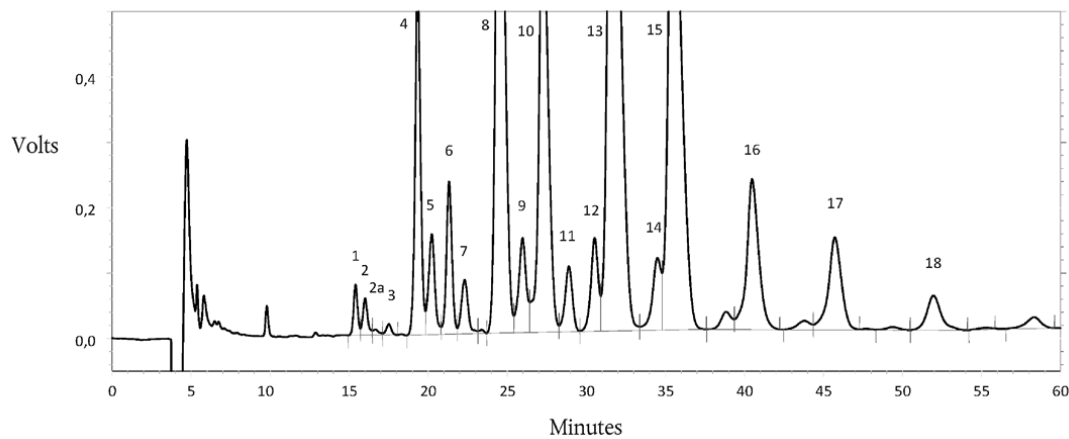


TABLE 1

Repeatability data of the determination of virgin olive oil TAGs by HPLC at a column temperature of 20 °C and using propionitrile as mobile phase

ECN	HPLC TAGs peaks	Sample 1		Sample 2		Sample 3		Sample 4		Sample 5		
		Mean (%)	RSD _r (%)	Mean (%)	RSD _r (%)	Mean (%)	RSD _r (%)	Mean (%)	RSD _r (%)	Mean (%)	RSD _r (%)	
42	1	LLL	0,02	7,23	0,066	5,18	0,095	4,1	0,113	0,95	0,34	1,05
	2	OLLn + PoLL	0,085	7,44	0,24	1,78	0,26	2,25	0,35	2,02	0,5	2,83
	3	PLLn	0,023	15,74	0,039	5,51	0,057	5,62	0,082	4,35	0,12	6,15
44	4	OLL	0,47	1,52	1,53	0,42	2,62	0,98	3,35	1,05	4,37	1,13

n = 3 replicates

RSD_r = Relative Standard Deviation of the repeatability

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

	5	OOLn + PoOL	1,07	2,01	1,54	0,46	1,61	0,71	1,72	1,07	1,77	2,4
	6	PLL + PoPoO	0,11	12,86	0,24	4,37	0,65	1,32	1,35	0,73	2,28	1,24
	7	POLn + PpoPo + PpoL	0,42	5,11	0,49	2,89	0,55	2,01	0,85	1,83	1,09	1,96
46	8	OOL + LnPP	6,72	0,63	8,79	0,31	11,21	0,42	13,25	0,33	15,24	0,23
	9	PoOO	1,24	2,86	1,49	0,95	1,63	0,85	2,12	0,45	2,52	0,56
	10	SLL + PLO	2,7	0,65	4,05	0,7	6,02	0,65	9,86	0,53	11,53	0,31
	11	PoOP + SpoL + SOLn + SpoPo	0,64	4,42	0,69	3,02	0,79	1,23	1,53	0,89	1,7	1,66
48	12+13	OOO + PLP+ PoPP	49,6	0,07	48,15	0,06	42,93	0,06	33,25	0,1	24,16	0,06
	14	SOL	0,82	1,72	0,92	1,56	1,05	1,32	1,25	1,05	1,6	1,77
	15	POO	22,75	0,25	21,8	0,2	21,05	0,3	20,36	0,35	20,17	0,14
50	16	POP	3,05	0,46	4,56	0,42	4,98	0,52	5,26	0,41	5,57	0,38
	17	SOO	6,87	0,21	5,56	0,33	4,86	0,43	4,12	0,72	3,09	0,69
	18	POS + SLS	1,73	1,23	1,65	1,1	1,54	0,99	1,49	1,1	1,41	1,0


n = 3 replicates
RSD_r = Relative Standard Deviation of the repeatability

Figure 2

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

GC chromatogram of fatty acid alkyl esters obtained from an olive-pomace oil by transesterification with a cold solution of KOH in methanol

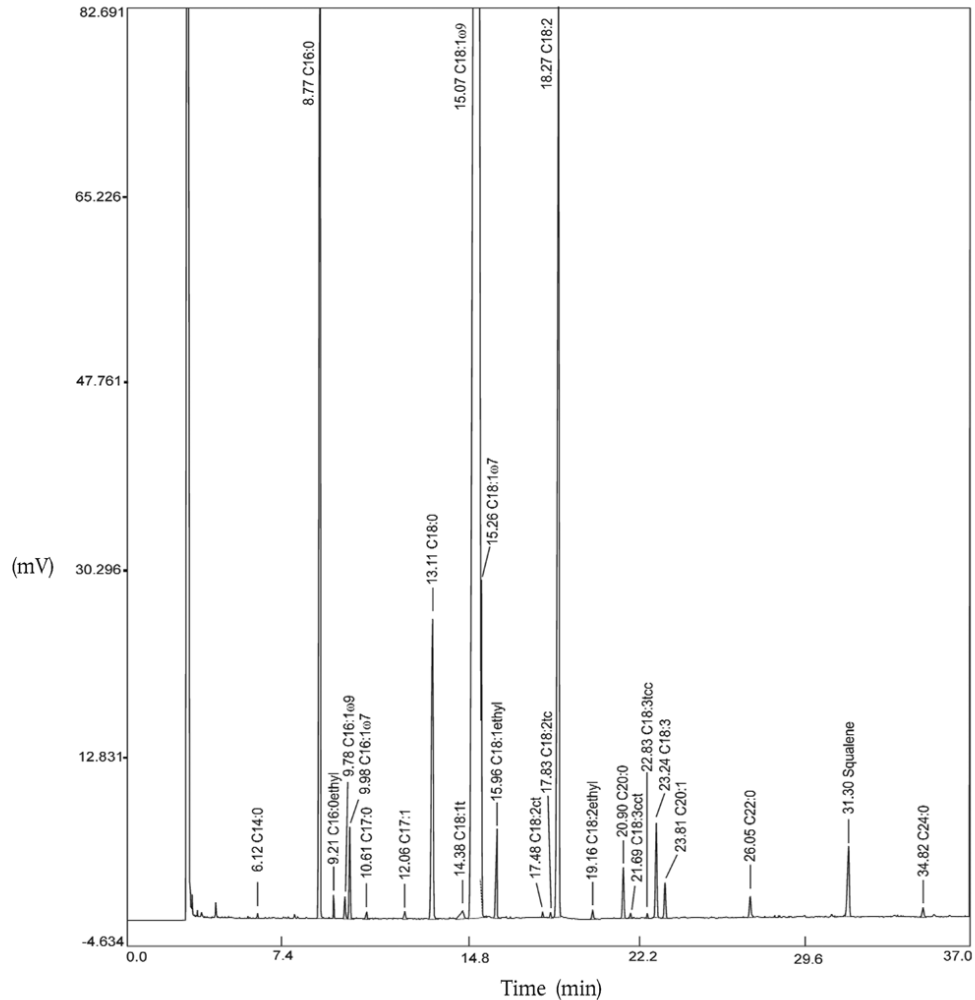

7. DETECTION OF EXTRANEEOUS OILS IN OLIVE OILS

The calculation method of the detection of extraneous oils in olive oils by means of a comparison of mathematical algorithms with a data base built from genuine olive oils is set out in Annex I to standard IOC/T.20/Doc. No 25.]

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

[^{F7} ANNEX XXI

RESULTS OF CONFORMITY CHECKS CARRIED OUT
 ON OLIVE OILS REFERRED TO IN ARTICLE 8(2)

Sample	Category	Country of origin	Place of inspection	Labelling				Chemical parameters	Organoleptic characteristics			Final conclusion
				Legal name of origin	Designation of origin	Standard information	Legibility		Medicinal	Flavour	Appearance	

a Internal market (mill, bottlers, retail stage), export, import.

b Each characteristic of olive oil set out in Annex I shall have a code.

c Conform/not conform.

d Not required for olive oil and pomace-oil.]

Status: Point in time view as at 01/03/2014.

Changes to legislation: There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations. (See end of Document for details)

- (1) OJ No 172, 30.9.1966, p. 3025/66.
- (2) OJ No L 353, 17.12.1990, p. 23.
- (3) OJ No L 128, 24.5.1977, p. 6.
- (4) OJ No L 166, 1.7.1988, p. 10.
- (5) [^{F3}OJ L 228, 1.9.2009, p. 3.]
- (6) [^{F2}Directive 2011/91/EU of the European Parliament and of the Council of 13 December 2011 on indications or marks identifying the lot to which a foodstuff belongs (OJ L 334, 16.12.2011, p. 1).]
- (7) [^{F10}After elution of the sterol esters the chromatogram trace must not show any significant peaks (triglycerides).]
- (8) [^{F2}They may refrain from tasting an oil when they notice any extremely intense negative attribute by direct olfactory means, in which case they shall record this exceptional circumstance in the profile sheet.]
- (9) [^{F19}In these cases in particular, a 95/5 by volume benzene/acetone eluent mixture must be used to obtain distinct band separation.]
- (10) [^{F14}After elution of the sterol esters, the chromatogram should not show any significant peaks (triacylglycerols).]

Textual Amendments

- F2** Substituted by Commission Implementing Regulation (EU) No 1348/2013 of 16 December 2013 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F3** Substituted by Commission Implementing Regulation (EU) No 299/2013 of 26 March 2013 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F10** Substituted by Commission Regulation (EC) No 702/2007 of 21 June 2007 amending Commission Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F14** Inserted by Commission Regulation (EU) No 61/2011 of 24 January 2011 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant methods of analysis.
- F19** Inserted by Commission Regulation (EC) No 796/2002 of 6 May 2002 amending Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-pomace oil and on the relevant methods of analysis and the additional notes in the Annex to Council Regulation (EEC) No 2658/87 on the tariff and statistical nomenclature and on the Common Customs Tariff.

Status:

Point in time view as at 01/03/2014.

Changes to legislation:

There are outstanding changes not yet made to Commission Regulation (EEC) No 2568/91. Any changes that have already been made to the legislation appear in the content and are referenced with annotations.